



# Diode Laser Modules and Pattern Generators



2011 Catalog

Superior Reliability & Performance

# How to Contact Us

## Contact Us by Phone

Coherent, Inc.  
(800) 527-3786 or (408) 764-4983

Benelux: +31 (30) 280 6060

China: +86 (10) 6280 0209

France: +33 (0)1 6985 5145

Germany: +49 (6071) 968333

Italy: +39 (02) 34 530 214

Japan: +81 (3) 5635 8700

Korea: +82 (2) 460 7900

United Kingdom: +44 (0) 1353 658833

## Contact Us by Email

Laser Products: [tech.sales@coherent.com](mailto:tech.sales@coherent.com)

Laser Measurement: [lmc.sales@coherent.com](mailto:lmc.sales@coherent.com)


Service: [product.support@coherent.com](mailto:product.support@coherent.com)



# COHERENT®

# Diode Laser Modules and Systems

## Introduction

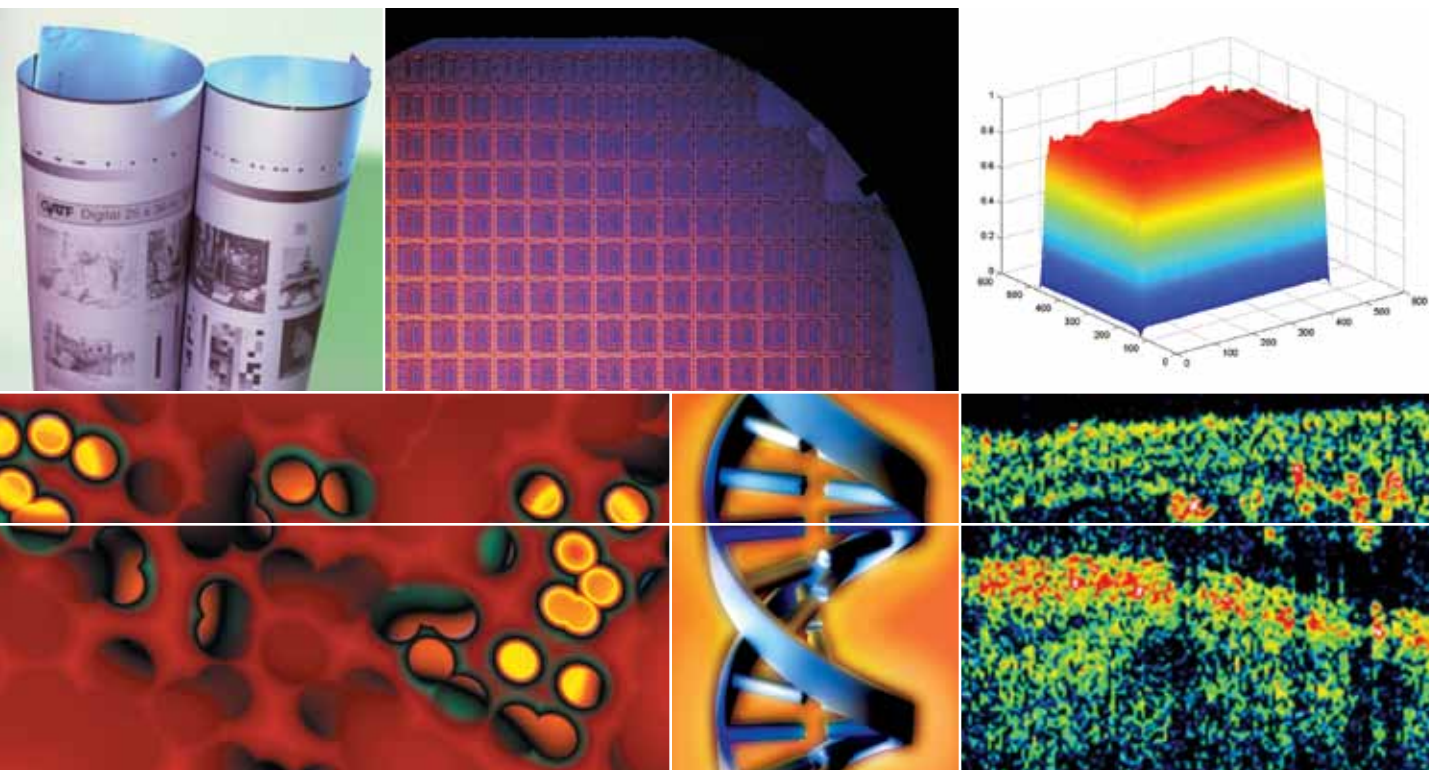

CUBE

Radius

Diode Laser  
Modules

Accessories

Coherent's Laser Diode Modules combine the technology of the diode laser with quality optics, sophisticated electronics and rugged modular packaging to provide a reliable laser source suitable to a variety of application needs.

### Applications:

Confocal Microscopy

Flow Cytometry

Medical Alignment

3D Contour Mapping

Genomics

High Speed Industrial Inspection

Homeland Security

Environmental Monitoring

Machine Vision

# Diode Laser Modules and Pattern Generators

## Table of Contents

### Introduction

|                                 |   |
|---------------------------------|---|
| The Lasiris Advantage           | 6 |
| Diode Laser Modules and Systems | 7 |

### Diode Laser Systems

|         |       |
|---------|-------|
| CUBE    | 8-19  |
| CUBE FP | 20-27 |
| Radius  | 28-30 |

### Diode Laser Modules

|                                    |       |
|------------------------------------|-------|
| ULN                                | 31-32 |
| Cross-Hair Line Generators         | 33    |
| Focusable Line Generator Modules   | 34-35 |
| Line Generator Diode Laser Modules | 36-37 |
| Miniature Diode Laser Modules      | 38-39 |

|                     |       |
|---------------------|-------|
| Application Primers | 40-43 |
|---------------------|-------|

### Lasiris Structured Light Lasers

|                            |       |
|----------------------------|-------|
| Lasiris SNF                | 44-47 |
| Lasiris Mini               | 48-51 |
| Lasiris DLS and DLSC       | 52-55 |
| Lasiris GLL                | 56-59 |
| Lasiris TMFL               | 60-63 |
| Lasiris ColdRay            | 64-67 |
| Lasiris PowerLine          | 68-71 |
| Lasiris Green PowerLine    | 72-75 |
| Lasiris Magnum II          | 76-79 |
| Lasiris Flat Top Projector | 80-83 |

# Diode Laser Modules and Pattern Generators

## Table of Contents

### Accessories

|                                           |       |
|-------------------------------------------|-------|
| Interference Filters                      | 84-85 |
| Diode LabLaser Universal Power Supply     | 86    |
| Universal Diode Laser Module Power Supply | 87    |
| Variable Angle Diode Laser Mounts         | 88    |
| FieldMaxII-TO                             | 89    |
| LaserCheck                                | 89    |
| High-Sensitivity Thermopile Sensors       | 90    |

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

# The Lasiris Advantage

## Introduction

### The Coherent Lasiris advantage

Coherent's world-renowned Lasiris™ laser's transforms the familiar laser dot into a wide range of structured light patterns including single and multiple laser lines with uniformity down to  $\pm 15\%$ . Straight laser lines are projected by allowing one dimension of light to fan out while maintaining tight control over the other, resulting in a uniform sheet-of-light.

What distinguishes the Lasiris family of laser's from conventional lasers (with cylindrical optics) is the evenness of the illumination pattern. Lasiris lasers incorporate a patented optical laser line generator that eliminates Gaussian distribution of the light, resulting in the most uniform laser lines on the market. The design provides superior quality light patterns while avoiding the intricacies of installation alignment and detector calibration.

### Increase the Performance and Reliability of your Application

#### Make full use of available power

The Lasiris non-gaussian line shows a sharp cut-off at the ends, enabling the detector to utilize the full length of the line of light instead of only the center region. All the power of a Lasiris laser is put to use.

### Simple no-nonsense calibration

Detector calibrations are easier because of the uniform intensity profile along a Lasiris line length.

### Industrial-grade protection

In addition to quality optical components, Lasiris lasers are protected from over-voltage, over-temperature, reverse-polarity, and ESD electrostatic discharge, complete with a full 2-year guarantee for selected lasers.

### Eliminate Laser Safety Problems


#### Eye-safe lasers

Since there is no hot spot and the laser energy is better spread along a Lasiris laser pattern, the safety rating of the structured light projector can be kept at the class II level for most model where other lasers cannot.

### Find the Solution that is Right for You

#### Over 48,000 off-the-shelf standard solutions

With so many standard possibilities as well as custom and OEM capabilities, whether you need a larger housing, a specific wavelength, a shorter cable, or a completely new type of laser, do not hesitate to test our skills!

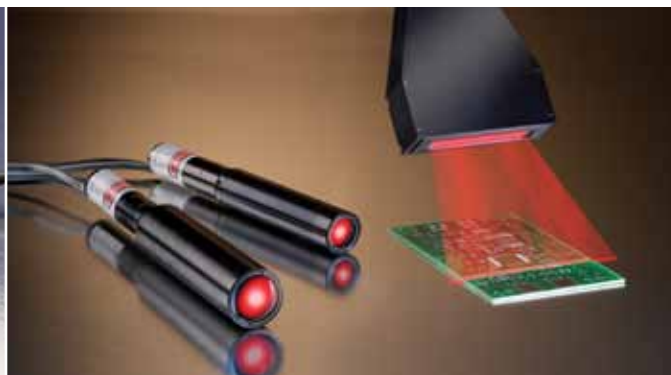

# Diode Laser Modules and Systems

## Reliable Solutions for all Your Photonics Needs

### Quality solutions

Based on size, reliability and cost, diode lasers from Coherent are the first choice for single-mode diode lasers (<100 mW) for scientific, OEM, medical or industrial applications. With the highest quality design and superior beam quality, our high-performance diode laser modules can meet your exact needs.

### Patented AlignLock technology

Only Coherent allows the laser beam to be simultaneously collimated and aligned. Other diode laser manufacturers can do one or the other, but not both. Patented AlignLock™ lets you know that the pointing stability and thermal/mechanical design is world class (US 5,111,476 and US 5,121,188).

### Rock-solid potting

Diode modules are engineered for protection from mechanical shock, optimized for thermal conductivity and vibrationally damped, so factory alignment is rock-solid and reliable.

### Electronic protection

Our diode laser modules are protected from electrostatic discharge (ESD) and are fully shielded with reverse polarity protection, meeting MIL-STD 1686B/Level 2 requirements, up to 4000 volts. With Radius™ and CUBE™ lasers, protection exceeds Level 4 ESD at 15,000 volts.

### Beam performance

Diode laser modules are smaller, require less power, generate less heat and have longer lifetimes than other laser sources. Our laser modules are configured as elliptical beams, circular beams or line generators with optional adjustable focus to increase or decrease the divergence.

### Temperature, wavelength and power

Temperature fluctuations can adversely affect diode lasers. For example, at 50°C typical red diode lifetimes are approximately 6,000 hours. However, when the diode is actively cooled, the expected lifetime becomes >30,000 hours. Both the Radius and CUBE systems offer an internal thermal-electric cooler (TEC) drive for complete thermal management and optimization of the diode laser's reliability and performance.

All Coherent diode modules feature constant power drive. Internal photodiodes provide for a closed loop power drive so that the output power is constant over temperature and lifetime. Beware of other modules with only constant-current drive, which have the disadvantage of continual power loss over the life of the laser. Coherent laser modules provide stable power—always.

We also offer the only ultra-low-noise (ULN) diode laser modules, which feature no mode-hop noise and no warm-up time, with RMS noise levels of  $\leq 0.06\%$  over temperature and lifetime of the module (bandwidth <10 MHz). ULN lasers are ideal for instrumentation and particle counting applications, providing the maximum system dynamic range so you can reliably count even the smallest particle.

### See our website

For product specification pages, product operation manuals, technical materials, competitive comparisons and prices, please visit the **Diode Modules** category at [www.Coherent.com](http://www.Coherent.com).

### Wavelength and Power Choices

Powers range from less than 1 mW up to 100 mW. Shorter wavelength blue/violet diodes provide output from 370 nm to 440 nm. For visual alignment systems, 635 nm is the wavelength of choice, while for instrumentation and measurement needs, wavelengths from 635 nm to 830 nm are available.

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

# CUBE

## Diode Laser Systems



CUBE™ lasers are high-performance, full-feature diode laser systems. Each contains a laser head, power supply, cable, and control box. The system is designed to be one of the easiest to set-up and use. The CUBE delivers the most power, stability and performance—in the smallest package—with the best value.

### A complete system

All CUBE electronics, optics and thermal management are built into the head. Every system includes a full-feature cable that lets you access all the features inside the CUBE. A compact 6V power supply, a control box for safety features, and USB software to drive the laser are also included.

### The CUBE

The CUBE offers you full control. First, the CUBE can work as a CW laser with a convenient auto-start mode. The CUBE also offers a digital input connector, located on the back panel, to drive the laser up to 150 MHz with a TTL source. Users can then command the CUBE to change the power level, either through the RS-232 or USB. For analog modulation the CUBE offers a separate analog input, located on the front panel of the control box, which allows variable power control up to 350 kHz.

For high dynamic range the digital modulation and analog power controls offer modulation depths greater than 250:1. For maximum flexibility in your application the CUBE also offers mode selection for constant-power or constant-current drive. The versatility of the CUBE includes CW, digital control, analog control and variable power operation.

With all the electronics inside the head, the control and instrumentation are available through the 15-pin connector on the laser's back panel. The interlock input, laser enable input, or auto-start mode provide control. Instrumentation on the CUBE gives you a laser-ready output, laser output power (actual laser power in an analog signal), baseplate temperature and general fault conditions. Every CUBE laser has a built-in photodiode that measures the laser final output power. The photodiode enables closed-loop power stability and allows you to know the laser power via the analog output or computer. CUBE lasers are also equipped with dual-temperature probes, which control thermal management and display both the diode and baseplate temperature. Overall, the CUBE offers complete control and instrumentation without the need for additional equipment.

# CUBE

## Diode Laser Systems

### Key Features

- Compact and modular OEM design
- Superior beam quality
- Minimal laser output RMS noise
- Complete CDRH safety features for scientific and laboratory use
- TEC temperature-stabilized and conductively cooled through baseplate
- Analog and digital modulation
- USB and RS-232

### Multiple CUBE Connections


CUBE lasers make it possible to configure multiple lasers into the same computer. Using the USB can allow up to 16 CUBE lasers to be connected on a single computer. The CUBE connection software will display each laser on its own control panel. The CUBE also has intelligence, allowing it to report which version (375 nm, 405 nm, etc.) is on the port, so you can easily drive multiple lasers simultaneously.

### The Software

The CUBE includes a full-feature software interface for Windows 2000 or Windows XP with multiple pages of status, control and configuration options. The home page includes an instrumentation dashboard for power, current, diode temperature and baseplate temperature. Buttons control laser on, mode, and power. A counter changes every one hundredth of an hour to show the actual 'on' time of the laser. The CUBE connection software makes it easy to use with help screens that feature descriptions, graphics and other assistance.


### Home Screen

The Home Screen allows for computer control of the CUBE. Control laser On/Off, set laser power and change operating modes.

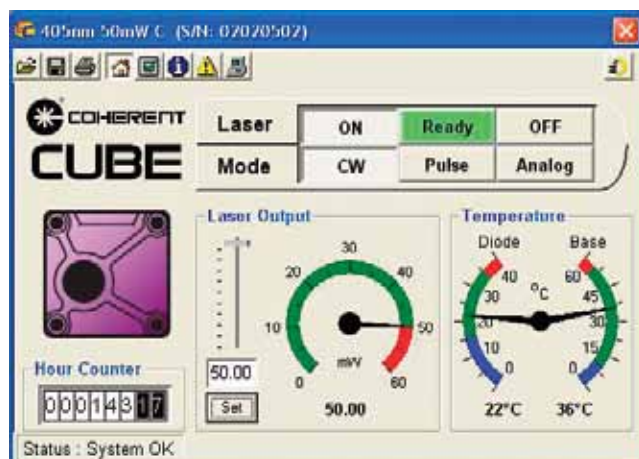

### Options Screen

The Options Screen allows for configuration choices on settings such as laser current display, 5-second delay, external analog power and laser diode temperature.

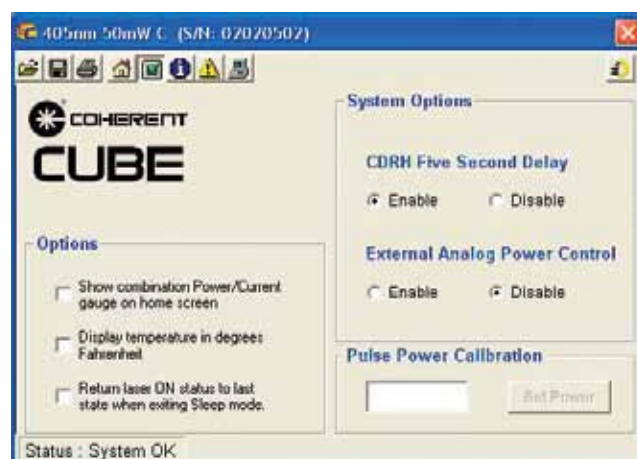

Diode Laser  
Systems


Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

# CUBE

## Diode Laser Systems

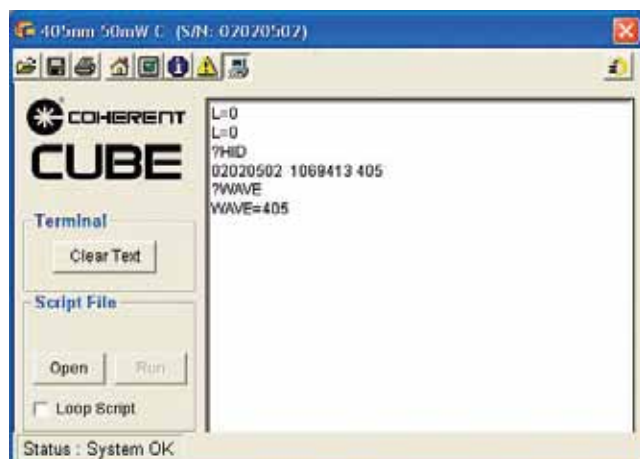


### Terminal Screen

The Terminal Screen allows for direct communication to the CUBE through the command line. Scripts can also be executed and set to loop. The Terminal Screen can accelerate the learning curve for laser control in embedded applications.

### Status Screen

The Status Screen shows the details of the laser from serial number to operating hours. A separate Fault Screen also shows any faults from an interlock fault to over-temperature.

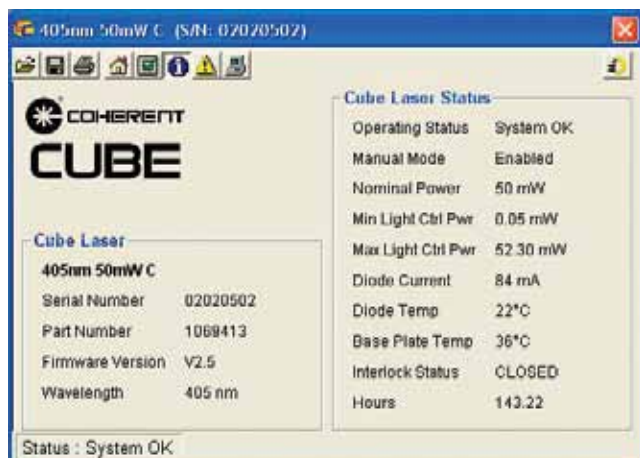

### Control and I/O

CUBE systems offer an extensive list of commands that can be executed through the RS-232 or USB. Please visit [www.coherent.com](http://www.coherent.com) and download the CUBE User Manual or the CUBE help screens for a full description of the commands and I/O.

### The Control Box

To comply with the CDRH safety requirements for a Class IIb laser, the CUBE includes both the shutter on the front of the laser and a separate control box. The control box includes the key switch that activates the laser, a five-second power-on delay, an interlock RCA jack on the back panel and a laser emission indicator. (For countries requiring a second laser emission indicator, a second lamp that connects to the cable I/O connector can be purchased separately [part no. 1079150].) The control box also offers a front panel BNC connection for the analog power control, which can use an external voltage source (0-5V) to change the CUBE's output power.

### The Cable

The CUBE's full-feature cable allows access and control to all its capabilities. The cable offers four separate lines: for power, I/O, control box, and RS-232. The I/O cable is labeled with the signal's name, pin number, and description so you can easily connect to the eight (8) separate laser control/instrumentation signals. For example, connect a voltmeter between pin numbers 2 and 7 and instantly the meter is now reporting laser power with 2V equivalent to 100% laser output power. Or connect the convenient Second Emission Lamp [part no. 1079150] as an additional safety measure.

The RS-232 cable can plug the laser directly into a computer for remote control and instrumentation.

The control box is conveniently connected to allow required safety features and analog power control from the front panel BNC. Also included with each CUBE is a USB cable, which allows the laser to plug directly into any Windows\* 2000 or Windows XP system.

### The Power Supply

The CUBE, which has been designed for easy integration into your system with a single power supply input, can accept from 4.8V to 6.5VDC.

# CUBE

## Diode Laser Systems

Every CUBE system comes standard with a compact 6VDC supply with a power switch and an LED. The power supply is equipped with an IEC320 universal input, which allows most computer-style power cords to be easily used with your country's power (U.S. -style power cord is included).

### Reliability

The CUBE is fully protected from electro-static discharge (ESD) to Level 4 standards and is also equipped with temperate shutdown systems to protect the laser diode. The diode and optics are temperature-controlled through a TEC to provide the optimal performance and maximum lifetime. Reliability is the priority—within every CUBE is a computer processor controlling all the aspects of the laser drive, thermal systems and protection measures.

### The Extras

For even simpler set-up and integration, a heat sink accessory provides a convenient 2.5" beam height and a universal bolt pattern for easy mounting—with the beam axis traveling directly down a row of holes on the laser table. For high-temperature applications a fan inside the heat sink provides additional cooling. Please note that the CUBE comes standard with the laser head, control box, cables, power supply, manual, software, and a U.S. power cord. The heat sink, however, is sold separately.

Overall, CUBE lasers can give you more powerful and efficient ways to use diode laser light in the fields of bioinstrumentation, particle counting, confocal microscopy, high-throughput screening, or general laboratory work.



Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

# CUBE

## Diode Laser Systems

| System Specifications | CUBE                                                  | 375-16C                        | 405-50C                                 | 405-100C |
|-----------------------|-------------------------------------------------------|--------------------------------|-----------------------------------------|----------|
|                       | Part Number                                           | 1112774                        | 1174298                                 | 1142279  |
|                       | Wavelength <sup>1</sup> (nm)                          | 375                            | 405                                     | 405      |
|                       | Output Power <sup>2</sup> (mW)                        | 16                             | 50                                      | 100      |
|                       | Beam Diameter at 1/e <sup>2</sup> (mm)(typical)       | 1.1                            | 1.4                                     | 1.4      |
|                       | Beam Divergence (mrad)(typical)                       | 0.5                            | 0.4                                     | 0.4      |
|                       | M <sup>2</sup> (ModeMaster beam quality) <sup>3</sup> | ≤1.3                           | ≤1.3                                    | ≤1.2     |
|                       | Beam Asymmetry <sup>3</sup>                           |                                | <1:1.2                                  |          |
|                       | Laser Drive Modes                                     | Digital, Analog and PC Control |                                         |          |
|                       | Digital Modulation                                    |                                |                                         |          |
|                       | Maximum Bandwidth (MHz)                               |                                | 150                                     |          |
|                       | Rise Time (10% to 90%) (nsec)                         |                                | <2                                      |          |
|                       | Fall Time (10% to 90%) (nsec)                         |                                | <2                                      |          |
|                       | Modulation Depth (extinction ratio)                   |                                | >1,000,000:1 at 0 Hz, >250:1 at 150 MHz |          |
|                       | Analog Modulation                                     |                                |                                         |          |
|                       | Maximum Bandwidth (kHz)                               | 70                             | 350                                     | 350      |
|                       | Rise or Fall Time (10% to 90%) (μsec)                 | <5                             | <1                                      | <1       |
|                       | Modulation Depth (extinction ratio)                   |                                | >10,000:1                               |          |
|                       | Laser Enable Control Input                            |                                |                                         |          |
|                       | Maximum Bandwidth (kHz)                               |                                | 130                                     |          |
|                       | Rise or Fall Time (10% to 90%) (μsec)                 |                                | <1                                      |          |
|                       | Modulation Depth (extinction ratio)                   |                                | ∞, complete On/Off                      |          |
|                       | RMS Noise                                             |                                |                                         |          |
|                       | 20 Hz to 10 MHz                                       | ≤0.15%                         |                                         | ≤0.1%    |
|                       | 10 MHz to 500 MHz                                     |                                | <1.0%                                   |          |
|                       | Long-Term Power Stability (8 hours)                   |                                | <±2%                                    |          |
|                       | Warm-Up Time (minutes)                                |                                | <5                                      |          |
|                       | Spatial Mode (far field)                              |                                | TEM <sub>00</sub>                       |          |
|                       | Pointing Stability (μrad/°C)                          |                                | <6                                      |          |
|                       | Polarization Ratio <sup>4</sup>                       |                                | Min. 100:1, Vert. ±5°                   |          |
|                       | Static Alignment Tolerances <sup>5</sup>              |                                |                                         |          |
|                       | Beam Position (mm)                                    |                                | <1                                      |          |
|                       | Beam Angle (mrad)                                     |                                | <5                                      |          |
|                       | CDRH Laser Classification                             |                                | Class IIIb                              |          |
|                       | ESD Protection <sup>6</sup>                           |                                | Level 4                                 |          |

<sup>1</sup> Laser-to-laser tolerance. All lasers ±5 nm except CUBE 640-40 with 635 to 642 nm range, CUBE 640-100 with 635 to 644 nm range, CUBE 660 with 652 to 665 nm range, and CUBE 685, 730 and 785 with ±10 nm.

<sup>2</sup> Output power is variable via RS-232, USB or analog interface from 1% to 100%. Specifications are valid for 100% power. Recommended power range is 70% to 100% power.

<sup>3</sup> M<sup>2</sup> and Beam Asymmetry measured with ModeMaster with 90/10 clip levels.

<sup>4</sup> Polarization measurement is max./min. power with polarizer at 32 in. from laser with LM-2 power head at 37 in. from laser.

<sup>5</sup> Static alignment tolerances are relative to the right bottom edge (in beam direction).

<sup>6</sup> Electro-Static Discharge Standard IEC 1000-4-2, 1995.

# CUBE

## Diode Laser Systems

| System Specifications                                 | CUBE | 445-40C                                 | 488-50C               |
|-------------------------------------------------------|------|-----------------------------------------|-----------------------|
| Part Number                                           |      | 1069418                                 | 1174040               |
| Wavelength <sup>1</sup> (nm)                          |      | 445                                     | 488                   |
| Output Power <sup>2</sup> (mW)                        |      | 40                                      | 50                    |
| Beam Diameter at 1/e <sup>2</sup> (mm)(typical)       |      | 1.1                                     | 1.3                   |
| Beam Divergence (mrad)(typical)                       |      |                                         | 0.6                   |
| M <sup>2</sup> (ModeMaster beam quality) <sup>3</sup> |      |                                         | ≤1.2                  |
| Beam Asymmetry <sup>3</sup>                           |      |                                         | <1:1.2                |
| Laser Drive Modes                                     |      | Digital, Analog and PC Control          |                       |
| Digital Modulation                                    |      |                                         |                       |
| Maximum Bandwidth (MHz)                               |      |                                         | 125                   |
| Rise Time (10% to 90%) (nsec)                         |      |                                         | <2                    |
| Fall Time (10% to 90%) (nsec)                         |      |                                         | <3                    |
| Modulation Depth (extinction ratio)                   |      | >1,000,000:1 at 0 Hz, >250:1 at 150 MHz |                       |
| Analog Modulation                                     |      |                                         |                       |
| Maximum Bandwidth (kHz)                               |      |                                         | 350                   |
| Rise or Fall Time (10% to 90%) (μsec)                 |      |                                         | <1                    |
| Modulation Depth (extinction ratio)                   |      |                                         | >10,000:1             |
| Laser Enable Control Input                            |      |                                         |                       |
| Maximum Bandwidth (kHz)                               |      |                                         | 130                   |
| Rise or Fall Time (10% to 90%) (μsec)                 |      |                                         | <1                    |
| Modulation Depth (extinction ratio)                   |      |                                         | ∞, complete On/Off    |
| RMS Noise                                             |      |                                         |                       |
| 20 Hz to 10 MHz                                       |      |                                         | ≤0.1%                 |
| 10 MHz to 500 MHz                                     |      |                                         | <1.0%                 |
| Long-Term Power Stability (8 hours)                   |      |                                         | <±2%                  |
| Warm-Up Time (minutes)                                |      |                                         | <5                    |
| Spatial Mode (far field)                              |      |                                         | TEM <sub>00</sub>     |
| Pointing Stability (μrad/°C)                          |      |                                         | <6                    |
| Polarization Ratio <sup>4</sup>                       |      |                                         | Min. 100:1, Vert. ±5° |
| Static Alignment Tolerances <sup>5</sup>              |      |                                         |                       |
| Beam Position (mm)                                    |      |                                         | <1                    |
| Beam Angle (mrad)                                     |      |                                         | <5                    |
| CDRH Laser Classification                             |      |                                         | Class IIIb            |
| ESD Protection <sup>6</sup>                           |      |                                         | Level 4               |

<sup>1</sup> Laser-to-laser tolerance. All lasers ±5 nm except CUBE 640-40 with 635 to 642 nm range, CUBE 640-100 with 635 to 644 nm range, CUBE 660 with 652 to 665 nm range, and CUBE 685, 730 and 785 with ±10 nm.

<sup>2</sup> Output power is variable via RS-232, USB or analog interface from 1% to 100%. Specifications are valid for 100% power. Recommended power range is 70% to 100% power.

<sup>3</sup> M<sup>2</sup> and Beam Asymmetry measured with ModeMaster with 90/10 clip levels.

<sup>4</sup> Polarization measurement is max./min. power with polarizer at 32 in. from laser with LM-2 power head at 37 in. from laser.

<sup>5</sup> Static alignment tolerances are relative to the right bottom edge (in beam direction).

<sup>6</sup> Electro-Static Discharge Standard IEC 1000-4-2, 1995.

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

# CUBE

## Diode Laser Systems

| System Specifications                                 | CUBE | 635-30E                                 | 640-40C               | 640-100C |
|-------------------------------------------------------|------|-----------------------------------------|-----------------------|----------|
| Part Number                                           |      | 1069410                                 | 1069417               | 1150205  |
| Wavelength <sup>1</sup> (nm)                          |      | 635                                     | 640                   | 640      |
| Output Power <sup>2</sup> (mW)                        |      | 30                                      | 40                    | 100      |
| Beam Diameter at 1/e <sup>2</sup> (mm)(typical)       |      | 1 x 3.8                                 | 1                     | 1        |
| Beam Divergence (mrad)(typical)                       |      | 0.8 x 0.3                               | 1                     | 1        |
| M <sup>2</sup> (ModeMaster beam quality) <sup>3</sup> |      |                                         | ≤1.2                  |          |
| Beam Asymmetry <sup>3</sup>                           |      | 1:3.8                                   | <1:1.2                | <1:1.2   |
| Laser Drive Modes                                     |      | Digital, Analog and PC Control          |                       |          |
| Digital Modulation                                    |      |                                         |                       |          |
| Maximum Bandwidth (MHz)                               |      |                                         | 150                   |          |
| Rise Time (10% to 90%) (nsec)                         |      |                                         | <2                    |          |
| Fall Time (10% to 90%) (nsec)                         |      |                                         | <2                    |          |
| Modulation Depth (extinction ratio)                   |      | >1,000,000:1 at 0 Hz, >250:1 at 150 MHz |                       |          |
| Analog Modulation                                     |      |                                         |                       |          |
| Maximum Bandwidth (kHz)                               |      |                                         | 350                   |          |
| Rise or Fall Time (10% to 90%) (μsec)                 |      |                                         | <1                    |          |
| Modulation Depth (extinction ratio)                   |      |                                         | >10,000:1             |          |
| Laser Enable Control Input                            |      |                                         |                       |          |
| Maximum Bandwidth (kHz)                               |      |                                         | 130                   |          |
| Rise or Fall Time (10% to 90%) (μsec)                 |      |                                         | <1                    |          |
| Modulation Depth (extinction ratio)                   |      |                                         | ∞, complete On/Off    |          |
| RMS Noise                                             |      |                                         |                       |          |
| 20 Hz to 10 MHz                                       |      |                                         | ≤0.1%                 |          |
| 10 MHz to 500 MHz                                     |      |                                         | <1.0%                 |          |
| Long-Term Power Stability (8 hours)                   |      |                                         | <±2%                  |          |
| Warm-Up Time (minutes)                                |      |                                         | <5                    |          |
| Spatial Mode (far field)                              |      |                                         | TEM <sub>00</sub>     |          |
| Pointing Stability (μrad/°C)                          |      |                                         | <6                    |          |
| Polarization Ratio <sup>4</sup>                       |      |                                         | Min. 100:1, Vert. ±5° |          |
| Static Alignment Tolerances <sup>5</sup>              |      |                                         |                       |          |
| Beam Position (mm)                                    |      |                                         | <1                    |          |
| Beam Angle (mrad)                                     |      |                                         | <5                    |          |
| CDRH Laser Classification                             |      |                                         | Class IIIb            |          |
| ESD Protection <sup>6</sup>                           |      |                                         | Level 4               |          |

<sup>1</sup> Laser-to-laser tolerance. All lasers ±5 nm except CUBE 640-40 with 635 to 642 nm range, CUBE 640-100 with 635 to 644 nm range, CUBE 660 with 652 to 665 nm range, and CUBE 685, 730 and 785 with ±10 nm.

<sup>2</sup> Output power is variable via RS-232, USB or analog interface from 1% to 100%. Specifications are valid for 100% power. Recommended power range is 70% to 100% power.

<sup>3</sup> M<sup>2</sup> and Beam Asymmetry measured with ModeMaster with 90/10 clip levels.

<sup>4</sup> Polarization measurement is max./min. power with polarizer at 32 in. from laser with LM-2 power head at 37 in. from laser.

<sup>5</sup> Static alignment tolerances are relative to the right bottom edge (in beam direction).

<sup>6</sup> Electro-Static Discharge Standard IEC 1000-4-2, 1995.

# CUBE

## Diode Laser Systems

| System Specifications | CUBE                                                  | 660-60C                                 | 660-100C              | 685-40C |
|-----------------------|-------------------------------------------------------|-----------------------------------------|-----------------------|---------|
|                       | Part Number                                           | 1094060                                 | 1130061               | 1184303 |
|                       | Wavelength <sup>1</sup> (nm)                          | 660                                     | 660                   | 685     |
|                       | Output Power <sup>2</sup> (mW)                        | 60                                      | 100                   | 40      |
|                       | Beam Diameter at 1/e <sup>2</sup> (mm)(typical)       | 1.4                                     | 1                     | 1       |
|                       | Beam Divergence (mrad)(typical)                       | 0.6                                     | 1                     | 1       |
|                       | M <sup>2</sup> (ModeMaster beam quality) <sup>3</sup> |                                         | ≤1.2                  |         |
|                       | Beam Asymmetry <sup>3</sup>                           |                                         | <1:1.2                |         |
|                       | Laser Drive Modes                                     | Digital, Analog and PC Control          |                       |         |
|                       | Digital Modulation                                    |                                         |                       |         |
|                       | Maximum Bandwidth (MHz)                               | 150                                     | 150                   | 150     |
|                       | Rise Time (10% to 90%) (nsec)                         | <2                                      | <2                    | <2      |
|                       | Fall Time (10% to 90%) (nsec)                         | <2                                      | <2                    | <2      |
|                       | Modulation Depth (extinction ratio)                   | >1,000,000:1 at 0 Hz, >250:1 at 150 MHz |                       |         |
|                       | Analog Modulation                                     |                                         |                       |         |
|                       | Maximum Bandwidth (kHz)                               |                                         | 350                   |         |
|                       | Rise or Fall Time (10% to 90%) (μsec)                 |                                         | <1                    |         |
|                       | Modulation Depth (extinction ratio)                   |                                         | >10,000:1             |         |
|                       | Laser Enable Control Input                            |                                         |                       |         |
|                       | Maximum Bandwidth (kHz)                               |                                         | 130                   |         |
|                       | Rise or Fall Time (10% to 90%) (μsec)                 |                                         | <1                    |         |
|                       | Modulation Depth (extinction ratio)                   |                                         | ∞, complete On/Off    |         |
|                       | RMS Noise                                             |                                         |                       |         |
|                       | 20 Hz to 10 MHz                                       |                                         | ≤0.1%                 |         |
|                       | 10 MHz to 500 MHz                                     |                                         | <1.0%                 |         |
|                       | Long-Term Power Stability (8 hours)                   |                                         | <±2%                  |         |
|                       | Warm-Up Time (minutes)                                |                                         | <5                    |         |
|                       | Spatial Mode (far field)                              |                                         | TEM <sub>00</sub>     |         |
|                       | Pointing Stability (μrad/°C)                          |                                         | <6                    |         |
|                       | Polarization Ratio <sup>4</sup>                       |                                         | Min. 100:1, Vert. ±5° |         |
|                       | Static Alignment Tolerances <sup>5</sup>              |                                         |                       |         |
|                       | Beam Position (mm)                                    |                                         | <1                    |         |
|                       | Beam Angle (mrad)                                     |                                         | <5                    |         |
|                       | CDRH Laser Classification                             |                                         | Class IIIb            |         |
|                       | ESD Protection <sup>6</sup>                           |                                         | Level 4               |         |

<sup>1</sup> Laser-to-laser tolerance. All lasers ±5 nm except CUBE 640-40 with 635 to 642 nm range, CUBE 640-100 with 635 to 644 nm range, CUBE 660 with 652 to 665 nm range, and CUBE 685, 730 and 785 with ±10 nm.

<sup>2</sup> Output power is variable via RS-232, USB or analog interface from 1% to 100%. Specifications are valid for 100% power. Recommended power range is 70% to 100% power.

<sup>3</sup> M<sup>2</sup> and Beam Asymmetry measured with ModeMaster with 90/10 clip levels.

<sup>4</sup> Polarization measurement is max./min. power with polarizer at 32 in. from laser with LM-2 power head at 37 in. from laser.

<sup>5</sup> Static alignment tolerances are relative to the right bottom edge (in beam direction).

<sup>6</sup> Electro-Static Discharge Standard IEC 1000-4-2, 1995.

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories

# CUBE

## Diode Laser Systems

| System Specifications                                 | CUBE | 730-30C                                 | 785-40C |
|-------------------------------------------------------|------|-----------------------------------------|---------|
| Part Number                                           |      | 1178855                                 | 1069416 |
| Wavelength <sup>1</sup> (nm)                          |      | 730                                     | 785     |
| Output Power <sup>2</sup> (mW)                        |      | 30                                      | 40      |
| Beam Diameter at 1/e <sup>2</sup> (mm)(typical)       |      | 1.2                                     | 1.2     |
| Beam Divergence (mrad)(typical)                       |      | 0.9                                     | 1       |
| M <sup>2</sup> (ModeMaster beam quality) <sup>3</sup> |      | ≤1.2                                    |         |
| Beam Asymmetry <sup>3</sup>                           |      | <1:1.2                                  |         |
| Laser Drive Modes                                     |      | Digital, Analog and PC Control          |         |
| Digital Modulation                                    |      |                                         |         |
| Maximum Bandwidth (MHz)                               |      | 150                                     |         |
| Rise Time (10% to 90%) (nsec)                         |      | <2                                      |         |
| Fall Time (10% to 90%) (nsec)                         |      | <2                                      |         |
| Modulation Depth (extinction ratio)                   |      | >1,000,000:1 at 0 Hz, >250:1 at 150 MHz |         |
| Analog Modulation                                     |      |                                         |         |
| Maximum Bandwidth (kHz)                               |      | 350                                     |         |
| Rise or Fall Time (10% to 90%) (μsec)                 |      | <1                                      |         |
| Modulation Depth (extinction ratio)                   |      | >10,000:1                               |         |
| Laser Enable Control Input                            |      |                                         |         |
| Maximum Bandwidth (kHz)                               |      | 130                                     |         |
| Rise or Fall Time (10% to 90%) (μsec)                 |      | <1                                      |         |
| Modulation Depth (extinction ratio)                   |      | ∞, complete On/Off                      |         |
| RMS Noise                                             |      |                                         |         |
| 20 Hz to 10 MHz                                       |      | ≤0.1%                                   |         |
| 10 MHz to 500 MHz                                     |      | <1.0%                                   |         |
| Long-Term Power Stability (8 hours)                   |      | <±2%                                    |         |
| Warm-Up Time (minutes)                                |      | <5                                      |         |
| Spatial Mode (far field)                              |      | TEM <sub>00</sub>                       |         |
| Pointing Stability (μrad/°C)                          |      | <6                                      |         |
| Polarization Ratio <sup>4</sup>                       |      | Min. 100:1, Vert. ±5°                   |         |
| Static Alignment Tolerances <sup>5</sup>              |      |                                         |         |
| Beam Position (mm)                                    |      | <1                                      |         |
| Beam Angle (mrad)                                     |      | <5                                      |         |
| CDRH Laser Classification                             |      | Class IIIb                              |         |
| ESD Protection <sup>6</sup>                           |      | Level 4                                 |         |

<sup>1</sup> Laser-to-laser tolerance. All lasers ±5 nm except CUBE 640-40 with 635 to 642 nm range, CUBE 640-100 with 635 to 644 nm range, CUBE 660 with 652 to 665 nm range, and CUBE 685, 730 and 785 with ±10 nm.

<sup>2</sup> Output power is variable via RS-232, USB or analog interface from 1% to 100%. Specifications are valid for 100% power. Recommended power range is 70% to 100% power.

<sup>3</sup> M<sup>2</sup> and Beam Asymmetry measured with ModeMaster with 90/10 clip levels.

<sup>4</sup> Polarization measurement is max./min. power with polarizer at 32 in. from laser with LM-2 power head at 37 in. from laser.

<sup>5</sup> Static alignment tolerances are relative to the right bottom edge (in beam direction).

<sup>6</sup> Electro-Static Discharge Standard IEC 1000-4-2, 1995.

# CUBE

## Diode Laser Systems

### Utility and Environmental Requirements

|                                                      |                             |
|------------------------------------------------------|-----------------------------|
| Laser Head Operating Voltage                         | +4.8 to 6.5 VDC             |
| Laser Head Operating Current                         | <2.5 amps                   |
| Baseplate Temperature Range <sup>1</sup>             | +10° to 50°C (50° to 122°F) |
| Maximum Heat Dissipation of Head (baseplate at 50°C) | 13W                         |
| Storage Temperature                                  | -20° to 60°C (-4° to 140°F) |
| Power Supply <sup>2</sup> (included)                 |                             |
| Operating Voltage <sup>3</sup>                       | 100 to 240 VAC, 50 to 60 Hz |
| Power Consumption                                    | <15W                        |


|                      | Dimensions (L x W x H)                 | Weight             |
|----------------------|----------------------------------------|--------------------|
| Laser Head           | 100 x 40 x 40 mm (3.9 x 1.6 x 1.6 in.) | 0.3 kg (0.7 lbs.)  |
| Control Box          | 138 x 62 x 26 mm (5.4 x 2.4 x 1.1 in.) | 0.12 kg (0.3 lbs.) |
| Power Supply         | 87 x 47 x 32 mm (3.4 x 1.9 x 1.3 in.)  | 0.22 kg (0.5 lbs.) |
| Head Cable           | 254 to 406 mm (10 to 18 in.)           | 0.08 kg (0.2 lbs.) |
| USB Cable            | 1.8m (6 ft.)                           | 0.1 kg (0.2 lbs.)  |
| Heat Sink (optional) | 127 x 66 x 44 mm (5 x 2.6 x 1.74 in.)  | 0.4 kg (0.9 lbs.)  |

<sup>1</sup> Non-condensing with laser diode TEC at set point of 22°C.

<sup>2</sup> On-off power switch and LED included in cord to laser head.

<sup>3</sup> System includes USA-type power cord, IEC320 input connection.

### CUBE Laser Head

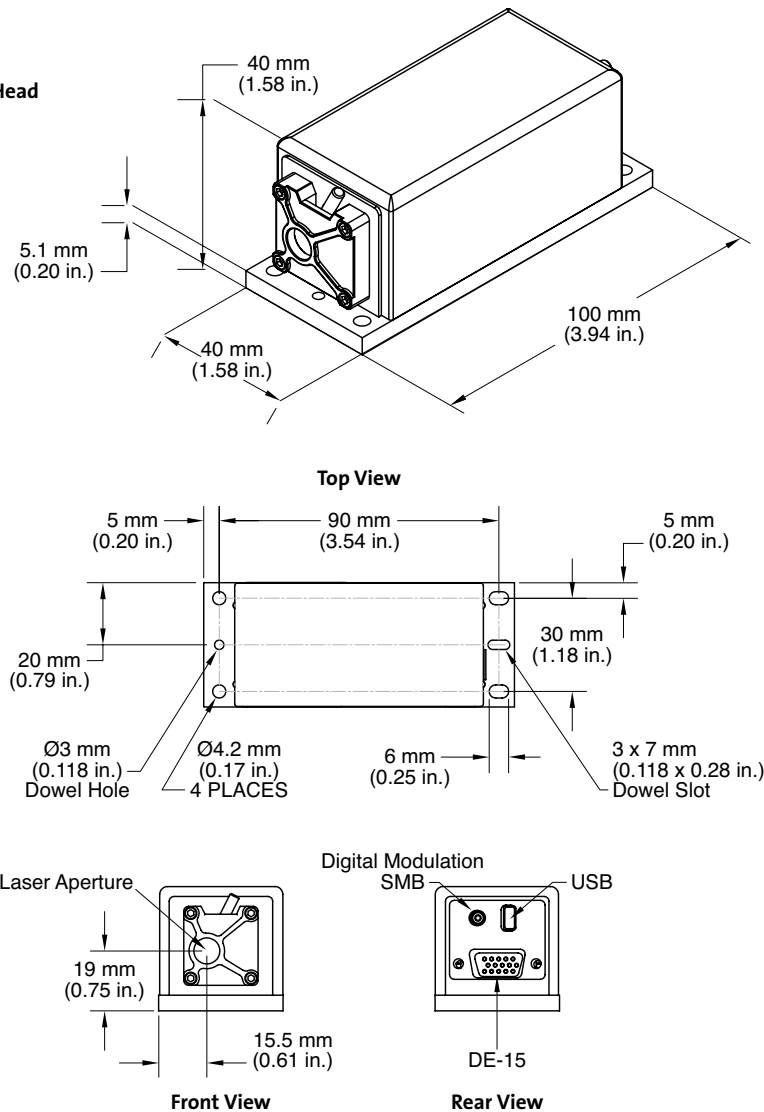

# CUBE


## Diode Laser Systems

|                         |           |                                           |                                                                                           |
|-------------------------|-----------|-------------------------------------------|-------------------------------------------------------------------------------------------|
| <b>CUBE Accessories</b> | 1073840   | CUBE Heat Sink Accessory                  | CUBE heat sink with fan, 127 x 66 x 44 mm                                                 |
|                         | 1079150   | CUBE 2nd Emission Lamp                    | CDRH lamp for laser emission. Connects to CUBE I/O cable.                                 |
|                         | 1072166   | *CUBE Interface Cable                     | CUBE to control box, I/O, RS-232 and power                                                |
|                         | 1072454   | *CUBE Power Supply, 6 VDC, 2.5 amps       | IEC320 input, 110V/220V, switch and LED                                                   |
|                         | 1108906   | *CUBE USB Cable, 6 ft. (1.8m)             | USB standard to mini for CUBE                                                             |
|                         | 1039966   | * Control Box, CDRH, Keyswitch, Interlock | Control box for CUBE lasers. Included box, key, interlock RCA and cable.                  |
|                         | 1080090   | Cable, Extension, DB-9 M/F, 6 ft. (1.8m)  | Extension cable to extend control box or RS-232. M-F, DB-9.                               |
|                         | 1040408   | * Plug, RCA                               | Interlock RCA plug for control box, non-shorted                                           |
|                         | 1079890   | *CUBE Manual and Software                 | CUBE user manual and software CD for Windows 2000** and Windows XP**                      |
|                         | 6005-0059 | *Power Cord, USA, 6 ft. (1.8m)            | Power cord, USA plug style to IEC320                                                      |
|                         | 1116779   | CUBE Heat Sink Right Angle                | CUBE right angle mount, 102 x 47 x 44 mm, same beam position with horizontal polarization |


\* Included free with every CUBE system.

\*\* Windows 2000 and Windows XP are registered trademarks of Microsoft Corporation.

**CUBE Heat Sink**

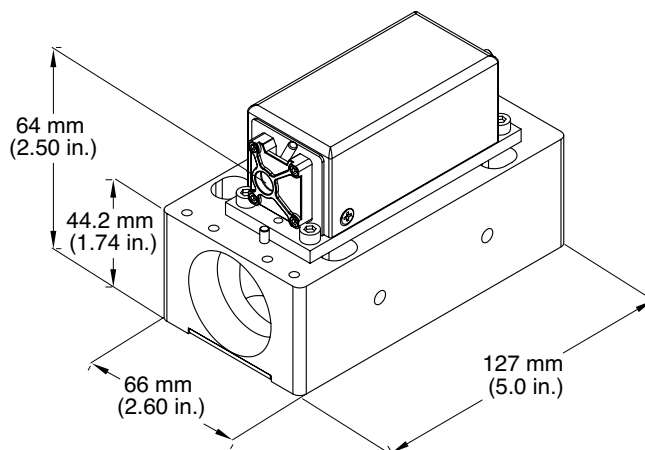

**Heat Sink Top View**

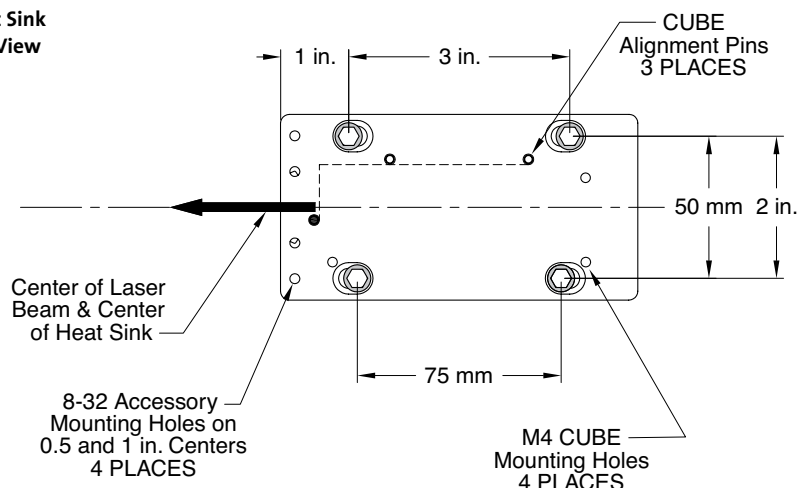

# CUBE

## Diode Laser Systems




Diode Laser  
Systems

Diode Laser  
Modules


Lasiris  
Structured  
Light Lasers

Accessories

**CUBE  
Power Supply**

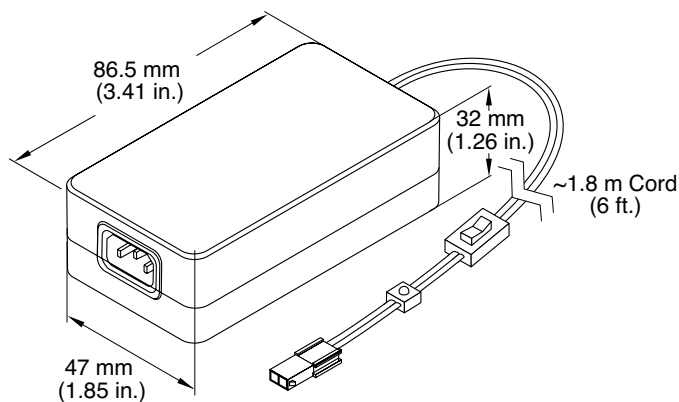

**CUBE  
Control Box**

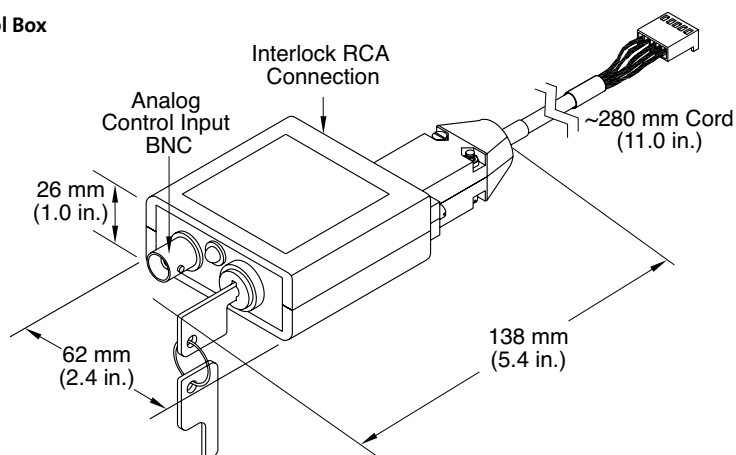

# CUBE FP


## Diode Laser System



### Key Features

- Permanent fiber attachment – constant power without adjustments
- Compact and modular OEM design
- Complete CDRH safety features for scientific and laboratory use
- TEC temperature-stabilized and conductively cooled through baseplate
- Analog and digital modulation
- USB and RS-232

### CUBE FP Laser Head

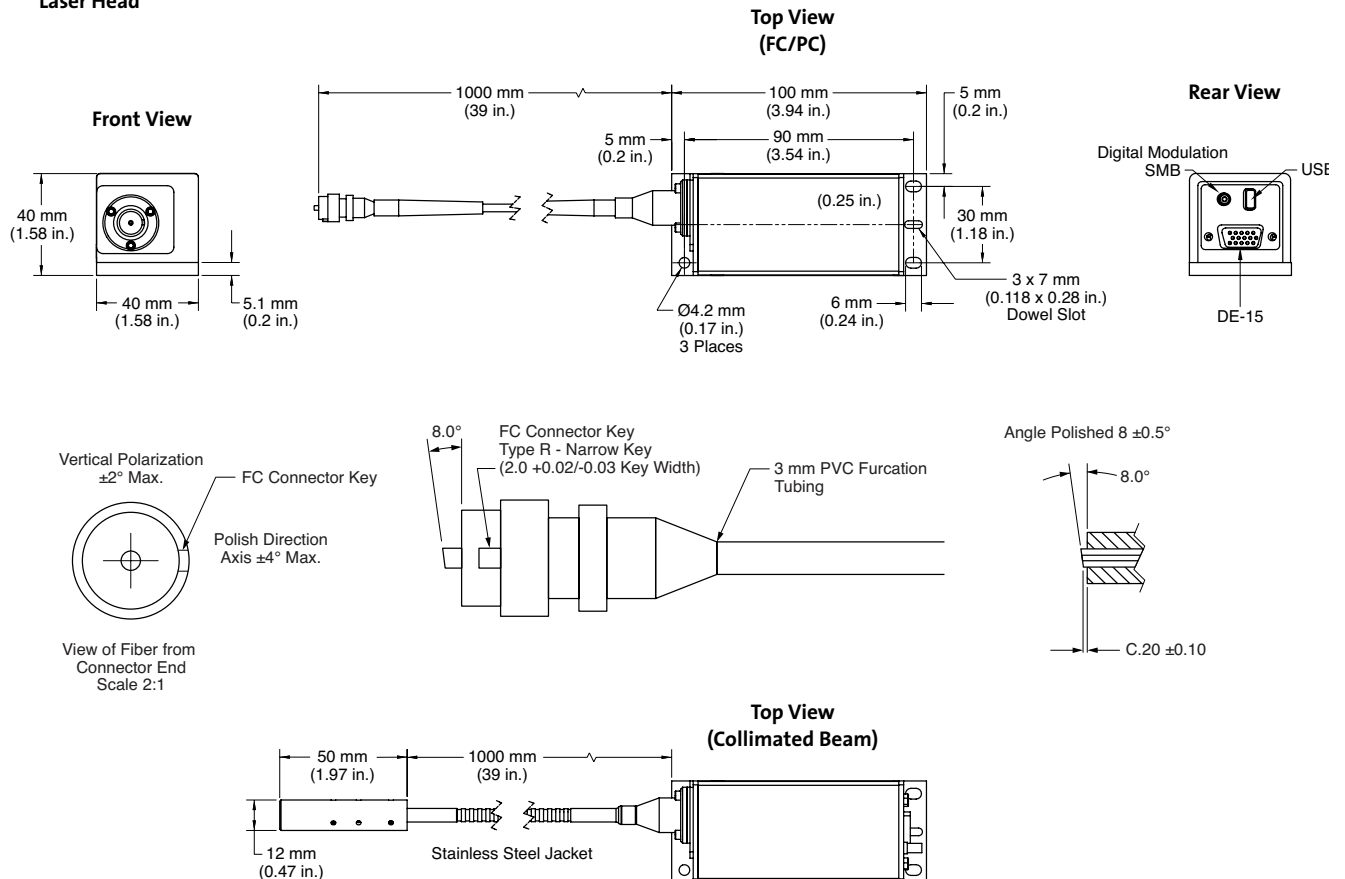

# CUBE FP

## Diode Laser System

| System Specifications                            | CUBE FP | 405-50FP                        | 405-50FP                                |
|--------------------------------------------------|---------|---------------------------------|-----------------------------------------|
| Part Number                                      |         | 1139602                         | 1139601                                 |
| Output from Fiber                                |         | FC/APC (8° angled) <sup>6</sup> | Collimated Beam                         |
| Fiber Cable Type                                 |         | 3 mm PVC                        | Stainless Steel                         |
| Fiber Cable Length (m)(minimum)                  |         | 1                               | 1                                       |
| Wavelength <sup>1</sup> (nm)                     |         | 405                             | 405                                     |
| Output Power <sup>2</sup> (mW)                   |         | 50                              | 50                                      |
| Beam Diameter (1/e <sup>2</sup> )(mm)(typical)   |         | -                               | 1                                       |
| Beam Divergence (mrad)(typical)                  |         | -                               | 1                                       |
| Fiber Numerical Aperture (NA)(1/e <sup>2</sup> ) |         | 0.05                            | -                                       |
| Mode-field Diameter (μm)(typical)                |         |                                 | 3.5                                     |
| M <sup>2</sup> (ModeMaster beam quality)         |         |                                 | ≤1.1                                    |
| Beam Asymmetry                                   |         |                                 | ≤1:1.1                                  |
| Laser Drive Modes                                |         | Digital, Analog and PC Control  |                                         |
| Digital Modulation                               |         |                                 |                                         |
| Maximum Bandwidth (MHz)                          |         |                                 | 150                                     |
| Rise Time (10% to 90%)(nsec)                     |         |                                 | <2                                      |
| Fall Time (10% to 90%)(nsec)                     |         |                                 | <2                                      |
| Modulation Depth (extinction ratio)              |         |                                 | >1,000,000:1 at 0 Hz, >250:1 at 150 MHz |
| Analog Modulation                                |         |                                 |                                         |
| Maximum Bandwidth (kHz)                          |         |                                 | 350                                     |
| Rise or Fall Time (10% to 90%)(μsec)             |         |                                 | <1                                      |
| Modulation Depth (extinction ratio)              |         |                                 | >10,000:1                               |
| Laser Enable Control Input                       |         |                                 |                                         |
| Maximum Bandwidth (kHz)                          |         |                                 | 130                                     |
| Rise or Fall Time (10% to 90%)(μsec)             |         |                                 | <1                                      |
| Modulation Depth (extinction ratio)              |         |                                 | ∞, complete On/Off                      |
| RMS Noise                                        |         |                                 |                                         |
| 20 Hz to 10 MHz                                  |         |                                 | <0.2%                                   |
| 10 MHz to 500 MHz                                |         |                                 | <1.0%                                   |
| Power Stability (8 hours, ±3°C)                  |         |                                 | ≤±2%                                    |
| Long-Term Output Power Average                   |         |                                 | ≤5%/1000 hrs.                           |
| Warm-Up Time (minutes)                           |         |                                 | <5                                      |
| Spatial Mode (far field)                         |         |                                 | TEM <sub>00</sub>                       |
| Pointing Stability (μrad/°C)                     |         | -                               | <6                                      |
| Polarization Ratio <sup>3</sup>                  |         |                                 | Minimum 100:1                           |
| Tolerances <sup>4</sup>                          |         |                                 |                                         |
| Beam Position (mm)                               |         | -                               | ≤0.5                                    |
| Beam Angle (mrad)                                |         | -                               | ≤15                                     |
| CDRH Laser Classification                        |         |                                 | Class IIIB                              |
| ESD Protection <sup>5</sup>                      |         |                                 | Level 4                                 |

<sup>1</sup> Laser-to-laser tolerance. CUBE 640-30 has center wavelength range of 635 nm to 642 nm, CUBE 640-75 with 635 nm to 644 nm range, CUBE 660 with 652 nm to 665 nm range, CUBE 405, 445 and 488 with ±5 nm range.

<sup>2</sup> Output power is variable via RS-232, USB or analog interface from 1% to 100%. Specifications are valid for 100% to 110% power. Recommended power range is 70% to 100% power. Factory setting 100% to 110% of power specification.

<sup>3</sup> Polarization measurement is max./min. power with polarizer at 32 cm from laser with LM-2 power head at 37 cm from laser.

<sup>4</sup> Static alignment tolerances are relative to center of output beam aperture.

<sup>5</sup> Electro-Static Discharge Standard IEC 1000-4-2, 1995.

<sup>6</sup> Fiber FC/APC connector output not compatible for patchcord-to-patchcord connection.

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

# CUBE FP

## Diode Laser System

| System Specifications                            | CUBE FP | 445-25FP                        | 445-25FP                                |
|--------------------------------------------------|---------|---------------------------------|-----------------------------------------|
| Part Number                                      |         | 1169225                         | 1169224                                 |
| Output from Fiber                                |         | FC/APC (8° angled) <sup>6</sup> | Collimated Beam                         |
| Fiber Cable Type                                 |         | 3 mm PVC                        | Stainless Steel                         |
| Fiber Cable Length (m)(minimum)                  |         | 1                               | 1                                       |
| Wavelength <sup>1</sup> (nm)                     |         | 445                             | 445                                     |
| Output Power <sup>2</sup> (mW)                   |         | 25                              | 25                                      |
| Beam Diameter (1/e <sup>2</sup> )(mm)(typical)   |         | -                               | 1                                       |
| Beam Divergence (mrad)(typical)                  |         | -                               | 1                                       |
| Fiber Numerical Aperture (NA)(1/e <sup>2</sup> ) |         | 0.05                            | -                                       |
| Mode-field Diameter (μm)(typical)                |         |                                 | 3.5                                     |
| M <sup>2</sup> (ModeMaster beam quality)         |         |                                 | ≤1.1                                    |
| Beam Asymmetry                                   |         |                                 | ≤1:1.1                                  |
| Laser Drive Modes                                |         | Digital, Analog and PC Control  |                                         |
| Digital Modulation                               |         |                                 |                                         |
| Maximum Bandwidth (MHz)                          |         |                                 | 125                                     |
| Rise Time (10% to 90%)(nsec)                     |         |                                 | <2                                      |
| Fall Time (10% to 90%)(nsec)                     |         |                                 | <3                                      |
| Modulation Depth (extinction ratio)              |         |                                 | >1,000,000:1 at 0 Hz, >250:1 at 150 MHz |
| Analog Modulation                                |         |                                 |                                         |
| Maximum Bandwidth (kHz)                          |         |                                 | 350                                     |
| Rise or Fall Time (10% to 90%)(μsec)             |         |                                 | <1                                      |
| Modulation Depth (extinction ratio)              |         |                                 | >10,000:1                               |
| Laser Enable Control Input                       |         |                                 |                                         |
| Maximum Bandwidth (kHz)                          |         |                                 | 130                                     |
| Rise or Fall Time (10% to 90%)(μsec)             |         |                                 | <1                                      |
| Modulation Depth (extinction ratio)              |         |                                 | ∞, complete On/Off                      |
| RMS Noise                                        |         |                                 |                                         |
| 20 Hz to 10 MHz                                  |         |                                 | <0.2%                                   |
| 10 MHz to 500 MHz                                |         |                                 | <1.0%                                   |
| Power Stability (8 hours, ±3°C)                  |         |                                 | ≤±2%                                    |
| Long-Term Output Power Average                   |         |                                 | ≤5%/1000 hrs.                           |
| Warm-Up Time (minutes)                           |         |                                 | <5                                      |
| Spatial Mode (far field)                         |         |                                 | TEM <sub>00</sub>                       |
| Pointing Stability (μrad/°C)                     |         | -                               | <6                                      |
| Polarization Ratio <sup>3</sup>                  |         |                                 | Minimum 100:1                           |
| Tolerances <sup>4</sup>                          |         |                                 |                                         |
| Beam Position (mm)                               |         | -                               | ≤0.5                                    |
| Beam Angle (mrad)                                |         | -                               | ≤15                                     |
| CDRH Laser Classification                        |         |                                 | Class IIb                               |
| ESD Protection <sup>5</sup>                      |         |                                 | Level 4                                 |

<sup>1</sup> Laser-to-laser tolerance. CUBE 640-30 has center wavelength range of 635 nm to 642 nm, CUBE 640-75 with 635 nm to 644 nm range, CUBE 660 with 652 nm to 665 nm range, CUBE 405, 445 and 488 with ±5 nm range.

<sup>2</sup> Output power is variable via RS-232, USB or analog interface from 1% to 100%. Specifications are valid for 100% to 110% power. Recommended power range is 70% to 100% power. Factory setting 100% to 110% of power specification.

<sup>3</sup> Polarization measurement is max./min. power with polarizer at 32 cm from laser with LM-2 power head at 37 cm from laser.

<sup>4</sup> Static alignment tolerances are relative to center of output beam aperture.

<sup>5</sup> Electro-Static Discharge Standard IEC 1000-4-2, 1995.

<sup>6</sup> Fiber FC/APC connector output not compatible for patchcord-to-patchcord connection.

# CUBE FP

## Diode Laser System

| System Specifications                            | CUBE FP | 488-30FP                        | 488-30FP                                |
|--------------------------------------------------|---------|---------------------------------|-----------------------------------------|
| Part Number                                      |         | 1169223                         | 1169222                                 |
| Output from Fiber                                |         | FC/APC (8° angled) <sup>6</sup> | Collimated Beam                         |
| Fiber Cable Type                                 |         | 3 mm PVC                        | Stainless Steel                         |
| Fiber Cable Length (m)(minimum)                  |         | 1                               | 1                                       |
| Wavelength <sup>1</sup> (nm)                     |         | 488                             | 488                                     |
| Output Power <sup>2</sup> (mW)                   |         | 30                              | 30                                      |
| Beam Diameter (1/e <sup>2</sup> )(mm)(typical)   |         | -                               | 1                                       |
| Beam Divergence (mrad)(typical)                  |         | -                               | 1                                       |
| Fiber Numerical Aperture (NA)(1/e <sup>2</sup> ) |         | 0.05                            | -                                       |
| Mode-field Diameter (μm)(typical)                |         |                                 | 3.5                                     |
| M <sup>2</sup> (ModeMaster beam quality)         |         |                                 | ≤1.1                                    |
| Beam Asymmetry                                   |         |                                 | ≤1:1.1                                  |
| Laser Drive Modes                                |         | Digital, Analog and PC Control  |                                         |
| Digital Modulation                               |         |                                 |                                         |
| Maximum Bandwidth (MHz)                          |         |                                 | 125                                     |
| Rise Time (10% to 90%)(nsec)                     |         |                                 | <2                                      |
| Fall Time (10% to 90%)(nsec)                     |         |                                 | <3                                      |
| Modulation Depth (extinction ratio)              |         |                                 | >1,000,000:1 at 0 Hz, >250:1 at 150 MHz |
| Analog Modulation                                |         |                                 |                                         |
| Maximum Bandwidth (kHz)                          |         |                                 | 350                                     |
| Rise or Fall Time (10% to 90%)(μsec)             |         |                                 | <1                                      |
| Modulation Depth (extinction ratio)              |         |                                 | >10,000:1                               |
| Laser Enable Control Input                       |         |                                 |                                         |
| Maximum Bandwidth (kHz)                          |         |                                 | 130                                     |
| Rise or Fall Time (10% to 90%)(μsec)             |         |                                 | <1                                      |
| Modulation Depth (extinction ratio)              |         |                                 | ∞, complete On/Off                      |
| RMS Noise                                        |         |                                 |                                         |
| 20 Hz to 10 MHz                                  |         |                                 | <0.2%                                   |
| 10 MHz to 500 MHz                                |         |                                 | <1.0%                                   |
| Power Stability (8 hours, ±3°C)                  |         |                                 | ≤±2%                                    |
| Long-Term Output Power Average                   |         |                                 | ≤5%/1000 hrs.                           |
| Warm-Up Time (minutes)                           |         |                                 | <5                                      |
| Spatial Mode (far field)                         |         |                                 | TEM <sub>00</sub>                       |
| Pointing Stability (μrad/°C)                     |         | -                               | <6                                      |
| Polarization Ratio <sup>3</sup>                  |         |                                 | Minimum 100:1                           |
| Tolerances <sup>4</sup>                          |         |                                 |                                         |
| Beam Position (mm)                               |         | -                               | ≤0.5                                    |
| Beam Angle (mrad)                                |         | -                               | ≤15                                     |
| CDRH Laser Classification                        |         |                                 | Class IIIB                              |
| ESD Protection <sup>5</sup>                      |         |                                 | Level 4                                 |

<sup>1</sup> Laser-to-laser tolerance. CUBE 640-30 has center wavelength range of 635 nm to 642 nm, CUBE 640-75 with 635 nm to 644 nm range, CUBE 660 with 652 nm to 665 nm range, CUBE 405, 445 and 488 with ±5 nm range.

<sup>2</sup> Output power is variable via RS-232, USB or analog interface from 1% to 100%. Specifications are valid for 100% to 110% power. Recommended power range is 70% to 100% power. Factory setting 100% to 110% of power specification.

<sup>3</sup> Polarization measurement is max./min. power with polarizer at 32 cm from laser with LM-2 power head at 37 cm from laser.

<sup>4</sup> Static alignment tolerances are relative to center of output beam aperture.

<sup>5</sup> Electro-Static Discharge Standard IEC 1000-4-2, 1995.

<sup>6</sup> Fiber FC/APC connector output not compatible for patchcord-to-patchcord connection.

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

# CUBE FP

## Diode Laser System

| System Specifications | CUBE FP                                          | 640-30FP                                | 640-30FP        | 640-75FP           | 640-75FP        |
|-----------------------|--------------------------------------------------|-----------------------------------------|-----------------|--------------------|-----------------|
|                       | Part Number                                      | 1139604                                 | 1139603         | 1169221            | 1169220         |
|                       | Output from Fiber                                | FC/APC (8° angled)                      | Collimated Beam | FC/APC (8° angled) | Collimated Beam |
|                       | Fiber Cable Type                                 | 3 mm PVC                                | Stainless Steel | 3 mm PVC           | Stainless Steel |
|                       | Fiber Cable Length (m)(minimum)                  | 1                                       | 1               | 1                  | 1               |
|                       | Wavelength <sup>1</sup> (nm)                     | 640                                     | 640             | 640                | 640             |
|                       | Output Power <sup>2</sup> (mW)                   | 30                                      | 30              | 75                 | 75              |
|                       | Beam Diameter (1/e <sup>2</sup> )(mm)(typical)   | -                                       | 1               | -                  | 1               |
|                       | Beam Divergence (mrad)(typical)                  | -                                       | 1               | -                  | 1               |
|                       | Fiber Numerical Aperture (NA)(1/e <sup>2</sup> ) | 0.09                                    | -               | 0.09               | -               |
|                       | Mode-field Diameter (μm)(typical)                | 4.5                                     | -               | 4.5                | -               |
|                       | M <sup>2</sup> (ModeMaster beam quality)         | ≤1.1                                    |                 |                    |                 |
|                       | Beam Asymmetry                                   | ≤1:1.1                                  |                 |                    |                 |
|                       | Laser Drive Modes                                | Digital, Analog and PC Control          |                 |                    |                 |
|                       | Digital Modulation                               |                                         |                 |                    |                 |
|                       | Maximum Bandwidth (MHz)                          | 150                                     |                 |                    |                 |
|                       | Rise Time (10% to 90%)(nsec)                     | <2                                      |                 |                    |                 |
|                       | Fall Time (10% to 90%)(nsec)                     | <2                                      |                 |                    |                 |
|                       | Modulation Depth (extinction ratio)              | >1,000,000:1 at 0 Hz, >250:1 at 150 MHz |                 |                    |                 |
|                       | Analog Modulation                                |                                         |                 |                    |                 |
|                       | Maximum Bandwidth (kHz)                          | 350                                     |                 |                    |                 |
|                       | Rise or Fall Time (10% to 90%)(μsec)             | <1                                      |                 |                    |                 |
|                       | Modulation Depth (extinction ratio)              | >10,000:1                               |                 |                    |                 |
|                       | Laser Enable Control Input                       |                                         |                 |                    |                 |
|                       | Maximum Bandwidth (kHz)                          | 130                                     |                 |                    |                 |
|                       | Rise or Fall Time (10% to 90%)(μsec)             | <1                                      |                 |                    |                 |
|                       | Modulation Depth (extinction ratio)              | ∞, complete On/Off                      |                 |                    |                 |
|                       | RMS Noise                                        |                                         |                 |                    |                 |
|                       | 20 Hz to 10 MHz                                  | <0.2%                                   |                 |                    |                 |
|                       | 10 MHz to 500 MHz                                | <1.0%                                   |                 |                    |                 |
|                       | Power Stability (8 hours, ±3°C)                  | ≤±2%                                    |                 |                    |                 |
|                       | Long-Term Output Power Average                   | ≤3%/1000 hrs.                           |                 |                    |                 |
|                       | Warm-Up Time (minutes)                           | <5                                      |                 |                    |                 |
|                       | Spatial Mode (far field)                         | TEM <sub>00</sub>                       |                 |                    |                 |
|                       | Pointing Stability (μrad/°C)                     | -                                       | <6              | -                  | <6              |
|                       | Polarization Ratio <sup>3</sup>                  | Minimum 100:1                           |                 |                    |                 |
|                       | Tolerances <sup>4</sup>                          |                                         |                 |                    |                 |
|                       | Beam Position (mm)                               | -                                       | ≤0.5            | -                  | ≤0.5            |
|                       | Beam Angle (mrad)                                | -                                       | ≤15             | -                  | ≤15             |
|                       | CDRH Laser Classification                        | Class IIb                               |                 |                    |                 |
|                       | ESD Protection <sup>5</sup>                      | Level 4                                 |                 |                    |                 |

<sup>1</sup> Laser-to-laser tolerance. CUBE 640-30 has center wavelength range of 635 nm to 642 nm, CUBE 640-75 with 635 nm to 644 nm range, CUBE 660 with 652 nm to 665 nm range, CUBE 405, 445 and 488 with ±5 nm range.

<sup>2</sup> Output power is variable via RS-232, USB or analog interface from 1% to 100%. Specifications are valid for 100% to 110% power. Recommended power range is 70% to 100% power. Factory setting 100% to 110% of power specification.

<sup>3</sup> Polarization measurement is max./min. power with polarizer at 32 cm from laser with LM-2 power head at 37 cm from laser.

<sup>4</sup> Static alignment tolerances are relative to center of output beam aperture.

<sup>5</sup> Electro-Static Discharge Standard IEC 1000-4-2, 1995.

# CUBE FP

## Diode Laser System

| System Specifications                            | CUBE FP | 660-75FP                       | 660-75FP                                |
|--------------------------------------------------|---------|--------------------------------|-----------------------------------------|
| Part Number                                      |         | 1151423                        | 1151412                                 |
| Output from Fiber                                |         | FC/APC (8° angled)             | Collimated Beam                         |
| Fiber Cable Type                                 |         | 3 mm PVC                       | Stainless Steel                         |
| Fiber Cable Length (m)(minimum)                  |         | 1                              | 1                                       |
| Wavelength <sup>1</sup> (nm)                     |         | 660                            | 660                                     |
| Output Power <sup>2</sup> (mW)                   |         | 75                             | 75                                      |
| Beam Diameter (1/e <sup>2</sup> )(mm)(typical)   |         | -                              | 1                                       |
| Beam Divergence (mrad)(typical)                  |         | -                              | 1                                       |
| Fiber Numerical Aperture (NA)(1/e <sup>2</sup> ) |         | 0.09                           | -                                       |
| Mode-field Diameter (μm)(typical)                |         |                                | 4.5                                     |
| M <sup>2</sup> (ModeMaster beam quality)         |         |                                | ≤1.1                                    |
| Beam Asymmetry                                   |         |                                | ≤1:1.1                                  |
| Laser Drive Modes                                |         | Digital, Analog and PC Control |                                         |
| Digital Modulation                               |         |                                |                                         |
| Maximum Bandwidth (MHz)                          |         |                                | 150                                     |
| Rise Time (10% to 90%)(nsec)                     |         |                                | <2                                      |
| Fall Time (10% to 90%)(nsec)                     |         |                                | <2                                      |
| Modulation Depth (extinction ratio)              |         |                                | >1,000,000:1 at 0 Hz, >250:1 at 150 MHz |
| Analog Modulation                                |         |                                |                                         |
| Maximum Bandwidth (kHz)                          |         |                                | 350                                     |
| Rise or Fall Time (10% to 90%)(μsec)             |         |                                | <1                                      |
| Modulation Depth (extinction ratio)              |         |                                | >10,000:1                               |
| Laser Enable Control Input                       |         |                                |                                         |
| Maximum Bandwidth (kHz)                          |         |                                | 130                                     |
| Rise or Fall Time (10% to 90%)(μsec)             |         |                                | <1                                      |
| Modulation Depth (extinction ratio)              |         |                                | ∞, complete On/Off                      |
| RMS Noise                                        |         |                                |                                         |
| 20 Hz to 10 MHz                                  |         |                                | <0.2%                                   |
| 10 MHz to 500 MHz                                |         |                                | <1.0%                                   |
| Power Stability (8 hours, ±3°C)                  |         |                                | ≤±2%                                    |
| Long-Term Output Power Average                   |         |                                | ≤3%/1000 hrs.                           |
| Warm-Up Time (minutes)                           |         |                                | <5                                      |
| Spatial Mode (far field)                         |         |                                | TEM <sub>00</sub>                       |
| Pointing Stability (μrad/°C)                     |         | -                              | <6                                      |
| Polarization Ratio <sup>3</sup>                  |         |                                | Minimum 80:1                            |
| Tolerances <sup>4</sup>                          |         |                                |                                         |
| Beam Position (mm)                               |         | -                              | ≤0.5                                    |
| Beam Angle (mrad)                                |         | -                              | ≤15                                     |
| CDRH Laser Classification                        |         |                                | Class IIIB                              |
| ESD Protection <sup>5</sup>                      |         |                                | Level 4                                 |

<sup>1</sup> Laser-to-laser tolerance. CUBE 640-30 has center wavelength range of 635 nm to 642 nm, CUBE 640-75 with 635 nm to 644 nm range, CUBE 660 with 652 nm to 665 nm range, CUBE 405, 445 and 488 with ±5 nm range.

<sup>2</sup> Output power is variable via RS-232, USB or analog interface from 1% to 100%. Specifications are valid for 100% to 110% power. Recommended power range is 70% to 100% power. Factory setting 100% to 110% of power specification.

<sup>3</sup> Polarization measurement is max./min. power with polarizer at 32 cm from laser with LM-2 power head at 37 cm from laser.

<sup>4</sup> Static alignment tolerances are relative to center of output beam aperture.

<sup>5</sup> Electro-Static Discharge Standard IEC 1000-4-2, 1995.

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories

# CUBE FP

## Diode Laser System

### Utility and Environmental Requirements

|                                                      |                             |
|------------------------------------------------------|-----------------------------|
| Laser Head Operating Voltage                         | +4.8 to 6.5 VDC             |
| Laser Head Operating Current                         | <2.5 amps                   |
| Baseplate Temperature Range <sup>1</sup>             | +10° to 50°C (50° to 122°F) |
| Maximum Heat Dissipation of Head (baseplate at 50°C) | 13W                         |
| Storage Temperature                                  | -20° to 60°C (-4° to 140°F) |
| Power Supply <sup>2</sup> (included)                 |                             |
| Operating Voltage <sup>3</sup>                       | 100 to 240 VAC, 50 to 60 Hz |
| Power Consumption                                    | <15W                        |

### Dimensions (L x W x H)

### Weight

|                                           |                                        |                    |
|-------------------------------------------|----------------------------------------|--------------------|
| Laser Head                                | 100 x 40 x 40 mm (3.9 x 1.6 x 1.6 in.) | 0.3 kg (0.7 lbs.)  |
| Laser Head with Minimum Fiber Bend Radius | 203 mm (8.0 in.)                       |                    |
| Fiber Tensile Load (max.)                 |                                        | 1 kg (2.2 lbs.)    |
| Fiber Minimum Bend Radius                 | 51 mm (2.0 in.)                        |                    |
| Control Box                               | 138 x 62 x 26 mm (5.4 x 2.4 x 1.1 in.) | 0.12 kg (0.3 lbs.) |
| Power Supply                              | 87 x 47 x 32 mm (3.4 x 1.9 x 1.3 in.)  | 0.22 kg (0.5 lbs.) |
| Head Cable                                | 254 to 406 mm (10 to 18 in.)           | 0.08 kg (0.2 lbs.) |
| USB Cable                                 | 1.8 m (6 ft.)                          | 0.1 kg (0.2 lbs.)  |
| Heat Sink (optional - sold separately)    | 127 x 66 x 44 mm (5 x 2.6 x 1.74 in.)  | 0.4 kg (0.9 lbs.)  |

<sup>1</sup> Non-condensing with laser diode TEC at set point of 22°C.

<sup>2</sup> On-off power switch and LED included in cord to laser head.

<sup>3</sup> System includes USA-type power cord, IEC320 input connection.



Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories

# CUBE FP


## Diode Laser System

|                         |           |                                           |                                                                          |
|-------------------------|-----------|-------------------------------------------|--------------------------------------------------------------------------|
| <b>CUBE Accessories</b> | 1073840   | CUBE Heat Sink Accessory                  | CUBE heat sink with fan, 127 x 66 x 44 mm                                |
|                         | 1079150   | CUBE 2nd Emission Lamp                    | CDRH lamp for laser emission. Connects to CUBE I/O cable.                |
|                         | 1072166   | *CUBE Interface Cable                     | CUBE to control box, I/O, RS-232 and power                               |
|                         | 1072454   | *CUBE Power Supply, 6 VDC, 2.5 amps       | IEC320 input, 110V/220V, switch and LED                                  |
|                         | 1108906   | *CUBE USB Cable, 6 ft. (1.8m)             | USB standard to mini for CUBE                                            |
|                         | 1039966   | * Control Box, CDRH, Keyswitch, Interlock | Control box for CUBE lasers. Included box, key, interlock RCA and cable. |
|                         | 1080090   | Cable, Extension, DB-9 M/F, 6 ft. (1.8m)  | Extension cable to extend control box or RS-232. M-F, DB-9.              |
|                         | 1040408   | * Plug, RCA                               | Interlock RCA plug for control box, non-shorted                          |
|                         | 1150553   | *CUBE FP Manual                           | CUBE FP user manual and software CD for Windows 2000** and Windows XP**  |
|                         | 6005-0059 | *Power Cord, USA, 6 ft. (1.8m)            | Power cord, USA plug style to IEC320                                     |


\* Included free with every CUBE system.

\*\* Windows 2000 and Windows XP are registered trademarks of Microsoft Corporation.


**CUBE  
Heat Sink**

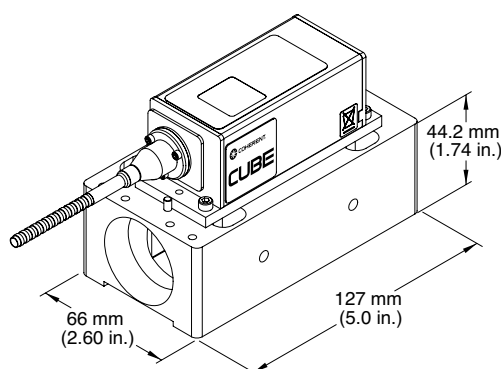


**Heat Sink  
Top View**

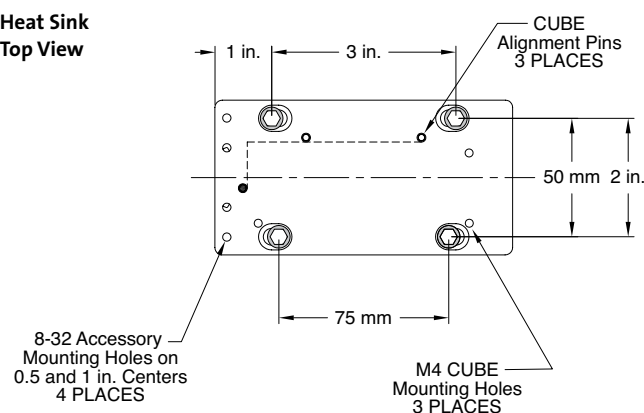

**CUBE  
Control Box**

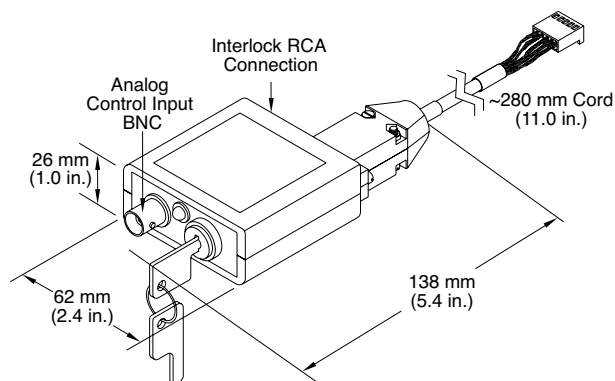

**CUBE  
Power Supply**

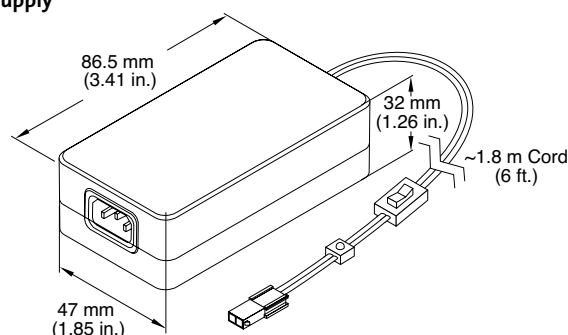

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

# Radius

## Diode Laser Systems

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories



Radius lasers are full-featured, high-performance, CW diode laser systems that include a laser head, power supply and optional control box. Designed for convenient mounting into a HeNe laser holder for fast and simple integration, they deliver power, stability and performance in a complete package with the best value.

### Adjustable Focus Radius

These lasers allow faster and easier system integration. When optimizing the beam for best focus or fiber launch, the adjustable focus eliminates the need for an expensive beam expander. Instead, a tool (included) adjusts the beam divergence through the front shutter port.

The optional control box—with a key-switch, 5-second power-on delay, emission indicator and interlock—allows the system to comply with CDRH safety regulations. The control box plugs into the side port of the Radius.

The side port can also be used to modulate the laser on-and-off as fast as 5 Hz (100 Hz for the Radius 635) through the interlock control input. This modulation can be used to extend the lifetime of the diode by enabling the laser to be 'on' only when needed—while the thermal system remains constant for improved stability.

Mounted on a large TEC, the diode laser and optics remain at a constant temperature of 22°C. A heat sink with an internal fan removes the heat. The closed-loop, constant temperature keeps the wavelength constant, minimizes mode-hops and maintains power stability.

### Extra Photodiode

The Radius 405-30EP features an extra photodiode that provides the best power stability over time.

### Circular Beam

The patented micro-lens technology and aspherical optics of the Radius 635-25 circularize the output of the diode laser to better than 80% with a typical  $M^2 < 1.2$ .

### Side Port

Using a side port, users can vary the diode temperature in order to change the center wavelength by  $\pm 1$  nm (635 nm diodes).

Varying the temperature adjustment also allows users to avoid modehops, thus providing the maximum coherence length.



# Radius

## Diode Laser Systems

### Key Features

- HeNe cylindrical-style package
- All-in-one design: Optics, electronics and TEC
- Adjustable beam option
- ESD protection
- Internal cooling fan
- CDRH Class IIIb safety compliance option
- Constant power CW

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

| System Specifications       | Radius                                                  | 405-50                                    | 635-25           | 635-25E         |
|-----------------------------|---------------------------------------------------------|-------------------------------------------|------------------|-----------------|
|                             | Wavelength (nm)                                         | 405 ±5                                    | 635 +7/-2        | 635 +7/-2       |
|                             | Power (mW)                                              | 50                                        | 25               | 25              |
|                             | Max. Modulation Frequency (Hz)                          | 5                                         | 100              | 100             |
|                             | RMS Noise (10 Hz to 10 MHz)                             | <0.5%                                     | <0.5%            | <0.5%           |
|                             | Power Stability (8 Hours)                               | <5%                                       | <2%              | <2%             |
|                             | Adjustable Beam Focus (Divergence)                      | Yes                                       | No               | Yes             |
|                             | Static Alignment Tolerances <sup>1</sup>                | <5 mrad (angle)                           | <15 mrad (angle) | <5 mrad (angle) |
|                             | Polarization (Horizontal)                               | >100:1                                    | >100:1           | >100:1          |
|                             | CDRH Class                                              | Class IIIb                                | Class IIIb       | Class IIIb      |
|                             | Spot Size (mm)(typical)                                 | 4.7 x 1.6                                 | 1                | 1.5 x 4.6       |
|                             | Divergence (mrad)(typical)                              | 0.2 x 0.3                                 | 1.3              | 0.5 x 0.2       |
|                             | Part Number for System with Control Box included        | 1064507                                   | 1051385          | —               |
|                             | Part Number for System without Control Box <sup>2</sup> | 1064506                                   | 1008053          | 1149819         |
| Laser Head Specifications   | Laser Drive                                             | CW, Constant Power                        |                  |                 |
|                             | Pointing Stability                                      | <6 µrad/°C                                |                  |                 |
|                             | Beam Spatial Mode (Far Field)                           | TEM <sub>00</sub>                         |                  |                 |
|                             | ESD Protection <sup>4</sup>                             | Level 4                                   |                  |                 |
|                             | Dimensions (Diameter x Length)                          | 44.5 x 218 mm (1.75 x 8.6 in.)            |                  |                 |
|                             | Optional Control Box (L x W x H)                        | 138 x 62 x 26 mm (5.4 x 2.4 x 1 in.)      |                  |                 |
|                             | Ambient Temperature Range for Operation                 | 10°C to 40°C (50°F to 104°F) <sup>5</sup> |                  |                 |
|                             | Storage Temperature                                     | -20°C to 60°C (-4°F to 140°F)             |                  |                 |
| Power Supply Specifications | Maximum Heat Load from Laser Head                       | 10W, typical 3W                           |                  |                 |
|                             | Power Supply included with System                       | Yes                                       |                  |                 |
|                             | Weight                                                  | 570 g (1.25 lbs.)                         |                  |                 |
|                             | Dimensions (L x W x H)                                  | 146 x 76 x 43 mm (5.8 x 3.0 x 1.7 in.)    |                  |                 |
|                             | Output to Laser Head                                    | +5VDC, -5VDC, +12VDC                      |                  |                 |
|                             | Cord Length to Laser Head                               | 120 cm (48 in.)                           |                  |                 |
|                             | AC Input Voltage                                        | 100 to 240VAC, 47 to 63 Hz <sup>6</sup>   |                  |                 |
|                             | AC Input Power                                          | <15W                                      |                  |                 |

<sup>1</sup> Static alignment tolerances are relative to the laser's outside diameter.

<sup>2</sup> Systems that do not include the control box do not comply with the Federal Regulations (21 CFR Subchapter J) as administered by the Center for Devices and Radiological Health. Purchaser acknowledges that his/her products must comply with these regulations before they can be sold to a customer.

<sup>3</sup> Modulation through side panel I/O port. Via TTL signal with 400 KOhm input impedance (VH>4V, VL < 3V);

TTL needs to be applied for minimum of 10 msec for Radius 635 and 100 msec for 375/405/440.

<sup>4</sup> Electro-Static Discharge Standard IEC 1000-4-2, 1995.

<sup>5</sup> Non-condensing.

<sup>6</sup> Shipped with USA power cord, IEC320 input connection.

# Radius


## Diode Laser Systems

### HeNe Comparison


Radius 635 specifications are equal to or better than most HeNe lasers in most applications. Radius lasers also cost less and come in a smaller package. Both power stability and pointing stability are superior with the Radius. Warm-up time and output power over the laser's life are also far superior when compared to the HeNe. For all the details, please see the Radius 635 Competitive Comparison, found on the Radius page of our website, [www.Coherent.com](http://www.Coherent.com).

|                           |         |                                           |                                                                            |
|---------------------------|---------|-------------------------------------------|----------------------------------------------------------------------------|
| <b>Radius Accessories</b> | 1039966 | Control Box, CDRH, Keyswitch, Interlock   | Control Box for Radius lasers. Included box, key, interlock RCA and cable. |
|                           | 1080090 | Cable, Extension, DB-9 M/F, 6 ft. (1.8 m) | Extension cable to extend control box. Male-Female, DB-9.                  |
|                           | 1040408 | Plug, RCA                                 | Interlock RCA plug for Control Box, non-shorted                            |
|                           | 1059300 | Power Supply, LDM, 110/220VAC             | Radius or ULN, Desktop, 43 Watt, Triple Output                             |


### Radius Laser Head

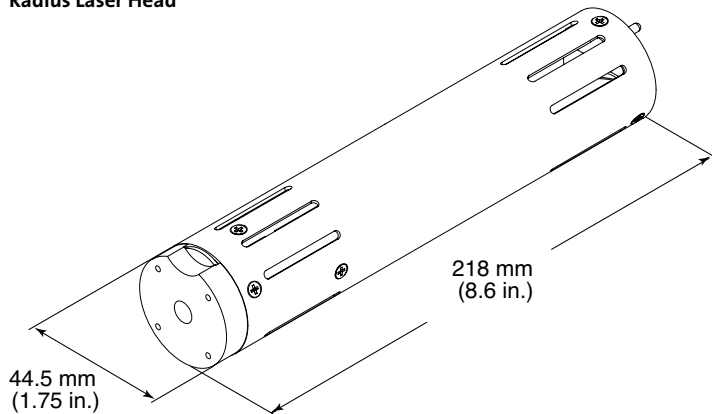

### Radius Power Supply

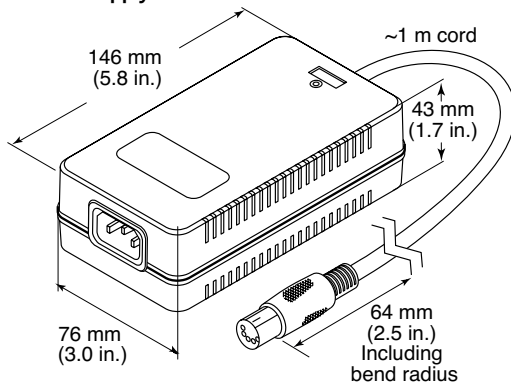

### Radius Control Box

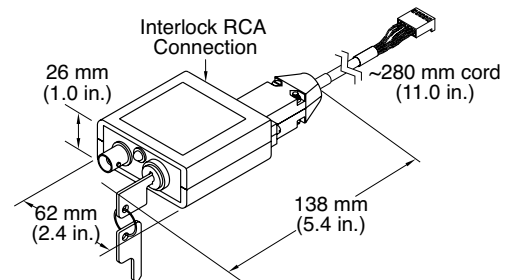

# ULN

## Diode Laser Modules



Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

ULN-Series diode laser modules are designed for applications that require particularly low noise or mode-hop, noise-free operation. Sophisticated drive electronics are used to ensure low noise output. Typical RMS noise is 0.06% or better for detector bandwidth from 10 Hz to 10 MHz, with practically no warm-up time.

MicroBlaze™ beam circulation technology produces a high-brightness, diffraction-limited circular beam that is fully collimated with low divergence.

Overall, these low-noise, compact LabLasers are ideal replacements for HeNe lasers, which require higher power and generate more heat.

Applications include genomics, particle counting, flow cytometry, optical storage, confocal microscopy, high throughput drug screening, and spectroscopy. They are also ideal for biomedical instrument design.

The ULN-Series features an alternative solution to the mode-hop problem. Rather than regulate temperature, a special modulation is used to force the laser into a multi-longitudinal mode.

While a normal laser will operate in a single-mode longitudinal or oscillate between two strong modes, the modulation creates several modes of lower

intensity. As the laser temperature increases, these modes move like a caterpillar across the wavelength spectrum. The movement does not allow abrupt changes, so the system operates as if there were no mode-hops.

The result is low RMS noise (~0.06%), which is stable over changes in temperature and the life of the diode laser module.

- ULN stabilizes both power and noise during temperature changes (note that wavelength varies with the standard temperature coefficient of the diode, ~0.2 nm/°C).
- ULN controls noise over the lifetime of the laser and will not develop mode-hop noise as the laser ages.
- ULN does not change the spatial qualities of the laser beam, so the beam can be focused to the same spot size, profile and shape.

The ULN laser is sold with a triple output power supply to provide the necessary +5V and +12V to drive the laser.

Use the LabLaser mount [part no. 0221-449-00] for convenient heat sinking and mounting with adjustable alignment.

# ULN

## Diode Laser Modules

### Key Features

- Ultra-low noise (ULN)
- RMS noise <0.06% for bandwidths of 10 Hz to 10 MHz
- Circular beam
- No mode-hop noise
- No warm-up period
- Wavelength 635 nm with 5 mW CW power
- Long lifetime


**ULN circular diode laser modules** create a 1 mm circular output beam using MicroBlaze patented technology. This is a technique that accurately positions a fiber microlens inside the diode laser cap. The microlens uniquely corrects the divergence of the fast axis. Once corrected, the beam shape changes from elliptical to circular, providing maximum power transmission and less astigmatism than other methods. The circular low-divergence beams have almost no astigmatism and a total wave-front error typically less than  $\lambda/10$ . The M2 beam quality factor is typically 1.2 or better. These diode laser modules produce a beam of such high quality that they can replace more expensive HeNe lasers in many applications.

### Compact, Rugged and Reliable

Our miniature laser modules are superb choices for many applications demanding affordable, compact lasers. Using AlignLock technology, all modules are optimized for thermal conductivity, protected from mechanical shock and deliver superior pointing accuracy and stability. These compact, rugged and reliable diode laser sources also have a sophisticated internal driver/controller circuit, which provides both surge and reverse polarity protection. The nominal operating temperature range is -10 to +40°C. Diode lifetime will be greater at lower temperatures.

| System Specifications | LabLaser™ | 635 nm 5 mW C ULN          |
|-----------------------|-----------|----------------------------|
| Wavelength (nm)       |           | 635 $\pm 7/-2$             |
| Power (mW)            |           | 5                          |
| CDRH Class            |           | Class IIIb, Non-Conforming |
| Spot Size (mm)        |           | 1                          |
| Divergence (mrad)     |           | 1                          |
| Laser Drive           |           | ULN, Ultra-Low Noise       |
| Power Supply          |           | Included                   |
| Connector             |           | 5-pin DIN                  |
| Part Number           |           | 31-0144-000                |

LabLaser

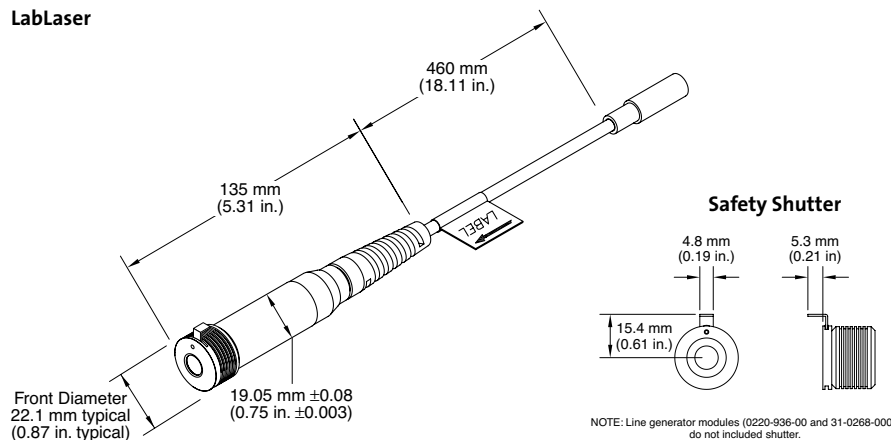

NOTE: Line generator modules (0220-936-00 and 31-0268-000) do not include shutter.

# Cross-Hair Line Generators




## Key Features

- Focusable VLM2-style package
- High visibility 635 nm
- Adjustable line width and focus
- Dual-element optic design
- High-quality diffractive optic

Cross-Hair Line Generators deliver two sharp, high-brightness, fan-angle lines at 635 nm. These adjustable focus cross-hair generators are ideal for long-life projection and industrial applications. The adjustable focus allows the width of the line to be focused, while the length of the line remains fixed. The line can typically be focused as small as  $<0.1$  mm at a distance of 10 cm or can be focused to as small as  $\sim 0.4$  mm at a distance of 1.0 meter.

The front knob allows for easy adjustment of the line width and focus location. During adjustment the optics do not rotate and focusing action is even and uniform while beam alignment is maintained<sup>3</sup>. A locking ring is used to secure the focus in any position. Use the VLM2 mount [part no. 0221-437-00] for convenient heat sinking and mounting with adjustable alignment. For lasers without a connector, order plug no.1057734.

## FVLM2

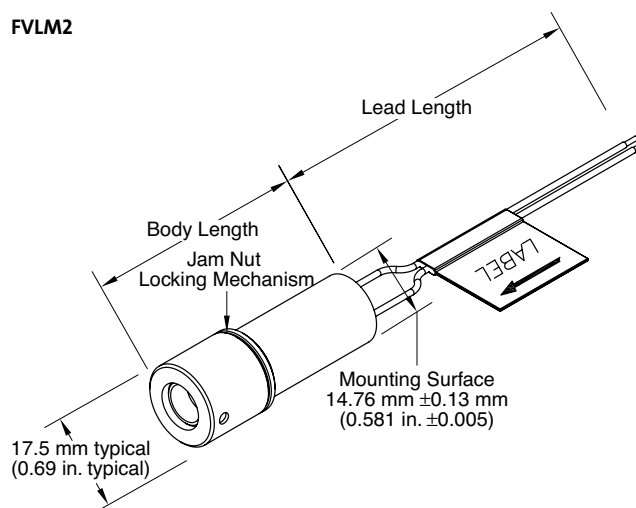

| System Specifications     | FVLM2                   | Cross-Hair 635 nm 0.6 mW | Cross-Hair 635 nm 1.7 mW |
|---------------------------|-------------------------|--------------------------|--------------------------|
| Wavelength (nm)           | 635 $\pm 7/-2$          | 635 $\pm 7/-2$           | 635 $\pm 7/-2$           |
| Power (mW)                | 0.55                    | 0.55                     | 1.7                      |
| CDRH Class                | Class II                | Class II                 | Class II                 |
| Spot Size (mm)            | 1 Adjustable            | 1 Adjustable             | 1 Adjustable             |
| Divergence (mrad)         | 24° x 1 mrad Adjustable | 24° x 1 mrad Adjustable  | 24° x 1 mrad Adjustable  |
| Laser Drive <sup>1</sup>  | CW                      | CW                       | CW                       |
| Power Supply <sup>2</sup> | 5 to 10VDC              | 5 to 10VDC               | 5 to 10VDC               |
| Body Length (mm)          | 48 to 55                | 48 to 55                 | 48 to 55                 |
| Lead Length (mm)          | 610                     | 610                      | 610                      |
| Connector                 | None                    | None                     | None                     |
| Part Number               | 0222-210-00             | 0222-209-00              | 0222-209-00              |

<sup>1</sup> Laser Drive is Continuous Wave (CW) with constant power.

<sup>2</sup> Power Supply: For modules with 5-10VDC it is recommended that the module be operated at 5V.

<sup>3</sup> Note the focus adjustment may slightly affect the beam-pointing direction. Due to the numerical aperture of the lens, slight variations in beam output power may occur as the lens is adjusted.

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

# Focusable Line Generator Modules

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories



Focusable diode laser modules combine the technology of high-quality diode lasers with superior optics, sophisticated electronics and rugged alignment packaging. These modules can also be used in many applications to reduce the requirements for external optics. The result is an overall system that is smaller, simpler and less costly.

The adjustable focus is a great benefit for many applications. For machine vision, particle measurements and inspection—and their inherent need for direct measurements—the beam can be focused at a distance for the smallest spot size to optimize the system resolution.

The front knob allows for easy adjustment. During adjustment the optics do not rotate and focusing action is even and uniform while beam alignment is maintained<sup>1</sup>. A locking ring is used to secure the focus in any position. Use the VLM2 mount [part no. 0221-437-00] for convenient heat sinking and mount-

ing with adjustable alignment. For lasers without a connector, order plug no. 1057734.

## Compact, Rugged and Reliable

Our miniature laser modules are superb choices for many applications demanding affordable, compact lasers. Using AlignLock technology, all modules are optimized for thermal conductivity, protected from mechanical shock and deliver superior pointing accuracy and stability. These compact, rugged and reliable diode laser sources also have a sophisticated internal driver/controller circuit, which provides both surge and reverse polarity protection. The nominal operating temperature range is -10 to +40°C. Diode lifetime will be greater at lower temperatures.

**Focusable line generators** deliver sharp lines with adjustable line width, while the length of the line remains fixed. The line can be typically focused as small as <0.1 mm at a distance of 10 cm or can be focused to as small as ~0.4 mm at a distance of 1 meter.

<sup>1</sup> The focus adjustment may slightly affect the beam-pointing direction. Due to the numerical aperture of the lens, slight variations in beam output power may occur as the lens is adjusted.

# Focusable Line Generator Modules

## Key Features


- Focusable VLM2 style package
- High visibility 635 nm
- Adjustable divergence and focus

## Advantages of the Focus Feature

The focus feature can be used to:

- Adjust the beam propagation so that it diverges to create a larger spot at the target for broader illumination.
- Adjust the beam so that it converges to create a specific focus spot size at a certain working distance.
- Adjust the beam so that the beam waist is at a particular point in space.
- Adjust the beam to move the focus spot location when an external objective lens is used. This is especially useful when matching multiple wavelengths into the same beam focus position.
- Adjust the beam at the laser to reduce the need for expensive and bulky adjustable beam expanders.

FVLM2

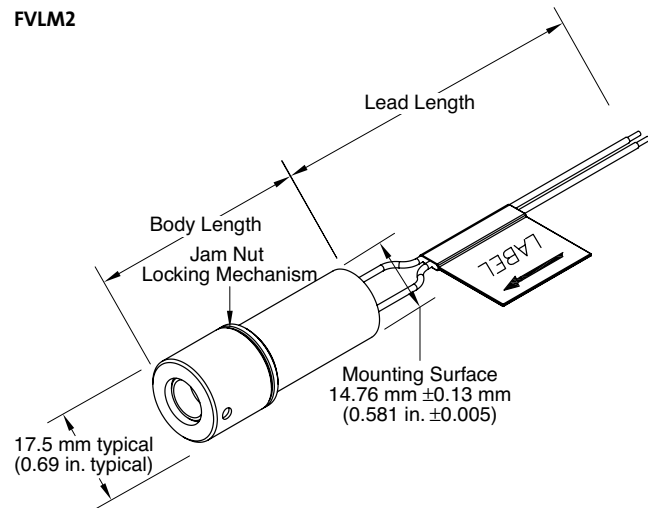

| System Specifications | FVLM2                     | 635 nm 1.8 mW 60°       |
|-----------------------|---------------------------|-------------------------|
|                       | Wavelength (nm)           | 635 +7/-2               |
|                       | Power (mW)                | 1.8                     |
|                       | CDRH Class                | Class II                |
|                       | Spot Size (mm)            | 1 Adjustable            |
|                       | Divergence (mrad)         | 60° x 1 mrad Adjustable |
|                       | Laser Drive               | CW                      |
|                       | Power Supply <sup>1</sup> | 5 to 10VDC              |
|                       | Body Length (mm)          | 48 to 54                |
|                       | Lead Length (mm)          | 610                     |
|                       | Connector                 | none                    |
|                       | Part Number               | 0220-846-00             |

<sup>1</sup> Power Supply: For modules with 5-10VDC it is recommended that the module be operated at 5V.

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories

# Line Generator Diode Laser Modules

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories



Line generators from Coherent project a clean, sharp line ~1 mm wide, with the line length represented by the specific fan angle. With fan angles of 40 or 60 degrees, a wavelength of 635 nm and various output powers available, line generators from Coherent offer the options you need for your particular application.

Use the appropriate VLM2 mount [part no. 0221-437-00]; or VLM3 mount [part no. 0221-448-00] for convenient heat sinking and mounting with adjustable alignment.

For lasers without a connector, order plug no. 1057734.


## Compact, Rugged and Reliable

Our miniature laser modules are superb choices for many applications demanding affordable, compact lasers. Using AlignLock technology, all modules are optimized for thermal conductivity, protected from mechanical shock and deliver superior pointing accuracy and stability. These compact, rugged and reliable diode laser sources also have a sophisticated internal driver/controller circuit, which provides both surge and reverse polarity protection. The nominal operating temperature range is -10 to +40°C. Diode lifetime will be greater at lower temperatures.

## Line Generator Line Length

The Coherent Line Generator projects a sharp line that is ~1 mm wide at the module face and ~1.5 mm wide at one meter. The length of a line is traditionally defined by specifying the visual fan angle, which is 85° for these modules. The 1/e<sup>2</sup> intensity fan angle is about 65°.

| Fan Angle<br>$\theta$ | Line Length L |           |         |         |
|-----------------------|---------------|-----------|---------|---------|
|                       | d = 0.25 m    | d = 0.5 m | d = 1 m | d = 2 m |
| 60°                   | 0.29 m        | 0.58 m    | 1.16 m  | 2.31 m  |
| 40°                   | 0.18 m        | 0.36 m    | 0.73 m  | 1.46 m  |

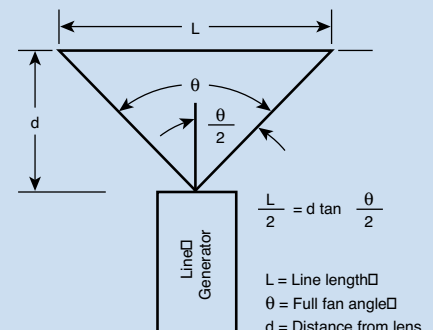

# Line Generator Diode Laser Modules

## Key Features



- VLM2 and VLM3-style packages
- Optional adjustable line width and focus
- 40° or 60° fan angles
- High brightness 635 nm

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

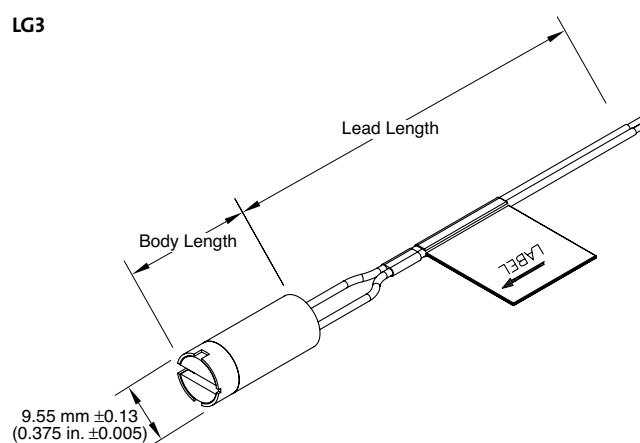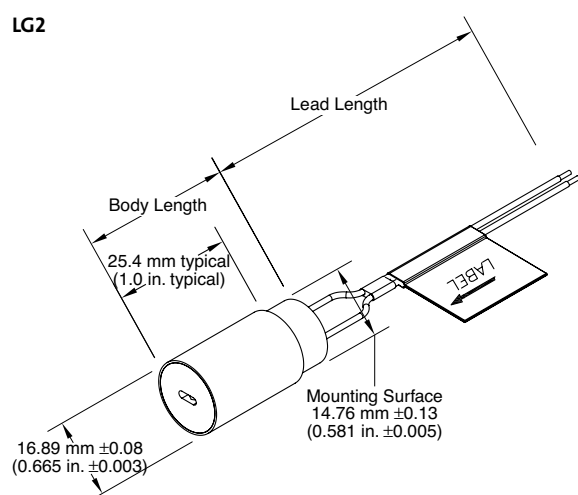

| System Specifications     | FVLM2 | LG3 635 nm 1.5 mW 40° | LG2 635 nm 2 mW 60° |
|---------------------------|-------|-----------------------|---------------------|
| Wavelength (nm)           |       | 635 +7/-2             | 635 +7/-2           |
| Power (mW)                |       | 1.5                   | 2                   |
| CDRH Class                |       | Class II              | Class II            |
| Spot Size (mm)            |       | 1                     | 1                   |
| Divergence (mrad)         |       | 40° x 1               | 60° x 1             |
| Laser Drive <sup>1</sup>  |       | CW                    | CW                  |
| Power Supply <sup>2</sup> |       | 3 to 6VDC             | 5 to 10VDC          |
| Body Length (mm)          |       | 24                    | 35                  |
| Lead Length (mm)          |       | 914                   | 914                 |
| Connector                 |       | none                  | none                |
| Part Number               |       | 0222-037-00           | 0220-934-00         |

<sup>1</sup> Laser Drive is Continuous Wave (CW) with constant power. For MVP option the laser includes a 3rd wire for Modulation and Variable Power Control

<sup>2</sup> Power Supply: For modules with 5-10VDC it is recommended that the module be operated at 5V. For modules with 3-6VDC it is recommended that the module be operated at 3.3V or 5V.


# Miniature Diode Laser Modules

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories

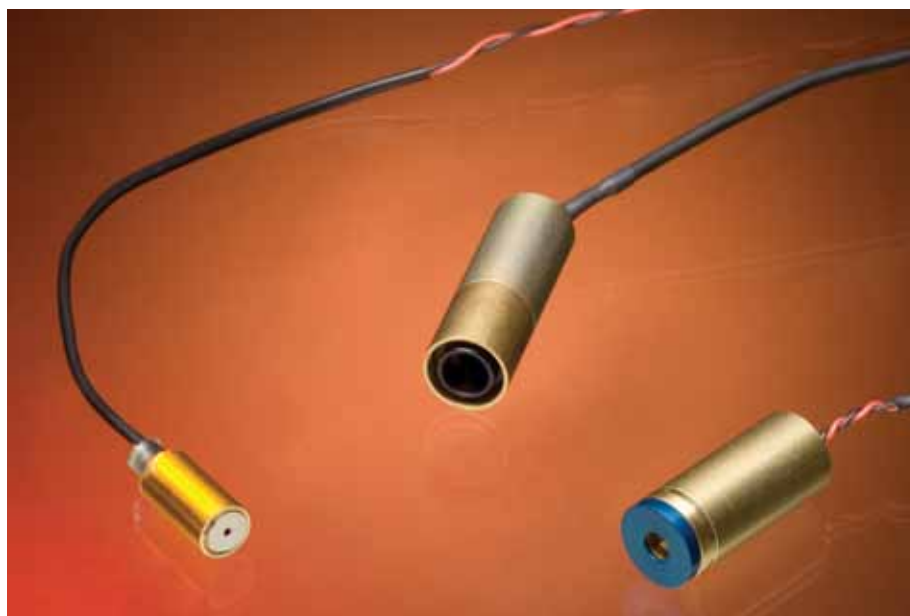

Circular and elliptical miniature diode laser modules provide the smallest package size and best value. Included inside every module are complete electronics to manage the constant power laser drive with a closed-light loop based on an internal photodiode. Also included in the electronics are protection circuits to provide a robust system in even the smallest of packages.

Miniature diode laser modules use the patented AlignLock technology to guarantee beam-pointing alignment of  $<4.3$  mrad. This allows for simpler and easier integration, as you can trust the unit-to-unit consistency of each beam.

Use the appropriate VLM2 mount [part no. 0221-437-00] or the VLM3 mount [part no. 0221-448-00] for convenient heat sinking and mounting with adjustable alignment.

For lasers without a connector, order plug no. 1057734.

**Miniature VLM2 circular diode laser modules** create a 1 mm circular output beam using MicroBlaze patented technology. This is a technique that accurately positions a fiber micro-lens inside the diode laser cap. The micro-lens uniquely corrects the divergence of the fast axis. Once corrected, the beam shape changes

from elliptical to circular, providing maximum power transmission and less astigmatism than other methods.

**Miniature elliptical diode laser modules** offer exceptional beam quality with a typical  $M2$  of 1.4 or better in a typical  $1 \times 3$  aspect ratio.

## Compact, Rugged and Reliable

Our miniature laser modules are superb choices for many applications demanding affordable, compact lasers. Using AlignLock technology, all modules are optimized for thermal conductivity, protected from mechanical shock and deliver superior pointing accuracy and stability. These compact, rugged and reliable diode laser sources also have a sophisticated internal driver/controller circuit, which provides both surge and reverse polarity protection. The nominal operating temperature range is  $-10$  to  $+40^{\circ}\text{C}$ . Diode lifetime will be greater at lower temperatures.

# Miniature Diode Laser Modules


## Key Features

- VLM2 and VLM3-style packages
- High visibility 635 nm, 650 nm, and 670 nm
- Circular or elliptical beams


## Circular 1 mm Beam

Circular diode laser modules produce a ~1 mm near-diffraction-limited beam. The circular lowdivergence beams have almost no astigmatism and a total wave-front error typically less than  $\lambda/10$ . The M2 beam quality factor is typically 1.2 or better. These diode laser modules produce a beam of such high quality that they can replace more expensive HeNe lasers in many applications. Circular diode laser modules are also smaller, require less power and generate less heat than traditional HeNe laser sources.

VLM2

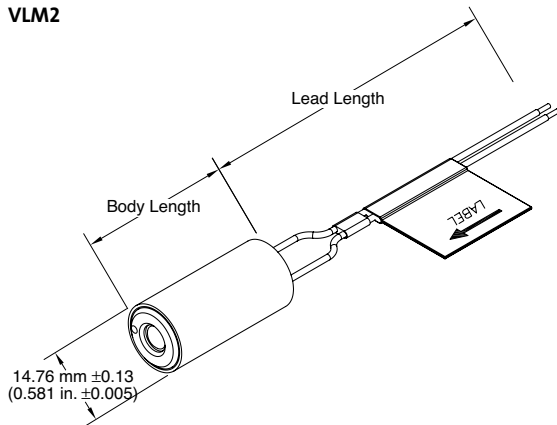

VLM3

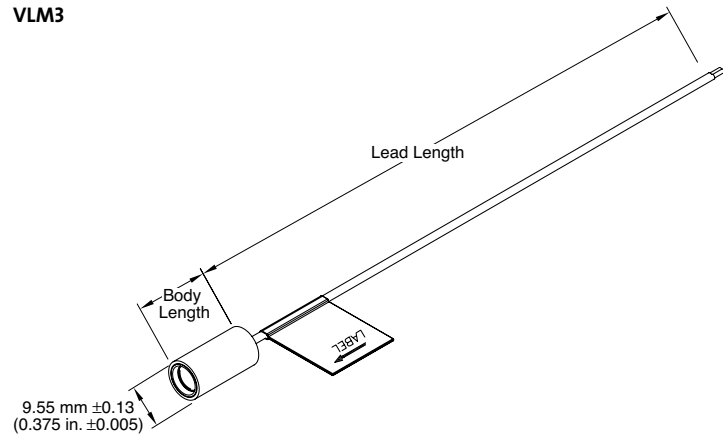

| System Specifications     | FVLM2 | VLM2 635 nm<br>4 mW E | VLM3 670 nm<br>4 mW E | VLM2 670 nm<br>7.5 mW E      | VLM2 635 nm<br>4 mW C 1 mm | VLM3 650 nm<br>1.6 mW C 1.3 mm |
|---------------------------|-------|-----------------------|-----------------------|------------------------------|----------------------------|--------------------------------|
| Wavelength (nm)           |       | 635 +7/-2             | 670 ±5                | 670 ±5                       | 635 +7/-2                  | 650 ±5                         |
| Power (mW)                |       | 4.2                   | 4.2                   | 7.5                          | 4.2                        | 1.6                            |
| CDRH Class                |       | Class IIIa            | Class IIIa            | Class IIIb<br>Non-Conforming | Class IIIa                 | Class IIIa<br>Non-Conforming   |
| Spot Size (mm)            |       | 5.8 x 1.7             | 4.6 x 1.3             | 3.4 x 1.3                    | 1.3                        | 1.3                            |
| Divergence (mrad)         |       | 0.2 x 0.6             | 0.2 x 0.7             | 0.3 x 0.8                    | 0.7                        | 0.8                            |
| Laser Drive <sup>1</sup>  |       | CW                    | CW                    | CW                           | CW                         | CW                             |
| Power Supply <sup>2</sup> |       | 5 to 10VDC            | 3 to 6VDC             | 5 to 10VDC                   | 5 to 10VDC                 | 5 to 6VDC                      |
| Body Length (mm)          |       | 32                    | 19                    | 28                           | 34                         | 26                             |
| Lead Length (mm)          |       | 914                   | 914                   | 305                          | 914                        | 914                            |
| Connector                 |       | none                  | none                  | none                         | none                       | none                           |
| Part Number               |       | 0221-534-00           | 0220-968-00           | 0223-038-00                  | 0222-200-00                | 0222-002-01                    |

<sup>1</sup> Laser Drive is Continuous Wave (CW) with constant power.

<sup>2</sup> Power Supply: For modules with 5-10VDC it is recommended that the module be operated at 5V. For modules with 3-6VDC it is recommended that the module be operated at 3.3V or 5V.

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers


Accessories

# Lasiris

## Application Primers

### What is Structured Light?

Structured light is the projection of a light pattern (plane, grid, or more complex shape) at a known angle onto an object. This technique can be very useful for imaging and acquiring dimensional information. The most often used light pattern is generated by fanning out a light beam into a sheet-of-light. When a sheet-of-light intersects with an object, a bright line of light can be seen on the surface of the object. By viewing this line of light from an angle, the observed distortions in the line can be translated into height variations.

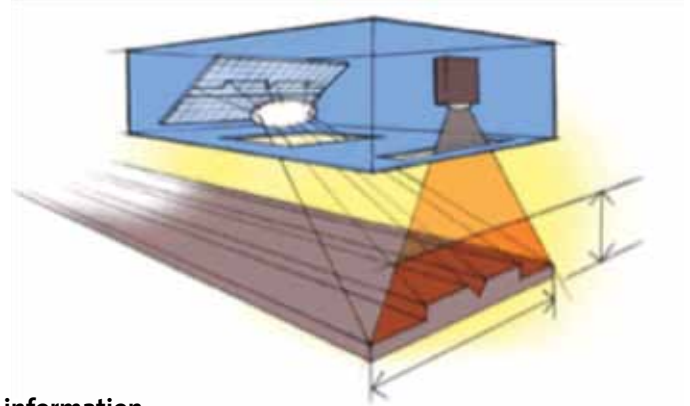

### Acquiring 2-dimensional information

Scanning the object with the light constructs 3-D information about the shape of the object. This is the basic principle behind depth perception for machines, or 3D machine vision. In this case, structured lighting is sometimes described as active triangulation.

Since structured lighting can be used to determine the shape of an object in machine vision applications, it can also help recognize and locate an object in an environment. These features make structured lighting useful in assembly lines implementing process control or quality control. These methods use structured lighting for alignment and inspection purposes. Structured light systems can

drastically transform a manufacturing plant by decreasing process variation, reducing production time, allowing automation of assembly lines, increasing precision, and generally decreasing overall cost. Although other types of light can be used for structured lighting, laser light is the best choice when precision and reliability are important.

### Structured Light Applications

Coherent's Lasiris uniform intensity laser projectors are especially useful for structured light applications, including Machine Vision, Inspection, and Alignment.

Machine vision is a combination of structured lighting, a detector, and a computer to precisely gather and analyze data. For example, it is used on robots as 3-D guiding systems to place or insert a part on a car, such as a windshield wiper or a door.

Structured light lasers used in inspection minimize process variation by drawing attention to parts that do not conform to specifications. They can pick out vegetables with blemishes from food processor lines or can ensure that the right colored capsule goes in the correct bottle in drug packaging lines.



Another laser application is alignment. In computer assembly, a laser system can help an operator determine if a computer chip is perfectly positioned on a circuit board.

Lasiris lasers can be useful for contour mapping of parts, surface defect detection, depth measurements, guidelines, edge detection, and alignment. Coherent has standard laser configurations available off-the-shelf for a wide variety of applications, but lasers can also be custom manufactured for OEM clients, designed for specific applications.

Companies that are interested in integrating laser technology for machine vision and industrial inspection into the manufacturing process should involve our engineers at the design stage.


# Lasiris

## Application Primers

### The Non-Gaussian Distribution

Most laser line generators on the market today use cylindrical optics to generate a line. We use a special Lasiris lens system that makes the lines non-gaussian (uniform) in intensity, even when the laser is off-axis.

The Gaussian or Non-Gaussian distribution refers to the distribution of power over the projected laser line. The phrase Gaussian distribution (also called Normal Distribution) is a statistical term that refers to a bell-shaped graph.

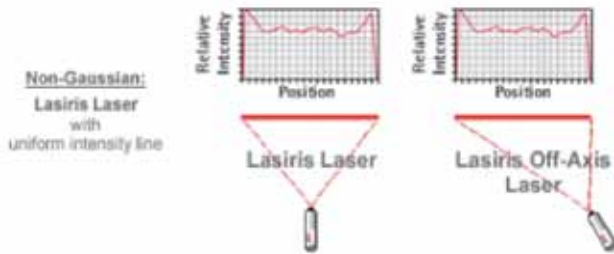


### Why Uniform Intensity Lines are Important


#### Non-Gaussian Lines are Efficient

The light intensity of a Gaussian line fades away towards the ends of the line, eventually falling below the threshold level of the detector and becoming invisible

to the system. Depending on the settings of the detector and the level of uniformity required by the application, as much as 50% of the available power can be lost.

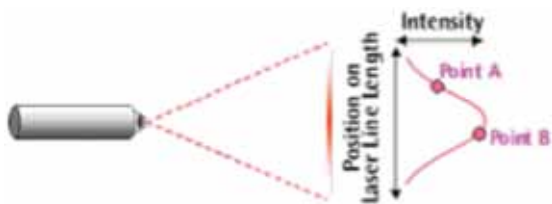


**The Competition** - The power is lowest at the tips of the projected laser line, and highest in the middle. Note how Points A and B have different intensities.

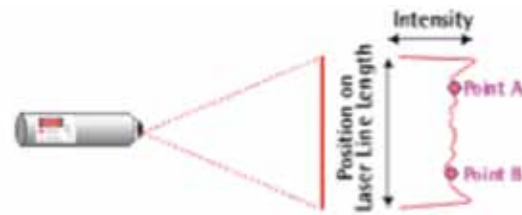

**Lasiris™ Laser** - A Lasiris line of light has a more uniform distribution of intensity. See how Points A and B have similar intensities.

#### Non-Gaussian Lines are Easy to Calibrate

Because the light intensity of gaussian lines is non-uniform, the calibration of CCDs can become very difficult. Separate calibrations must be made for pixels in the bright central area and for those in the transition area. The low intensity area cannot contribute to the calibration because it is invisible to the system.

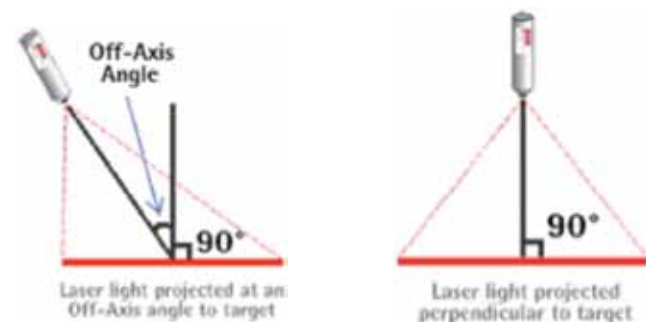

**An Off-Axis laser is a laser that is not perpendicular to the surface on which it is projecting laser light.**

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories

# Lasiris

## Application Primers

Lasiris Lasers can easily be adapted to an Off-Axis angle while maintaining uniform line intensity simply slide the line generating optic head slightly off center, see

figure below, turn the laser ON and move optical head up and down until the projected line is uniform, then tighten the set-screw.



### Diffractive Optics

The following equation is used to calculate the diffraction angle of the pattern (lines or dots).

$$d \sin(\theta_m) = m\lambda \text{ or } \theta_m = \sin^{-1}\left(\frac{m\lambda}{d}\right)$$

where  $d$  = DOE period ( $\mu\text{m}$ )

$\theta$  = diffraction angle of the  $m^{\text{th}}$  beam ( $^\circ$ )

$\lambda$  = wavelength ( $\mu\text{m}$ )

Single Line



Crosshair



Dot Line



7x7 Dot Matrix



Parallel Lines



4x4 Grid



Single Circle




7 Concentric Circles




$\theta_m$  is the diffraction angle measured from the center to the  $m^{\text{th}}$  beam. The interbeam angle is the angle between two neighboring beams (or orders) starting from the central beam (order 0) with the positive orders on one side (+1, +2, +3...) and the negative orders on the other side (-1, -2, -3...). The spreading angle,  $\theta_s$ , between the two outmost beams (or lines, when a line generating lens is added) is twice the angle of the higher order beam.

For example, for a 5-line pattern at 670 nm with a period  $d = 25 \mu\text{m}$ , the spreading angle  $\theta_s$  is twice the diffraction angle of the second order beam.

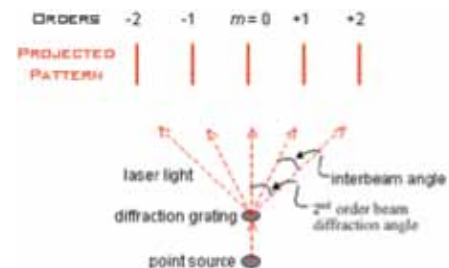

$$\theta_5 = 2\theta_m = 2\theta_2 \sin^{-1}\left(\frac{2 \cdot 670 \mu\text{m}}{25 \mu\text{m}}\right) = \theta_5 = 6.14^\circ$$

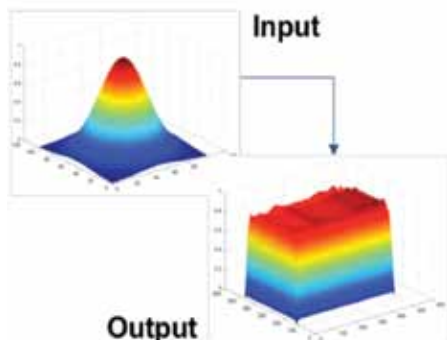

### Beam Shaping Optics (Flat Top)


Optical system that convert a Gaussian laser beam intensity profile into a "box" profile which is based on two divergent aspheric cylinders that purposely aberrates the beam to create a uniform far field irradiance pattern. Intensity from the center of the Gaussian is redistributed to its edges extremely efficiently. 1 or both axis can be shaped. This "Flat Top" is a more efficient use of optical power. Sharp edges are often required as this results in no residual light outside the region of interest.

# Lasiris

## Application Primers

Spatial uniformity directly correlates to system performance for Light Induce Fluorescence and scatter type


applications. This improves signal to noise especially for submicron LIF type applications.

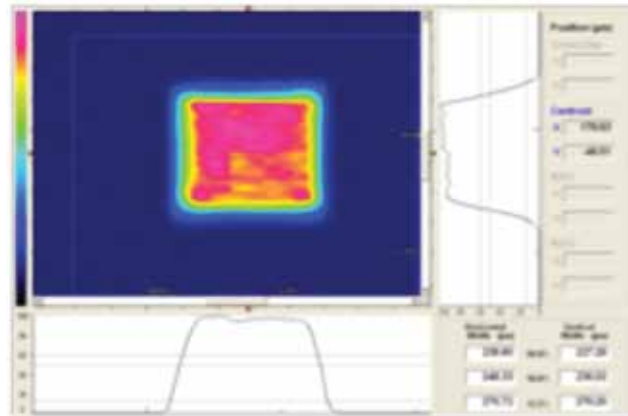

### Easier system alignment and more tolerant to system drift

The Flat top can be offered as either stand alone diverging flat top where the end user expands or focus the beam or as a complete focused system.

FA = Fan angle and is the divergence of the Powell, angles range from 0.25 up to 120 degrees  
E = Beam expander expansion power  
f = Focal length  
WD = Working distance  
 $\Phi_i$  = input beam diameter @ 1/e<sup>2</sup>

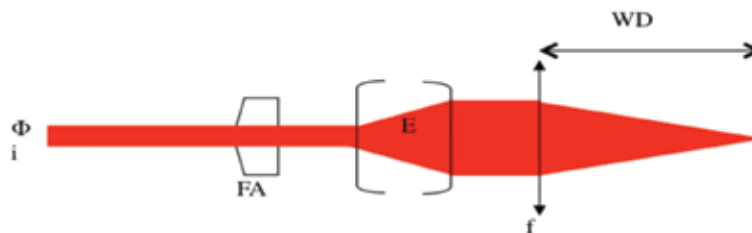

Dimension of the flat top is directly proportional to the fan angle of the Powell and the effective focal length of the focusing system, irrespective of distance between the Powell and the focussing system.

$$FTS \text{ (flat top size)} = \frac{2 \tan\left(\frac{FA}{2}\right) f}{E}$$

Flat top uniformity is defined by Coefficient of variation, CV, which is the standard deviation divided by the mean. We typically offer CV's of under 1 percent

$$CV = \frac{\text{StdDev}}{\text{Mean}}$$


The Depth of Uniformity is proportional to the size of the flat top and the effective focal length of the system and indirectly proportional to the input beam diameter

$$DOU = \frac{FTS * f}{2 \Phi_i}$$

$\Phi_i$  = distance where variation change is ~25%

The steepness ration is defined by the flat top size divided by the diffraction limited spot size without the shaper.

$$R \text{ (steepness ratio)} = \frac{\sin\left(\frac{FA}{2}\right) H \Phi}{2 \lambda}$$

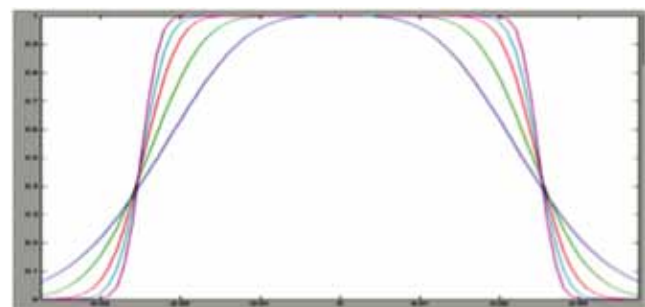

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories

# Lasiris SNF

## Structured Light Pattern Generating Laser


Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers




### Standard SNF Laser

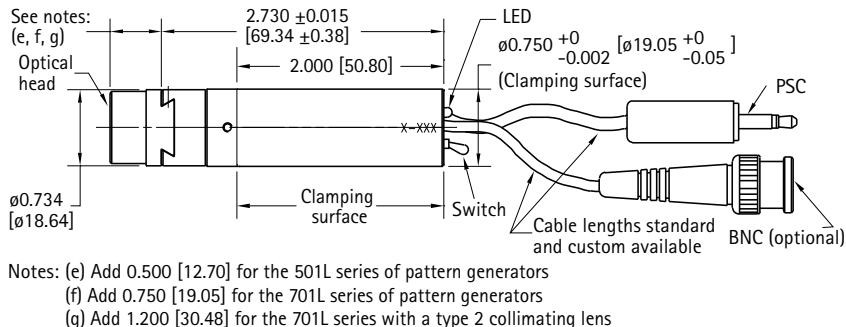

Top View

### SNF Laser with Separate Electronics Option

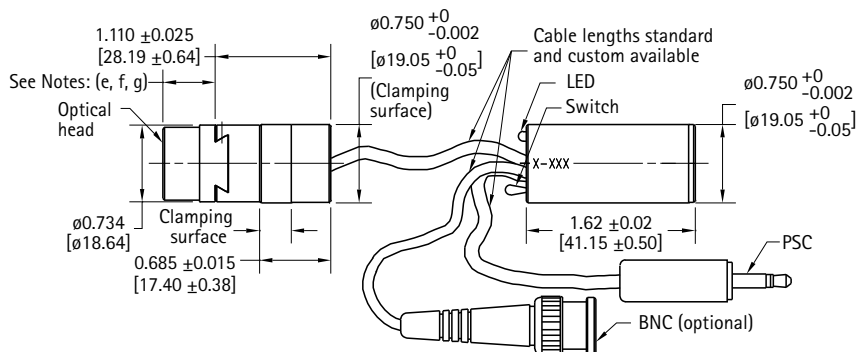

Top View

# Lasiris SNF

## Structured Light Pattern Generating Laser

### Key Features

- 635 nm to 1550 nm
- Power up to 200 mW
- Uniformity down to  $\pm 15\%$  for laser line generators
- Wide range of patterns and fan angles
- Focusable
- ESD, over-temperature, over-voltage, and reverse-polarity protection
- Rugged, shock and vibration resistant design
- 1D and 2D flat top compatible

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

### Ordering Information

SNF lasers are covered under a 2-year warranty (parts and labor). To order, use the following code: **MODEL-(Fast/Slow axis)-(Pattern)-(Interbeam Angle)-(Wavelength)(Electronic Option)-Power-Fan Angle-Custom Options**. E.g., **SNF-501L-635S-10-10**. Additional options can be specified at time of order. Contact us for more details.

|                              | 501L or 701L | 503L or 703L      | 505L or 705L | 509L or 709L | 511L or 711L |
|------------------------------|--------------|-------------------|--------------|--------------|--------------|
| Pattern                      | 1 line       | 3 lines           | 5 lines      | 9 lines      | 11 lines     |
| Interbeam Angle <sup>1</sup> | -            | 1.5°, 5.0°, 11.7° | 0.23°, 1.55° | 0.11°, 0.07° | 1.5°         |

|                              | 515L or 715L | 519L or 719L | 533L or 733L | 565L or 765L | 599L or 799L |
|------------------------------|--------------|--------------|--------------|--------------|--------------|
| Pattern                      | 15 lines     | 19 lines     | 33 lines     | 65 lines     | 99 lines     |
| Interbeam Angle <sup>1</sup> | 2.3°         | 0.77°        | 0.09°, 0.38° | 0.41°        | 0.149°       |

|                              | 501S     | 504G     | 501H      | 501C         | 507C                 |
|------------------------------|----------|----------|-----------|--------------|----------------------|
| Pattern                      | 1 square | 4x4 grid | crosshair | 1 circle     | 7 concentric circles |
| Interbeam Angle <sup>1</sup> | 2.9°     | 2.44°    | -         | 0.77°, 11.4° | 0.77°                |

|                              | 507X           | 519X             | Custom |
|------------------------------|----------------|------------------|--------|
| Pattern                      | 7x7 dot matrix | 19x19 dot matrix |        |
| Interbeam Angle <sup>1</sup> | 1.9°           | 0.77°            |        |

|                  |                                                                 |                      |        |                |          |              |
|------------------|-----------------------------------------------------------------|----------------------|--------|----------------|----------|--------------|
| Wavelength (nm)  | 635                                                             | 660                  | 690    | 785            | 810      | 830          |
| Diode Power (mW) | 1, 5, 10, 15, 35                                                | 1, 5, 10, 20, 35, 50 | 20, 35 | 20, 35, 75, 90 | 100, 200 | 30, 100, 150 |
| Fan Angle        | 1° <sup>2</sup> , 5°, 10°, 15°, 20°, 30°, 45°, 60°, 75°, custom |                      |        |                |          |              |

<sup>1</sup> The interbeam angle is specified at 670 nm. The angle is proportional to wavelength

<sup>2</sup> Not standard for crosshair projectors.


# Lasiris SNF

## Structured Light Pattern Generating Laser

### Uniform Intensity

Conventional laser line patterns are often generated by cylindrical optics that produce a Gaussian line profile with a bright center and fading ends. Lasiris patented beam shaping optics spread the light into an evenly illuminated line. The result is a crisp, uniform line with sharp ends.

### Line Intensity Profile Along Line Length

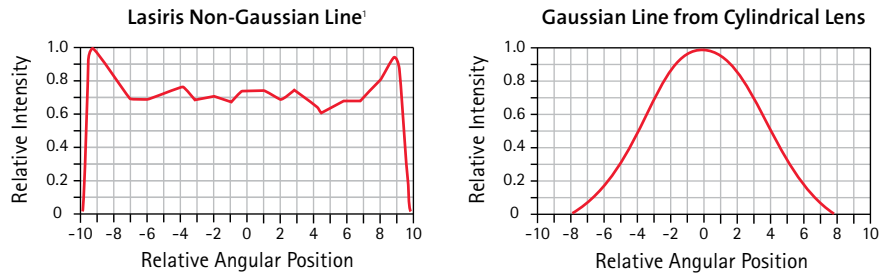


Relative Intensity vs. Angular Position Along Line Length

### Focusing Performance

The following figures show the typical focusing and depth-of-field performance of the SNF laser line generator. The focus charts indicate the minimum line thickness (at  $1/e^2$ ) achievable for a specific projection distance. The depth-of-field is defined as twice the distance over which the thickness of the line has increased by a factor of  $\sqrt{2}$ .

### Focusing and Depth-of-Field Performance

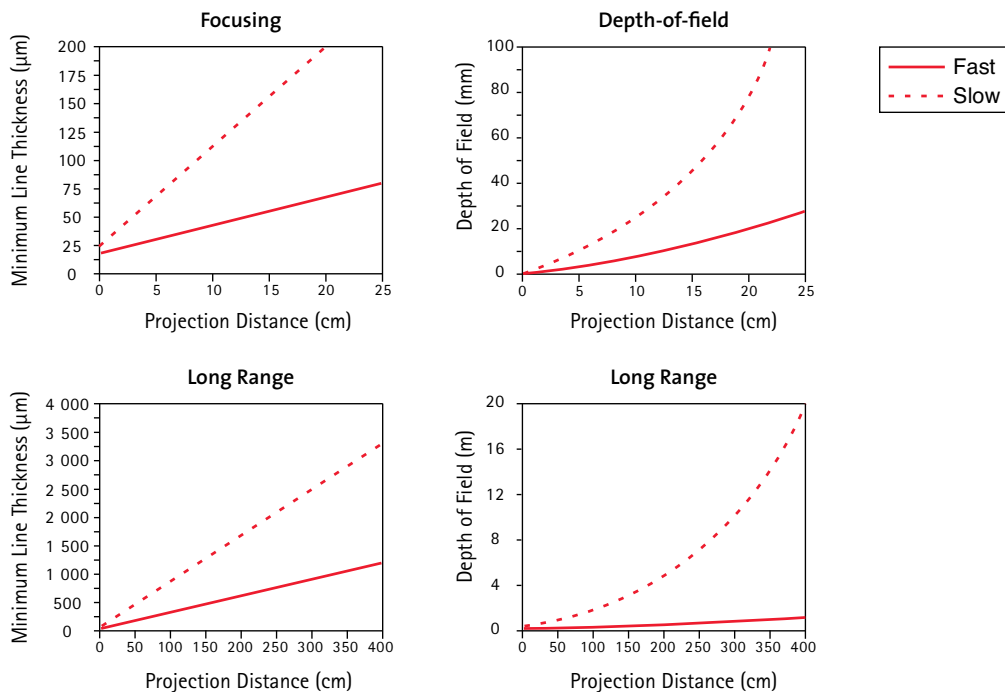


<sup>1</sup> Typical profile.

# Lasiris SNF

## Structured Light Pattern Generating Laser

|                                     |                                                      |                                                           |
|-------------------------------------|------------------------------------------------------|-----------------------------------------------------------|
| <b>Optical Specifications</b>       | Wavelength (nm)                                      | 635 to 1550                                               |
|                                     | Output Power (mW)                                    | 1 to 200                                                  |
|                                     | Intensity Distribution <sup>1</sup>                  | Uniform (non-Gaussian) lengthwise, Gaussian widthwise     |
|                                     | Fan Angles <sup>2</sup> (°)                          | 1 to 120 and custom                                       |
|                                     | Line Uniformity <sup>2</sup>                         | 30% uniformity over 80% of the pattern <sup>1</sup>       |
|                                     | Bore Sighting <sup>4</sup> (mrad)                    | <3                                                        |
| <b>Mechanical Specifications</b>    | Weight (g)                                           | 65 (with separate electronic: 70)                         |
|                                     | Housing Material                                     | Black anodized aluminum                                   |
| <b>Environmental Specifications</b> | Operating Temperature (°C)                           | -10 to +48 with bracket                                   |
|                                     | Wavelength Drift (nm/°C)                             | 0.25 typical                                              |
|                                     | Storage Temperature (°C)                             | -20 to 60                                                 |
|                                     | Over-Temperature Protection (°C)                     | up to 48                                                  |
| <b>Electrical Specifications</b>    | Power Supply Voltage (VDC)                           | 4.8 to 6.5 (6 optimal)<br>Optional 9, 12, 24, 115/220 VAC |
|                                     | Connector Type                                       | Male phono-jack 3.5 mm Ø or custom                        |
|                                     | Slow Start Time Delay (msec)                         | <250 for most models                                      |
|                                     | Reverse-Polarity Protection, Over-Voltage Protection |                                                           |
|                                     | Beam Modulation                                      | Via optional BNC connector                                |
|                                     | ESD Protection (kV)                                  | Discharge: Contact = 4; Indirect = 4; Air = 8             |
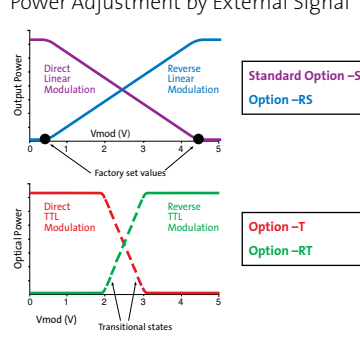
|                |                                                                                     |                                                                                                                                                                                                                                                                                                                                                                                                                                                                           |
|----------------|-------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>Options</b> | Power Options                                                                       |                                                                                                                                                                                                                                                                                                                                                                                                                                                                           |
|                | Power Adjustment Potentiometer                                                      | Laser power can be easily changed by adjusting an optional built-in potentiometer with a small screwdriver. Code "P".                                                                                                                                                                                                                                                                                                                                                     |
|                | Power Adjustment by External Signal                                                 |                                                                                                                                                                                                                                                                                                                                                                                                                                                                           |
|                |  | <p>Standard modulation<br/>S (synchro) or RS (reverse synchro): DC to 10 kHz<br/>S: Input voltage = 0VDC → laser "ON"<br/>RS: Input voltage = 5VDC → laser "ON"</p> <p>T (TTL) or RT (Reverse TTL): DC to 10 kHz<br/>Transition 2 to 3 VDC<br/>T: Input voltage = 0VDC → laser "ON"<br/>RT: Input voltage = 5V → laser "ON"</p> <p>"Fast" Modulation<br/>FS (fast synchro), FRS (fast reverse synchro) or<br/>FT (fast TTL), FRT (fast reverse TTL)<br/>DC to 300 kHz</p> |
|                | Input Impedance (kΩ)                                                                | >1                                                                                                                                                                                                                                                                                                                                                                                                                                                                        |
|                | Rise/Fall Time (μs)                                                                 | <10 for S and RS, <0.8 for FS and FRS                                                                                                                                                                                                                                                                                                                                                                                                                                     |
|                | Separate Electronic Option                                                          | The electronics of the laser can be separated from the main unit. See diagram "SNF laser with separate electronics option" on first page. Code "SD".                                                                                                                                                                                                                                                                                                                      |

|                           |                                                                                     |                                                                                     |                                                                                     |                                                                                     |                                                                                       |
|---------------------------|-------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------|
| <b>Available Patterns</b> | Single Line                                                                         | Crosshair <sup>3</sup>                                                              | Dot Line                                                                            | 7x7 Dot Matrix                                                                      | Single Dot                                                                            |
|                           |  |  |  |  |  |
|                           | Parallel Lines                                                                      | 4x4 Grid                                                                            | Single Circle                                                                       | 7 Concentric Circles                                                                |                                                                                       |
|                           |  |  |  |  |                                                                                       |

<sup>1</sup> Improve line uniformity to ±25, 20, or 15%. Request this option at the time of order.

<sup>2</sup> For line generators.

<sup>3</sup> Lasiris crosshair projectors have a single optical component, unlike conventional crosshairs that are formed either by using two lasers or by splitting and recombining one beam to form a cross.

<sup>4</sup> Define boresight method.

# Lasiris Mini

## Compact Structured Light Pattern Generating Laser

Diode Laser  
Systems


Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers




### Mechanical Specifications

#### Standard Mini Laser

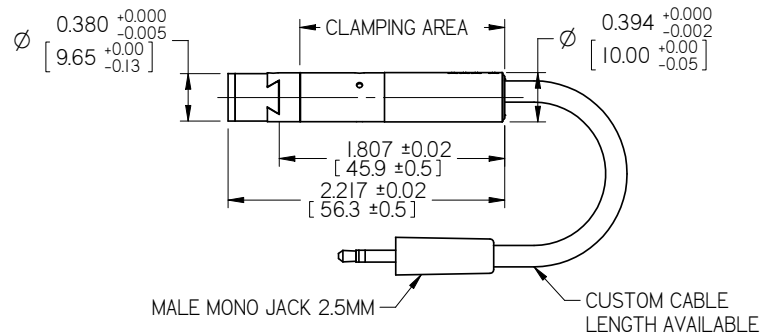

#### Mini Laser with Separate Electronics Option

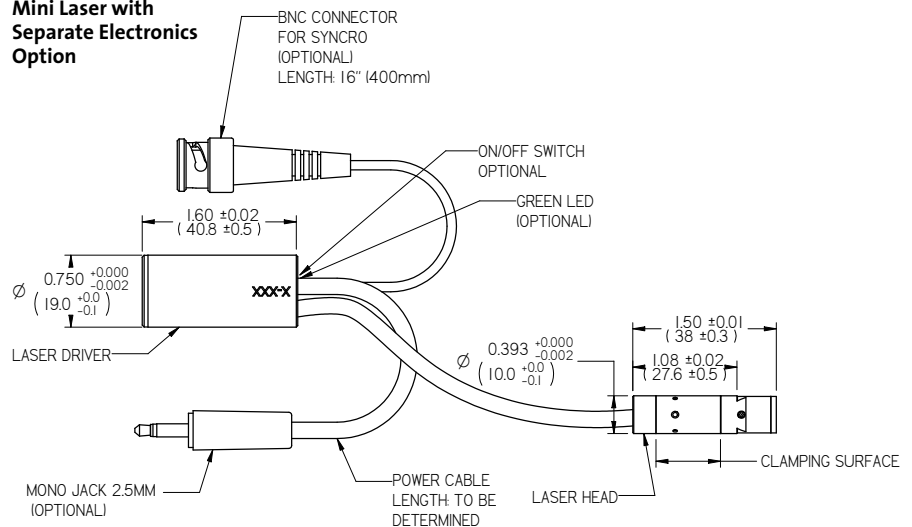

# Lasiris Mini

## Compact Structured Light Pattern Generating Laser

### Key Features

- Compact and lightweight
- 635 nm to 830 nm
- Wide range of powers, patterns, and fan angles
- Uniform intensity distribution with line generators
- Focusable
- ESD, over-temperature, and reverse-polarity protection
- 1D and 2D flat top compatible

### Ordering Information

Mini lasers are covered under a 2-year warranty (parts and labor). To order, use the following code: **MODEL-(Fast/Slow axis)-(Pattern)-(Interbeam Angle)-(Wavelength)(Electronic Option)-Power-Fan Angle-Custom Options**. E.g., **MINI-501L-635S-10-10**. Additional options can be specified at time of order. Contact us for more details.

|                  | 501L or 701L                                            | 503L or 703L      | 505L or 705L   | 509L or 709L     | 511L or 711L |        |
|------------------|---------------------------------------------------------|-------------------|----------------|------------------|--------------|--------|
| Pattern¹         | 1 line                                                  | 3 lines           | 5 lines        | 9 lines          | 11 lines     |        |
| Interbeam Angle² | -                                                       | 1.5°, 5.0°, 11.7° | 0.23°, 1.55°   | 0.11°, 0.07°     | 1.5°         |        |
|                  |                                                         |                   |                |                  |              |        |
|                  | 515L or 715L                                            | 519L or 719L      | 533L or 733L   | 599L or 799L     | 501S         |        |
| Pattern¹         | 15 lines                                                | 19 lines          | 33 lines       | 99 lines         | 1 square     |        |
| Interbeam Angle² | 2.3°                                                    | 0.77°             | 0.09°, 0.38°   | 0.149°           | 2.9°         |        |
|                  |                                                         |                   |                |                  |              |        |
|                  | 504G                                                    | 501H              | 507X           | 519X             | Custom       |        |
| Pattern¹         | 4x4 grid                                                | crosshair         | 7x7 dot matrix | 19x19 dot matrix |              |        |
| Interbeam Angle² | 2.44°                                                   | -                 | 1.9°           | 0.77°            |              |        |
|                  |                                                         |                   |                |                  |              |        |
| Wavelength (nm)  | 635                                                     | 660               | 690            | 785              | 830          | Custom |
| Diode Power (mW) | 1,5,10                                                  | 1,5,10,20,35      | 20,35          | 20,35            | 30           |        |
| Fan Angle        | 1°, 5°, 10°, 15°, 20°, 30°, 45°, 60°, 75°, 120°, custom |                   |                |                  |              |        |

<sup>1</sup> Line patterns are also available as dots -- add "D" nstead of "L" in the order code.

<sup>2</sup> At 660 nm.

<sup>3</sup> Not standard for crosshair projector.

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories


# Lasiris Mini

## Compact Structured Light Pattern Generating Laser

### Uniform Intensity

Conventional laser line patterns are often generated by cylindrical optics that produce a Gaussian line profile with a bright center and fading ends. Lasiris patented beam shaping optics spread the light into an evenly illuminated line. The result is a crisp, uniform line with sharp ends.

### Line Intensity Profile Along Line Length

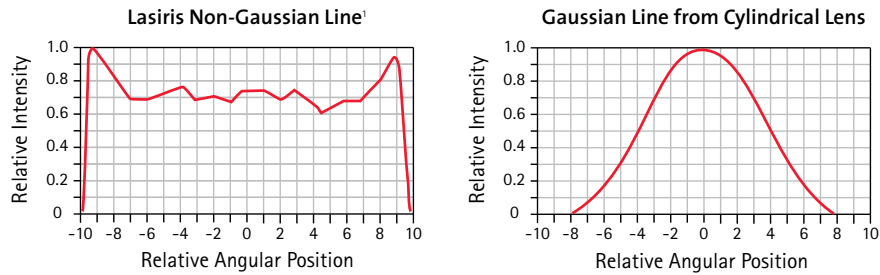


Relative Intensity vs. Angular Position Along Line Length

### Focusing Performance

The following figures show the typical focusing and depth-of-field performance of the SNF laser line generator. The focus charts indicate the minimum line thickness (at  $1/e^2$ ) achievable for a specific projection distance. The depth-of-field is defined as twice the distance over which the thickness of the line has increased by a factor of  $\sqrt{2}$ .

### Focusing and Depth-of-Field Performance

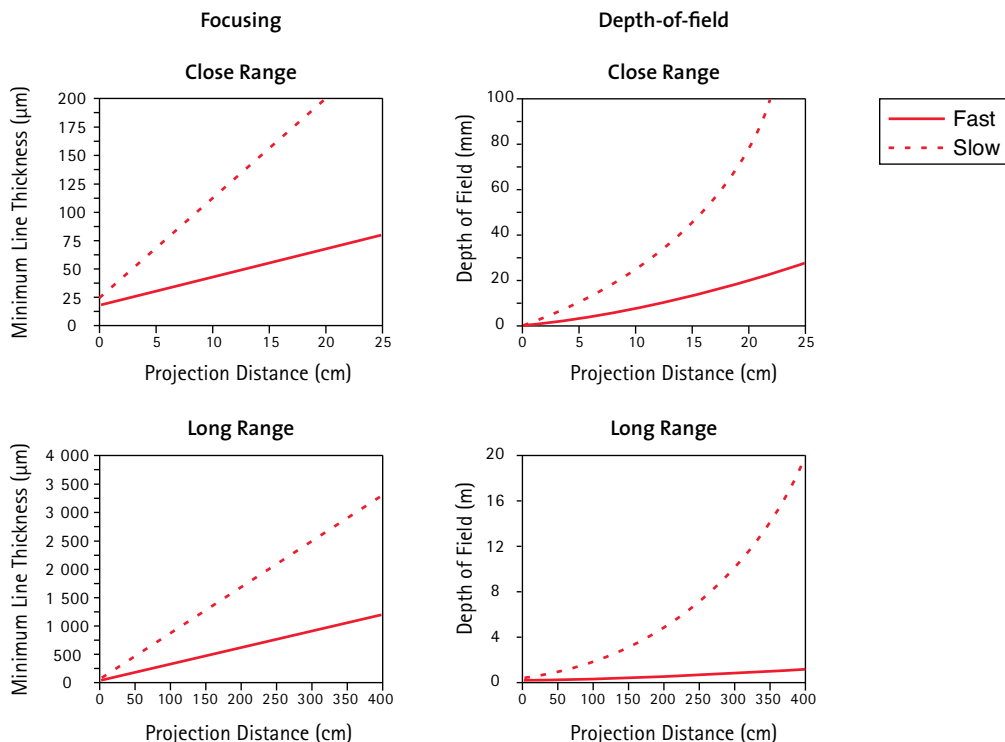

<sup>1</sup> Typical profile.

# Lasiris Mini

## Compact Structured Light Pattern Generating Laser

|                                     |                                  |                                                                  |
|-------------------------------------|----------------------------------|------------------------------------------------------------------|
| <b>Optical Specifications</b>       | Wavelength (nm)                  | 635 to 830                                                       |
|                                     | Output Power (mW)                | 1 to 40                                                          |
|                                     | Intensity Distribution           | Uniform (non-Gaussian) lengthwise, Gaussian widthwise            |
|                                     | Fan Angles (°)                   | 1 to 120                                                         |
|                                     | Line Uniformity                  | ±30% uniformity over 80% of the pattern                          |
|                                     | Bore Sighting (mrad)             | <3                                                               |
| <b>Mechanical Specifications</b>    | Weight (g)                       | 12                                                               |
|                                     | Housing Material                 | Black anodized aluminum                                          |
| <b>Environmental Specifications</b> | Operating Temperature (°C)       | -10 to +48 with bracket                                          |
|                                     | Wavelength Drift (nm/°C)         | 0.25 typical                                                     |
|                                     | Storage Temperature (°C)         | -20 to 60                                                        |
|                                     | Over-Temperature Protection (°C) | up to 48                                                         |
| <b>Electrical Specifications</b>    | Power Supply Voltage (VDC)       | 5 ±5%<br>6 with optional separate driver<br>Optional 115/220 VAC |
|                                     | Connector Type                   | Male phono-jack 2.5 mm Ø or custom                               |
|                                     | Slow Start Time Delay (msec)     | 250                                                              |
|                                     | Reverse-Polarity Protection      |                                                                  |
|                                     | ESD Protection (kV)              | Discharge: Contact = 4; Indirect =4; Air =8                      |


### Options

#### Power Options

Power Adjustment Potentiometer

Laser power can be easily changed by adjusting an optional built-in potentiometer with a small screwdriver. Code "P".

Power Adjustment by External Signal

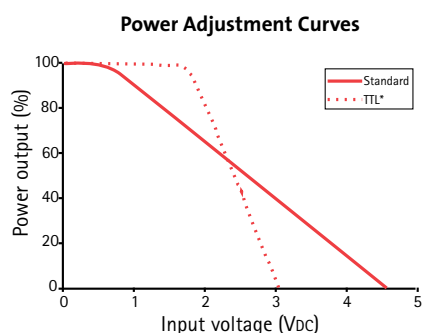

Standard Modulation (Available only on SD configuration)

S (syncro) or RS (reverse syncro): DC to 10 kHz

Linear for amplitude 1.0 to 4.0 Vdc

S: Input Voltage = 0VDC → laser "ON"

RS: Input Voltage = 5VDC → laser "ON"

T (TTL): DC to 10 kHz

Transition 2 to 3 VDC

T: Input Voltage = 0VDC → laser "ON"

"Fast" Modulation

FS (fast syncro), FRS (fast reverse syncro) or FT (fast TTL)

DC to 300 kHz

Input Impedance (kΩ)

>1

Rise/Fall Time (μs)

<10

Separate Electronic Option

The electronics of the laser can be separated from the main unit. See diagram "Mini laser with separate electronics option" on first page. Code "SD".

### Available Patterns

Single Line



Crosshair<sup>1</sup>



Dot Line



7x7 Dot Matrix



Parallel Lines



Single Square



Single Dot



<sup>1</sup> Lasiris crosshair projectors have a single optical component, unlike conventional crosshairs that are formed either by using two lasers or by splitting and recombining one beam to form a cross.

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

# Lasiris DLS and DLSC

## Elliptical and Circular Spot Projectors

Diode Laser  
Systems


Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers




### Mechanical Specifications

### Standard DLS/C Lasers

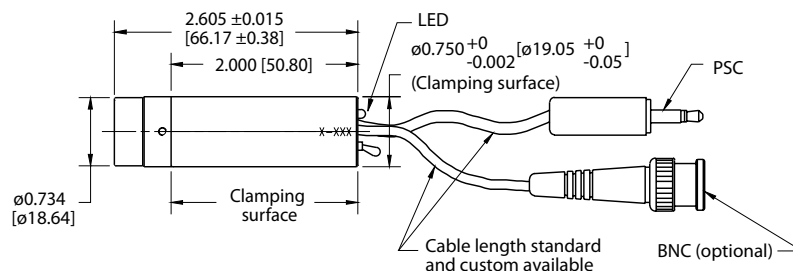

### DLS/C Laser with Separate Electronics Option

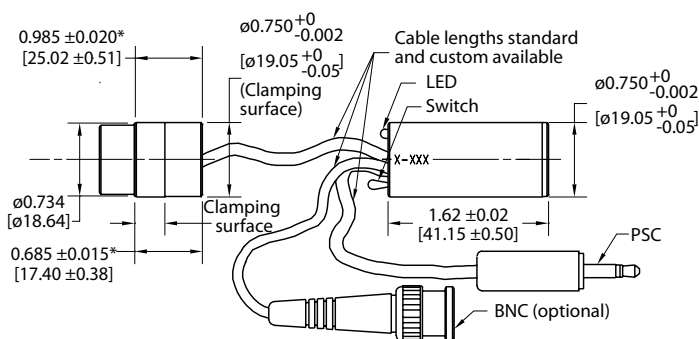

\*Add 0.325 in. (8.26 mm) for DLSC lasers.  
Some combinations of specifications may result in different dimensions.

# Lasiris DLS and DLSC

## Elliptical and Circular Spot Projectors

### Key Features

- Elliptical or circular beam
- Wide range of wavelengths and powers
- Focusable
- ESD, over-temperature, over-voltage, and reverse-polarity protection
- Rugged, shock and vibration resistant design

### Ordering Information

Lasers are covered under warranty (parts & labor) for 2 years (DLS) and 1 year (DLSC). To order, use the following code: **MODEL-(Fast/Slow)-Wavelength(Electronic Option)-Power-Custom Options**. E.g., **DLS-635T-1-SD**. Contact us for more details.

| Wavelength(nm)        | 635              | 660                       | 690    | 785        | 830          |
|-----------------------|------------------|---------------------------|--------|------------|--------------|
| DLS Diode Power (mW)  | 1, 5, 10, 15, 35 | 1, 5, 10, 20, 35, 50, 100 | 20, 35 | 20, 35, 75 | 30, 100, 150 |
| DLSC Diode Power (mW) | 1, 5, 10, 15, 35 | 1, 5, 20, 35, 50          | 20, 35 | 20, 35, 75 | 100, 150     |

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories

# Lasiris DLS and DLSC


## Elliptical and Circular Spot Projectors

### Focusing Performance

The DLSC typically deviates from circularity by no more than 1.2 to 1 ratio of the major to the minor axis. For the DLS, the spot is elliptical (1:3 ratio of axes) because the laser diode has a rectangular aperture. The following figures show the typical focusing and depth-of-field (at  $1/e^2$ ) performance of the DLS and DLSC spot projectors. The solid lines on the graphs apply to the major axis of the DLS and the dotted lines apply to the DLSC and to the minor axis of the DLS.

The focus charts indicate the minimum spot thickness (at  $1/e^2$ ) achievable for a specific projection distance. The depth-of-field is defined as twice the distance over which the spot thickness has increased by a factor of  $\sqrt{2}$ .


### Focusing and Depth-of-Field Performance

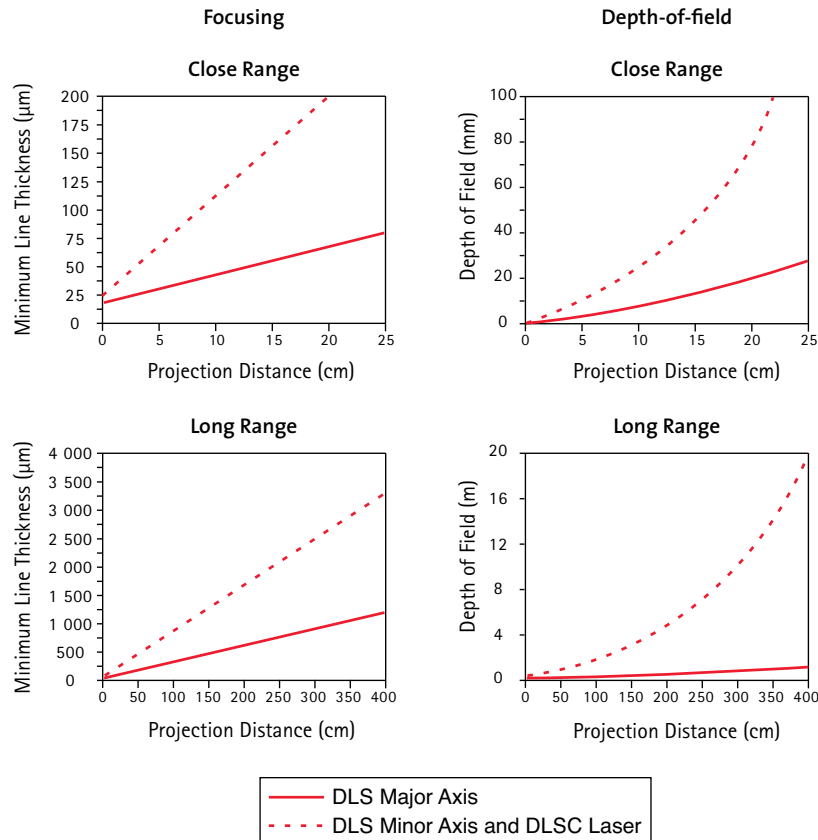

The smallest achievable beam width can be determined using these focus charts. However, the laser can also be adjusted to project a wider spot at a given projection distance, or collimated for minimum divergence. The focus is adjustable with our C-Thru Focus key, or the focus can be set at the factory. By specifying the desired spot diameter and working distance, the laser can be preset to your precise requirements. A collimating focal lens is available on request.

# Lasiris DLS and DLSC

## Elliptical and Circular Spot Projectors
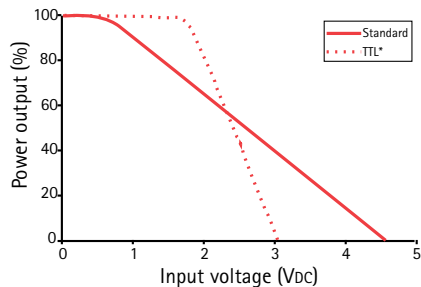
|                                     |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                   |                                                                                                                                                    |
|-------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>Optical Specifications</b>       | Wavelength (nm)                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                   | 635 to 1550 (varies with model)                                                                                                                    |
|                                     | Output Power (mW)                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                 | 1 to 150 (varies with model)                                                                                                                       |
|                                     | Spot Diameter                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     | See focus charts                                                                                                                                   |
|                                     | Bore Sighting (mrad)                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                              | <3                                                                                                                                                 |
| <b>Mechanical Specifications</b>    | Weight (g)                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                        | DLS: 54, DLSC: 60                                                                                                                                  |
|                                     | Housing Material                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                  | Black anodized aluminum                                                                                                                            |
| <b>Environmental Specifications</b> | Operating Temperature (°C)                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                        | -10 to +48 with bracket                                                                                                                            |
|                                     | Wavelength Drift (nm/°C)                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                          | 0.25 typical                                                                                                                                       |
|                                     | Over-Temperature Protection (°C)                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                  | Built-in: 48                                                                                                                                       |
| <b>Electrical Specifications</b>    | Power Supply Voltage (VDC)                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                        | 5 to 6<br>Optional 9, 12, 24, 115/220 VAC                                                                                                          |
|                                     | Connector Type                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                    | Male phono-jack 3.5 mm Ø or custom                                                                                                                 |
|                                     | Slow Start Time Delay (µsec)                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      | 250 (for most models)                                                                                                                              |
|                                     | Reverse-Polarity Protection, Over-Voltage Protection                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                              |                                                                                                                                                    |
| <b>Options</b>                      | Power Options                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     |                                                                                                                                                    |
|                                     | <p>Power Adjustment Potentiometer</p> <p>Power Adjustment by External Signal</p> <div data-bbox="332 995 750 1323"> <p><b>Power Adjustment Curves</b></p>  <p>The graph plots Power output (%) on the y-axis (0 to 100) against Input voltage (Vdc) on the x-axis (0 to 5). Two curves are shown: a solid red line for 'Standard' and a dotted red line for 'TTL*'. The 'Standard' curve starts at 100% at 0V, remains flat until ~0.5V, then decreases linearly to 0% at 5V. The 'TTL*' curve starts at 100% at 0V, remains flat until ~1.5V, then decreases linearly to 0% at 3V.</p> </div> |                                                                                                                                                    |
|                                     | <p>Laser power can be easily changed by adjusting an optional built-in potentiometer with a small screwdriver. Code "P".</p> <p>Standard modulation<br/>S (synchro) or RS (reverse synchro): DC to 10 kHz<br/>S: Input voltage = 0VDC → laser "ON"<br/>RS: Input voltage = 5VDC → laser "ON"</p> <p>T (TTL) or RT (Reverse TTL): DC to 10 kHz<br/>Transition 2 to 3 VDC<br/>T: Input voltage = 0VDC → laser "ON"<br/>RT: Input voltage = 5V → laser "ON"</p> <p>"Fast" Modulation<br/>FS (fast synchro), FRS (fast reverse synchro) or<br/>FT (fast TTL), FRT (fast reverse TTL)<br/>DC to 300 kHz</p>                                                                            |                                                                                                                                                    |
|                                     | Input Impedance (kΩ)                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                              | >1                                                                                                                                                 |
|                                     | Rise/Fall Time (µs)                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                               | 10 for 10 KHz, 0.8 for 300 kHz                                                                                                                     |
|                                     | Separate Electronic Option                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                        | Electronics of the laser can be separated from the main unit. See diagram "DLS/C laser with separate electronics option" in first page. Code "SD". |

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

# Lasiris GLL

Green Structured Light Pattern Generating Laser


Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers



## Mechanical Specifications

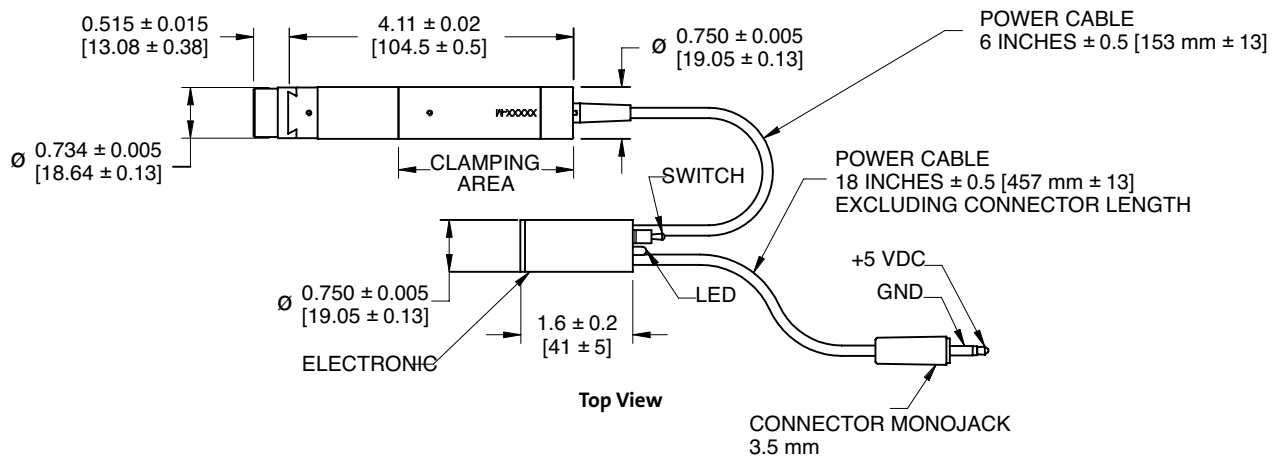

# Lasiris GLL

## Green Structured Light Pattern Generating Laser

### Key Features

- High visibility, high contrast laser
- Stable output power
- Line generators that have a uniform intensity distribution
- Many patterns available
- Focusable
- ESD, over-voltage, and reverse-polarity protection
- Uniformity down to  $\pm 15\%$  for line generators

### Ordering Information

Green lasers are covered under a warranty of 1 year (parts and labor). To order, use this code: **MODEL-(Fast/Slow axis) (Pattern)-(Interbeam Angle)-("532" for Wavelength)(Electronic Option)-Power-Fan Angle-Custom Options**. E.g., **GLL-503L(1.2°)-532-1-5-SD**.

|                              | 501L or 501D                                                    | 503L or 503D      | 505L or 505D      | 509L or 509D      | 511L or 511D        |
|------------------------------|-----------------------------------------------------------------|-------------------|-------------------|-------------------|---------------------|
| Pattern                      | 1 line or 1 dot                                                 | 3 lines or 3 dots | 5 lines or 5 dots | 9 lines or 9 dots | 11 lines or 11 dots |
| Interbeam Angle <sup>1</sup> | -                                                               | 1.2°, 4.0°, 9.2°  | 0.18°, 1.23°      | 0.09°, 0.06°      | 1.2°                |
| Diode Power (mW)             | 1, 5, 10                                                        |                   |                   |                   |                     |
| Fan Angle                    | 1° <sup>2</sup> , 5°, 10°, 15°, 20°, 30°, 45°, 60°, 75°, custom |                   |                   |                   |                     |

|                              | 515L or 5151D                                                   | 519L or 519D        | 533L or 533D        | 599L or 599D        | 501S     |
|------------------------------|-----------------------------------------------------------------|---------------------|---------------------|---------------------|----------|
| Pattern                      | 15 line or 15 dots                                              | 19 lines or 19 dots | 33 lines or 33 dots | 99 lines or 99 dots | 1 square |
| Interbeam Angle <sup>1</sup> | 1.8°                                                            | 0.61°               | 0.07°, 0.30°        | 0.118°              | 2.3°     |
| Diode Power (mW)             | 1, 5, 10                                                        |                     |                     |                     |          |
| Fan Angle                    | 1° <sup>2</sup> , 5°, 10°, 15°, 20°, 30°, 45°, 60°, 75°, custom |                     |                     |                     |          |

|                              | 504G                                                            | 501H      | 501C        | 507C                 | 507X           |
|------------------------------|-----------------------------------------------------------------|-----------|-------------|----------------------|----------------|
| Pattern                      | 4x4 grid                                                        | crosshair | 1 circle    | 7 concentric circles | 7x7 dot matrix |
| Interbeam Angle <sup>1</sup> | 1.93°                                                           | -         | 0.61°, 9.0° | 0.61°                | 1.5°           |
| Diode Power (mW)             | 1, 5, 10                                                        |           |             |                      |                |
| Fan Angle                    | 1° <sup>2</sup> , 5°, 10°, 15°, 20°, 30°, 45°, 60°, 75°, custom |           |             |                      |                |

|                              | 501H                                                            | Custom |
|------------------------------|-----------------------------------------------------------------|--------|
| Pattern                      | 19x19 dot matrix                                                |        |
| Interbeam Angle <sup>1</sup> | 0.61°                                                           |        |
| Diode Power (mW)             | 1, 5, 10                                                        |        |
| Fan Angle                    | 1° <sup>2</sup> , 5°, 10°, 15°, 20°, 30°, 45°, 60°, 75°, custom |        |

<sup>1</sup> At 532 nm.

<sup>2</sup> Not standard for crosshair projectors.

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories


# Lasiris GLL

## Green Structured Light Pattern Generating Laser

### Uniform Intensity

Laser line patterns are often generated by cylindrical optics that produce a Gaussian line profile with a bright center and fading ends. Lasiris patented optics spread the light into an evenly illuminated line with sharp ends.

### Line Intensity Profile Along Line Length

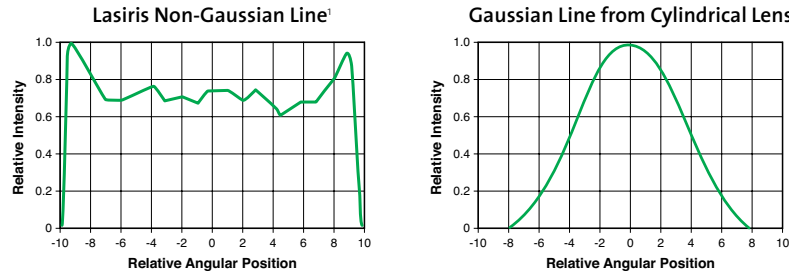


Relative Intensity vs. Angular Position Along Line Length

### Focusing Performance

The following figures show the typical focusing and depth-of-field performance of the Green laser projector. The focus charts indicate the minimum line thickness (at  $1/e^2$ ) achievable for a specific projection distance. The depth-of-field is defined as twice the distance over which the thickness of the line has increased by a factor of  $\sqrt{2}$ .

### Focusing and Depth-of-Field Performance

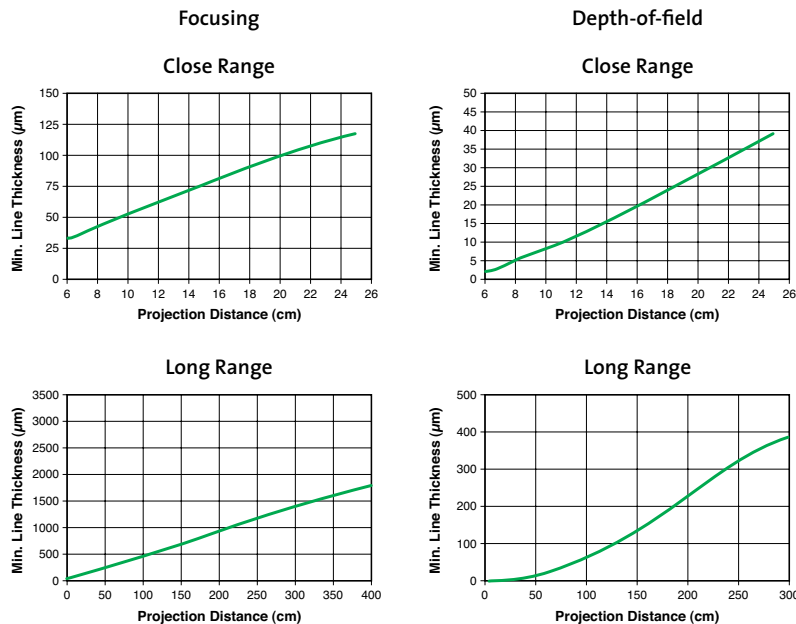

<sup>1</sup> Typical profile.

# Lasiris GLL

## Green Structured Light Pattern Generating Laser

|                                     |                                             |                                                       |
|-------------------------------------|---------------------------------------------|-------------------------------------------------------|
| <b>Optical Specifications</b>       | Wavelength (nm)                             | 532 ±1                                                |
|                                     | Output Power (mW)                           | 1, 5, 10                                              |
|                                     | Mode                                        | TEM <sub>00</sub> (Gaussian profile)                  |
|                                     | Intensity Distribution <sup>1</sup>         | Uniform (non-Gaussian) lengthwise, Gaussian widthwise |
|                                     | Fan Angles (°)                              | 1 to 120 and custom                                   |
|                                     | Line Uniformity                             | 30% uniformity over 80% of the pattern <sup>1</sup>   |
|                                     | M <sup>2</sup>                              | <1.5                                                  |
|                                     | Power Stability Over Temperature Range (%)  | 10                                                    |
|                                     | Power Stability at Constant Temperature (%) | 5                                                     |
|                                     | Factory Set Focus (m)                       | 1                                                     |
|                                     | Line Thickness (focus)                      | See focusing charts                                   |
| <b>Mechanical Specifications</b>    | Bore Sighting (mrad)                        | <3 over focus range                                   |
|                                     | Weight (g)                                  | 100                                                   |
| <b>Environmental Specifications</b> | Housing Material                            | Black anodized aluminum                               |
|                                     | Operating Temperature (°C)                  | +15 to +35                                            |
| <b>Electrical Specifications</b>    | Storage Temperature (°C)                    | -40 to +50                                            |
|                                     | Wavelength Drift (nm/°C)                    | 0.25 typical                                          |
|                                     | ESD Protection (kV)                         | Discharge: Contact = 4; Indirect = 4; Air = 8         |
|                                     | Power Supply Voltage (VDC)                  | 5 ±1                                                  |
| <b>Options</b>                      | Connector Type                              | Male phono-jack 3.5 mm                                |
|                                     | Operating Mode                              | Constant average power                                |
|                                     | Over-Voltage Protection (V)                 | Up to 12                                              |
|                                     | Reverse-Polarity Protection                 |                                                       |


**Options**

Power Options

Power Adjustment Potentiometer

Laser power can be easily changed by adjusting an optional built-in potentiometer with a small screwdriver. Code "P".

Power Adjustment by External Signal

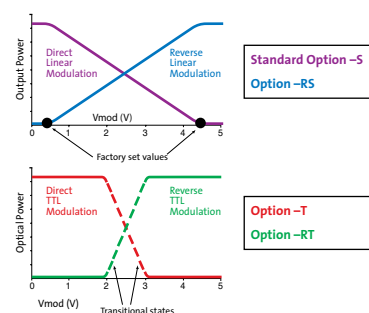


Standard modulation  
S (synchro) or RS (reverse synchro): DC to 10 kHz  
S: Input voltage = 0VDC → laser "ON"  
RS: Input voltage = 5VDC → laser "ON"

T (TTL) or RT (Reverse TTL): DC to 10 kHz  
Transition 2 to 3 VDC  
T: Input voltage = 0VDC → laser "ON"  
RT: Input voltage = 5V → laser "ON"

"Fast" Modulation  
FS (fast synchro), FRS (fast reverse synchro) or  
FT (fast TTL), FRT (fast reverse TTL)  
DC to 300 kHz

### Available Patterns<sup>2</sup>

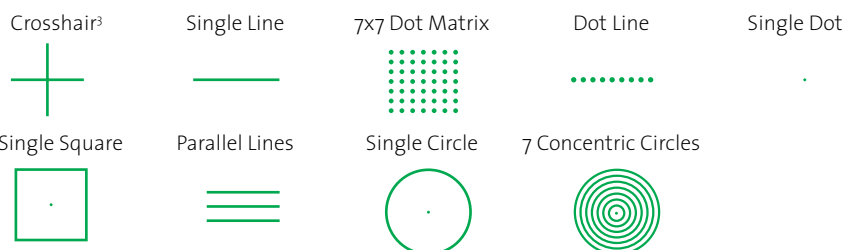

<sup>1</sup> Improve line uniformity to ±25, 20, or 15%. Request this option at the time of order.

<sup>2</sup> Other patterns available.

<sup>3</sup> Lasiris crosshair projectors have a single optical component, unlike conventional crosshairs that are formed either by using two lasers or by splitting and recombining one beam to form a cross.

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories



# Lasiris TMFL

## Micro-Focus Telecentric Laser Line Pattern Generating Laser

### Key Features

- Constant line width along entire line length
- Extremely thin lines down to 9.0  $\mu\text{m}$  at  $1/e^2$
- Line length up to 60 mm
- Highly efficient optical system
- Uniform intensity distribution
- ESD, over-temperature, over-voltage, and reverse-polarity protection

### Ordering Information

TMFL lasers 660 nm are covered under a 2-year warranty (parts & labor). To order, select from the specifications below. **Code: TMFL-Wavelength-Diode Power-Electronics Option-Working Distance-Line Length.** E.g., **TMFL-660-100-P-55-10, TMFL-660-100-FT-125-30.** Call us or visit our website for updates and other specifications.

|                       |                      |            |            |
|-----------------------|----------------------|------------|------------|
| Wavelength (nm)       | 660                  |            |            |
| Diode Power (mW)      | 100 <sup>1</sup>     |            |            |
| Electronics Option    | N/A, P, T, S, FS, FT |            |            |
| Working Distance (mm) | 55                   | 120        | 125        |
| Line Length (mm)      | 5, 7, 10             | 10, 20, 30 | 30, 45, 60 |

<sup>1</sup> Also available in 670 nm - 500 mW, 810 nm - 1W. Power can be adjusted down to 20 mW by potentiometer as a standard option.

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories


# Lasiris TMFL

## Micro-Focus Telecentric Laser Line Pattern Generating Laser

### Uniform Intensity

Conventional laser line patterns are often generated by cylindrical optics that produce a Gaussian line profile with a bright center and fading ends. Lasiris patented beam shaping optics spread the light into an evenly illuminated line. The result is a crisp, uniform line with sharp ends.

### Line Intensity Profile Along Line Length

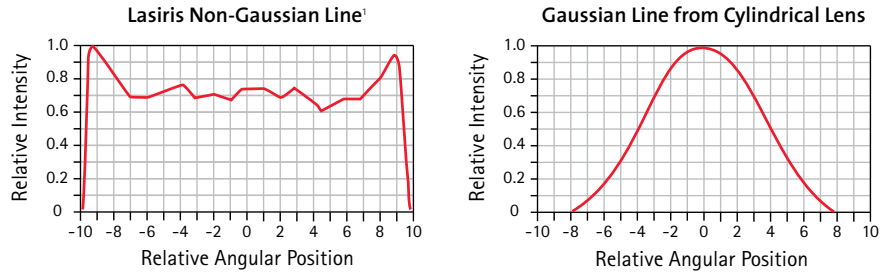

Relative Intensity vs. Angular Position Along Line Length


### Line Thickness and Depth-of-Field

TMFL lasers feature factory-set focusing for specific working distances. The table below shows the typical thickness (at  $1/e^2$ ) and depth-of-field performance of a 660 nm wavelength TMFL. The depth-of-field is defined as twice the distance over which the thickness of the line has increased by a factor of  $\sqrt{2}$ .

Main characteristic of the TMFL laser is its ability to maintain a constant line thickness along the entire line length when projected on a flat plane.

|                              | TMFL 55   | TMFL 120  | TMFL 125  |
|------------------------------|-----------|-----------|-----------|
| Depth-of-Field (mm)          | 0.2       | 1         | 1         |
| Working Distance (mm)        | 55        | 120       | 125       |
| Line Width ( $\mu\text{m}$ ) | 9         | 25        | 25        |
| Line Length² (mm)            | 5 to 10   | 10 to 30  | 30 to 60  |
| Line Uniformity (%)          | $\leq 25$ | $\leq 25$ | $\leq 25$ |

### Line Thickness vs. Position on Line

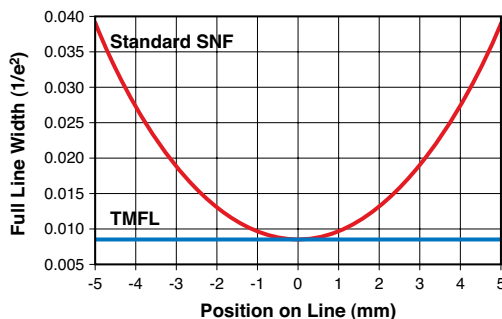



The "Thickness vs Position" graph shows how the line thickness of the SNF laser increases farther away from the center of the line. The TMFL is built with an optical system that creates a constant line thickness when the laser is projected on a flat plane.

<sup>1</sup> Typical profile.

<sup>2</sup> Line length can be customized.

# Lasiris TMFL

## Micro-Focus Telecentric Laser Line Pattern Generating Laser
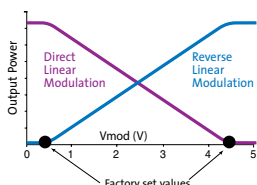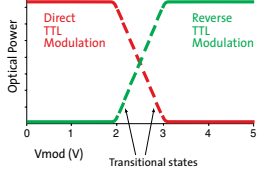
|                                     |                                                                                     |                                                                                                                                                                              |
|-------------------------------------|-------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>Optical Specifications</b>       | Wavelength (nm)                                                                     | 660                                                                                                                                                                          |
|                                     | Output Power (mW)                                                                   | 100                                                                                                                                                                          |
|                                     | Intensity Distribution                                                              | Uniform (non-Gaussian) for line                                                                                                                                              |
|                                     | Bore Sighting (mrad)                                                                | <3                                                                                                                                                                           |
| <b>Mechanical Specifications</b>    | Housing Material                                                                    | Black anodized aluminum                                                                                                                                                      |
| <b>Environmental Specifications</b> | Operating Temperature (°C)                                                          | -10 to +48                                                                                                                                                                   |
|                                     | Wavelength Drift (nm/°C)                                                            | 0.25 typical                                                                                                                                                                 |
|                                     | Thermal Stability                                                                   | Less than 5% change in focus over multiple temperature cycling from 10°C to 30°C                                                                                             |
|                                     | Over-Temperature Protection (°C)                                                    | Built-in: 48                                                                                                                                                                 |
| <b>Electrical Specifications</b>    | Power Supply Voltage (VDC)                                                          | 4.8 to 6.5 (6 optimal)<br>Optional 9, 12, 24, 115/220 VAC                                                                                                                    |
|                                     | Connector Type                                                                      | Male phono-jack 3.5 mm Ø or custom                                                                                                                                           |
|                                     | Slow Start Time Delay (msec)                                                        | <250 for most models                                                                                                                                                         |
|                                     | Reverse-Polarity Protection, Over-Voltage Protection                                |                                                                                                                                                                              |
|                                     | Beam Modulation                                                                     | Via optional BNC connector                                                                                                                                                   |
|                                     | ESD Protection (kV)                                                                 | Discharge: Contact = 4; Indirect = 4; Air = 8                                                                                                                                |
| <b>Options</b>                      | Power Options                                                                       |                                                                                                                                                                              |
|                                     | Power Adjustment Potentiometer                                                      | Laser power can be easily changed by adjusting an optional built-in potentiometer with a small screwdriver. Code "P".                                                        |
|                                     | Power Adjustment by External Signal                                                 |                                                                                                                                                                              |
|                                     |  | <p>Standard modulation</p> <p>S (synchro) or RS (reverse synchro): DC to 10 kHz</p> <p>S: Input voltage = 0VDC → laser "ON"</p> <p>RS: Input voltage = 5VDC → laser "ON"</p> |
|                                     |  | <p>T (TTL) or RT (Reverse TTL): DC to 10 kHz</p> <p>Transition 2 to 3 VDC</p> <p>T: Input voltage = 0VDC → laser "ON"</p> <p>RT: Input voltage = 5V → laser "ON"</p>         |
|                                     | Input Impedance (kΩ)                                                                | >1                                                                                                                                                                           |
|                                     | Rise/Fall Time (μs)                                                                 | 10 for 10 KHz, 0.8 for 300 KHz                                                                                                                                               |
|                                     | Separate Electronic Option                                                          | The electronics of the laser can be separated from the main unit. Code "SD".                                                                                                 |

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories


# Lasiris ColdRay

Compact Thermoelectrically Cooled Laser


Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

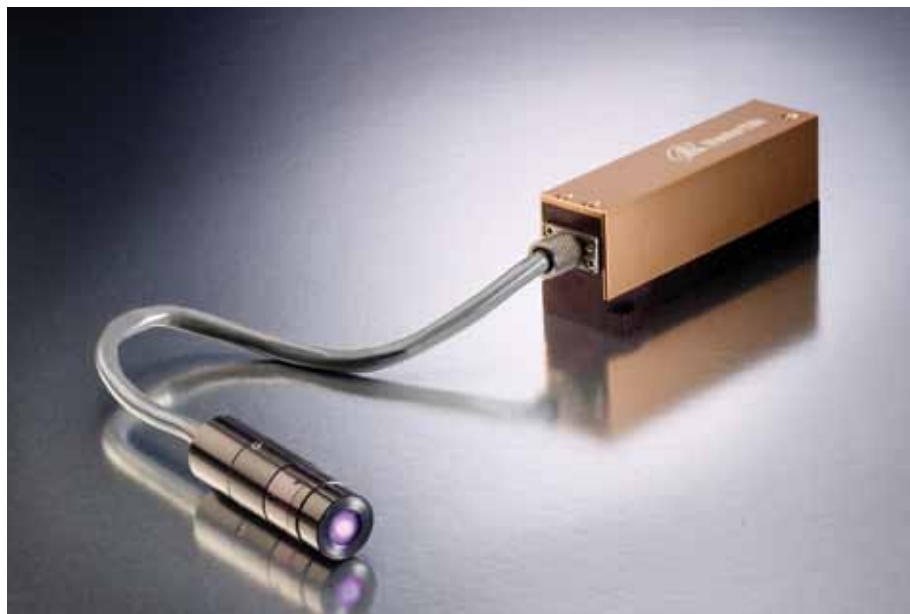

## Mechanical Specifications

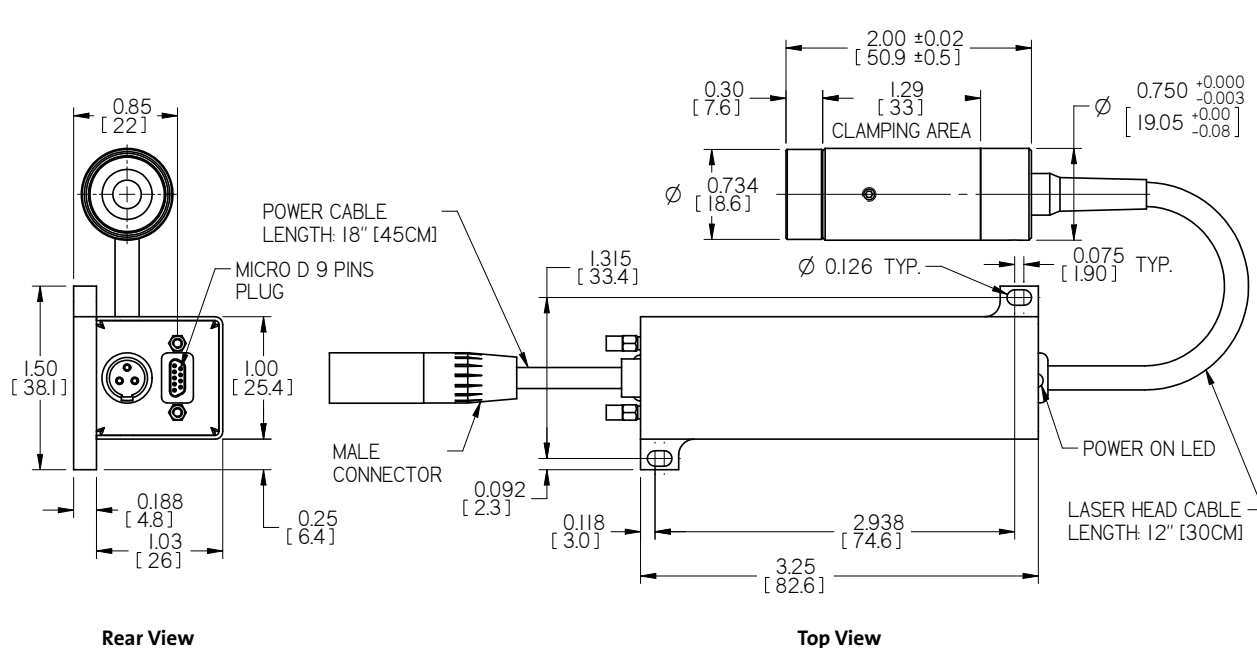

# Lasiris ColdRay

## Compact Thermoelectrically Cooled Laser

### Key Features

- 375 nm, 405 nm, 445 nm, 473 nm, 640 nm, 660 nm
- Excellent wavelength and power stability
- Excellent pointing stability
- Compatible with standard 1D and 2D flat top projector
- Extended lifetime
- Analog and digital modulation
- Circular models available
- Focusable

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

### Ordering Information

To order, use the following code: **MODEL-(Fast/Slow axis)(Pattern)-(Interbeam Angle)-(Wavelength)(Electronic Option)-Power-Fan Angle-Custom Options**. E.g., **CR-501L-405S-60-20**. Call us or visit our website for available optical patterns such as Flat-Top and diffractive optics.


|                    |              |                                                                               |     |     |               |                               |
|--------------------|--------------|-------------------------------------------------------------------------------|-----|-----|---------------|-------------------------------|
| Wavelength (nm)    | 375          | 405                                                                           | 445 | 473 | 640           | 660                           |
| Diode Power (mW)   | 20           | 5, 10, 20 <sup>1</sup> , 30, 40, 60 <sup>1</sup> , 80, 105 <sup>1</sup> , 120 | 50  | 20  | 5, 10, 35, 80 | 35, 50, 100, 150 <sup>1</sup> |
| Electronics Option | S, T, RS, RT |                                                                               |     |     |               |                               |

<sup>1</sup> No internal photodiode.


Accessories

### Compatible with Micro Flat Top Generators


Flat Top 2D

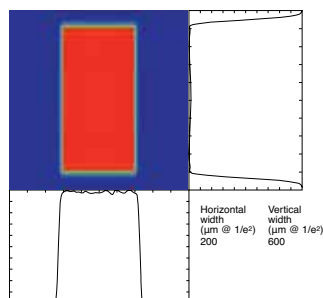


Flat Top Matrix 2D

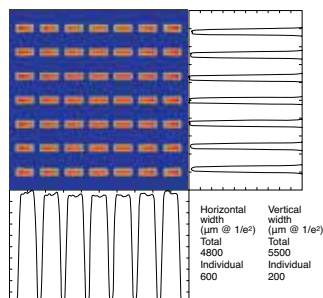

Flat Top 1D

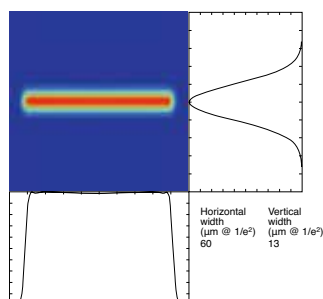

Flat Top Matrix 1D

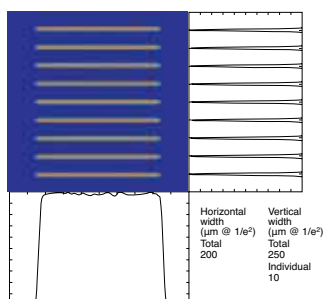


# Lasiris ColdRay

## Compact Thermoelectrically Cooled Laser

### Focusing Performance

The following figures show the typical focusing and depth-of-field performance of the ColdRay laser beam. The focus charts indicate the minimum line thickness (at  $1/e^2$ ) achievable for a specific projection distance. The depth-of-field is defined as twice the distance over which the thickness of the line has increased by a factor of  $\sqrt{2}$ .

### Focusing and Depth-of-Field Performance

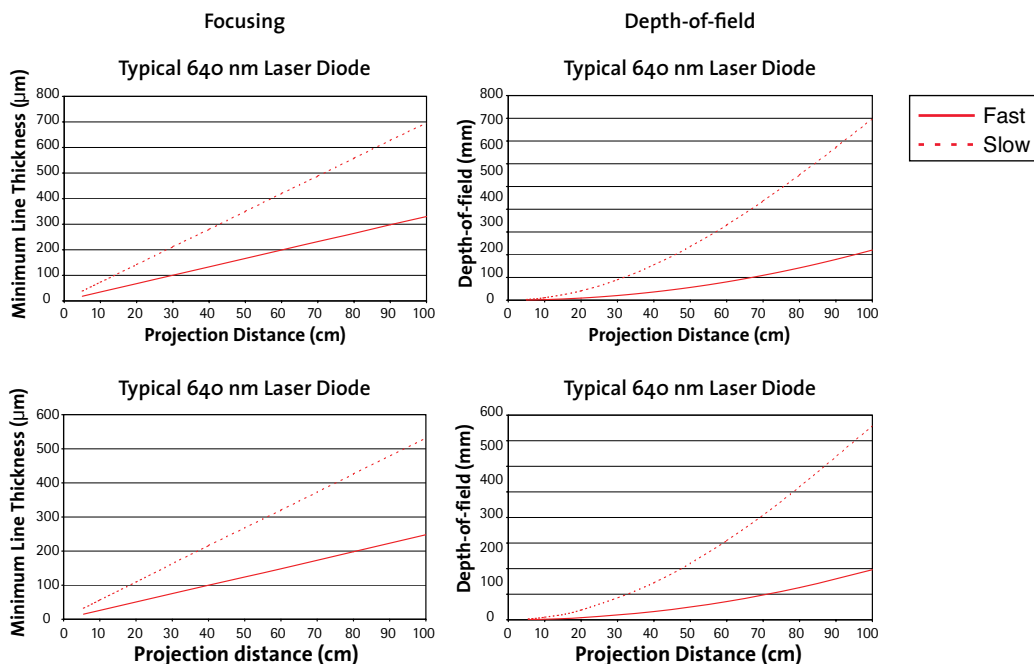

### Wavelength-Dependent Parameters

|                                                       |                   |                   |                   |                   |                   |                   |
|-------------------------------------------------------|-------------------|-------------------|-------------------|-------------------|-------------------|-------------------|
| Wavelength <sup>1</sup> (nm)                          | 375               | 405               | 445               | 473               | 640               | 660               |
| Wavelength Tolerance (nm)                             | +5/-5             | +8/-5             | +5/-5             | +5/-5             | +5/-5             | +10/-10           |
| Maximum Laser Diode Power (mW)                        | 20                | 120               | 50                | 20                | 80                | 150               |
| Output Power after Collimating Optics                 | 15                | 90                | 35                | 15                | 60                | 75                |
| Coll. Beam Dia. ( $1/e^2$ ) at Exit (mm) <sup>2</sup> | 2.6 x 1.0         | 2.3 x 1.1         | 2.6 x 1.0         | 2.8 x 1.1         | 2.8 x 1.3         | 2.2 x 1.3         |
| Collimated Beam Divergence (mrad) <sup>3</sup>        | <0.6              | <0.6              | <0.7              | <0.7              | <0.8              | <0.8              |
| M <sup>2</sup> (ModeMaster x,y)                       | ≤1.2              | ≤1.2              | ≤1.2              | ≤1.2              | ≤1.2              | ≤1.2              |
| Transverse Mode                                       | TEM <sub>00</sub> | TEM <sub>00</sub> | TEM <sub>00</sub> | TEM <sub>00</sub> | TEM <sub>00</sub> | TEM <sub>00</sub> |
| Polarization Extinction Ratio <sup>4</sup> (dB)       | >17               | >17               | >17               | >17               | >17               | >17               |
| Polarization Axis                                     | TE                | TE                | TE                | TE                | TE                | TE                |
| Warm-Up Time (sec.)                                   | 1                 | 1                 | 1                 | 1                 | 1                 | 1                 |

<sup>1</sup> Other wavelengths available. Please call for more details.

<sup>2</sup> Dependent on laser diode choice.

<sup>3</sup> Full angle. Dependent on laser diode choice.

<sup>4</sup> Typical.

# Lasiris ColdRay

## Compact Thermoelectrically Cooled Laser

|                                     |                                                       |                                                |
|-------------------------------------|-------------------------------------------------------|------------------------------------------------|
| <b>Optical Specifications</b>       | Power Stability in Operation 1 hr, 25°C (%)           | $\leq \pm 0.5$                                 |
|                                     | Power Stability in Operation 24 hrs, 25°C (%)         | $\leq \pm 1$                                   |
|                                     | Pointing Stability ( $\mu\text{rad}/^\circ\text{C}$ ) | $< 10$                                         |
|                                     | Noise                                                 |                                                |
|                                     | RMS, 20 Hz to 20 MHz (%)                              | $< 0.5$                                        |
| <b>Environmental Specifications</b> | P-P, 20 Hz to 20 MHz (%)                              | $< 1$                                          |
|                                     | Operating Temperature ( $^\circ\text{C}$ )            | 0 to 40                                        |
|                                     | Storage Temperature ( $^\circ\text{C}$ )              | -20 to 60                                      |
|                                     | Operating Pressure and Humidity                       | Atmospheric/non-condensing                     |
|                                     | Over-Temperature Protection, ESD Protection           |                                                |
| <b>Electrical Specifications</b>    | Input Voltage/Current (V)                             | 5 $\pm$ 10% at 2A max (typical 1A)             |
|                                     | Modulation                                            |                                                |
|                                     | Rise/fall time ( $\mu\text{s}$ )                      | 5 (faster response time available upon demand) |
|                                     | Reverse-Polarity Protection, Over-Voltage Protection  |                                                |
|                                     | Input Impedance ( $k\Omega$ )                         | $> 1$                                          |
| <b>Options</b>                      | Slow Start Time Delay (msec)                          | $< 250$ (most models)                          |


### Options

#### Power Options


Power Adjustment Potentiometer

Laser power can be easily changed by adjusting an optional built-in potentiometer with a small screwdriver. Code "P".

#### Power Adjustment by External Signal

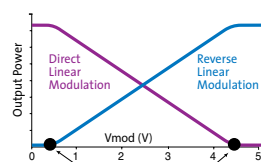

Standard Option -S  
Option -RS

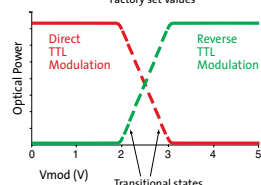

Option -T  
Option -RT

#### Standard modulation

S (synchro) or RS (reverse synchro): DC to 10 kHz

S: Input voltage = 0VDC  $\rightarrow$  laser "ON"

RS: Input voltage = 5VDC  $\rightarrow$  laser "ON"

T (TTL) or RT (Reverse TTL): DC to 10 kHz

Transition 2 to 3 VDC

T: Input voltage = 0VDC  $\rightarrow$  laser "ON"

RT: Input voltage = 5V  $\rightarrow$  laser "ON"


#### "Fast" Modulation

FS (fast synchro), FRS (fast reverse synchro) or


FT (fast TTL), FRT (fast reverse TTL)

DC to 300 kHz

### DB-9 Pinout and Back Panel Diagram

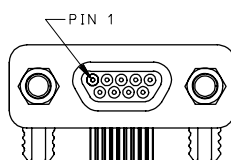

PIN 1: TEC Control Voltage (Black\*)  
PIN 2: Reserved for Calibration (Brown\*)  
PIN 3: Not Connected (Red\*)  
PIN 4: Current/Power Monitoring Output (Orange\*)  
PIN 5: Ground (Yellow\*)  
PIN 6: TEC Current Monitoring (Green\*)  
PIN 7: Modulation Input (Blue\*)  
PIN 8: Reserved for Calibration (Violet\*)  
PIN 9: Diode Temperature Monitoring (Gray\*)  
\*Wire color on optional Micro DB-9 cable

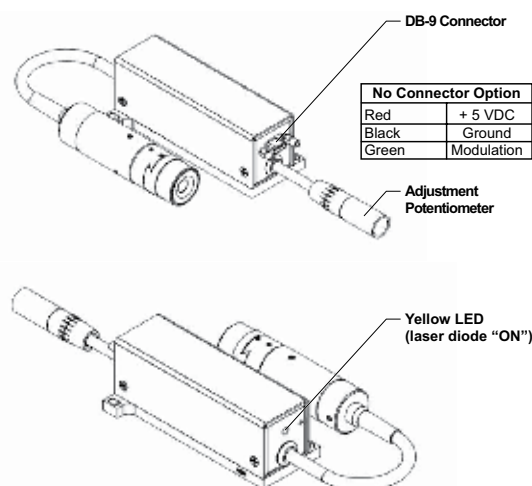

| No Connector Option |            |
|---------------------|------------|
| Red                 | + 5 VDC    |
| Black               | Ground     |
| Green               | Modulation |

# Lasiris PowerLine

Compact, High Power Laser Line Generator

Diode Laser  
Systems



Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers



## Mechanical Specifications

Top View

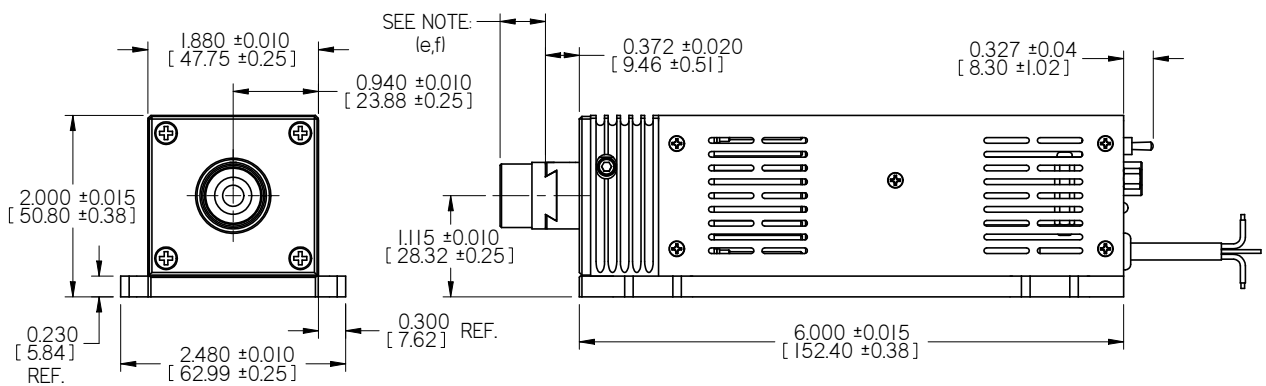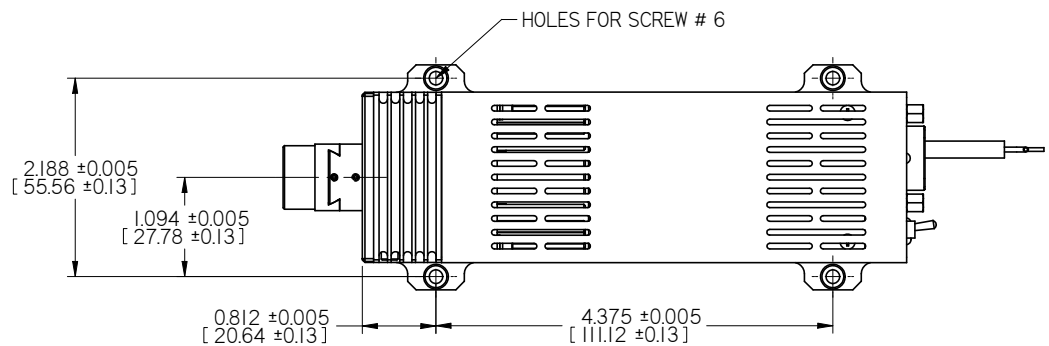

Front View

Side View

# Lasiris PowerLine

## Compact, High Power Laser Line Generator

### Key Features

- High power laser in a compact package
- Uniform intensity distribution for laser line generators
- Focusable
- High pointing stability
- Over-voltage, reverse-polarity, over-heating and ESD protection
- Compatible with standard pattern generators

### Ordering Information

For all PowerLine series, the warranty period shall be one (1) year. To order, use the following code: **MODEL-(Fast/Slow axis)-(Pattern)-(Interbeam Angle)-(Wavelength)(Electronic Option)-Power-Fan Angle-Custom Options**. E.g., **PL-501L-670T-500-45°**. Note that the projected fan angle may be less than the lens fan angle.

|                              | PowerLine UV                                        | PowerLine Violet          | PowerLine Blue            | PowerLine Red         | PowerLine IR              |
|------------------------------|-----------------------------------------------------|---------------------------|---------------------------|-----------------------|---------------------------|
| Wavelength <sup>1</sup> (nm) | 375                                                 | 405                       | 445                       | 670                   | 810                       |
| Diode Power (mW)             | 150                                                 | 600                       | 500                       | 500                   | 1000, 2000                |
| Beam Power (mW)              | 112                                                 | 450                       | 375                       | 375                   | 750, 1500                 |
| Electrical Power             | 12/24VDC, 3A <sup>2</sup>                           | 12/24VDC, 3A <sup>2</sup> | 12/24VDC, 3A <sup>2</sup> | 5VDC, 2A <sup>2</sup> | 5VDC, 3A, 4A <sup>2</sup> |
| Lens Fan Angle               | 10°, 15°, 20°, 30°, 40°, 45°, 55°, 60°, 75°, custom |                           |                           |                       |                           |

<sup>1</sup> ±10 nm.

<sup>2</sup> At ambient temperature of 23°C.

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories


# Lasiris PowerLine

## Compact, High Power Laser Line Generator

### Uniform Intensity

Conventional laser line patterns are often generated by cylindrical optics that produce a Gaussian line profile with a bright center and fading ends. Lasiris patented beam shaping optics spread the light into an evenly illuminated line. The result is a crisp, uniform line with sharp ends.

### Line Intensity Profile Along Line Length

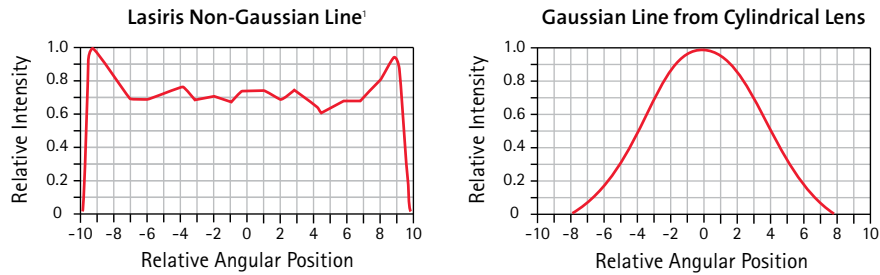


Relative Intensity vs. Angular Position Along Line Length

### Focusing Performance

The following figures show the typical focusing and depth-of-field performance (at  $1/e^2$ ) performance. Lasiris PowerLine lasers are focusable and can be adjusted by the user to produce a focused line at any projection distance. In addition, the line can be collimated so that its thickness remains fairly constant over a long projection distance.

### Focusing and Depth-of-Field Performance

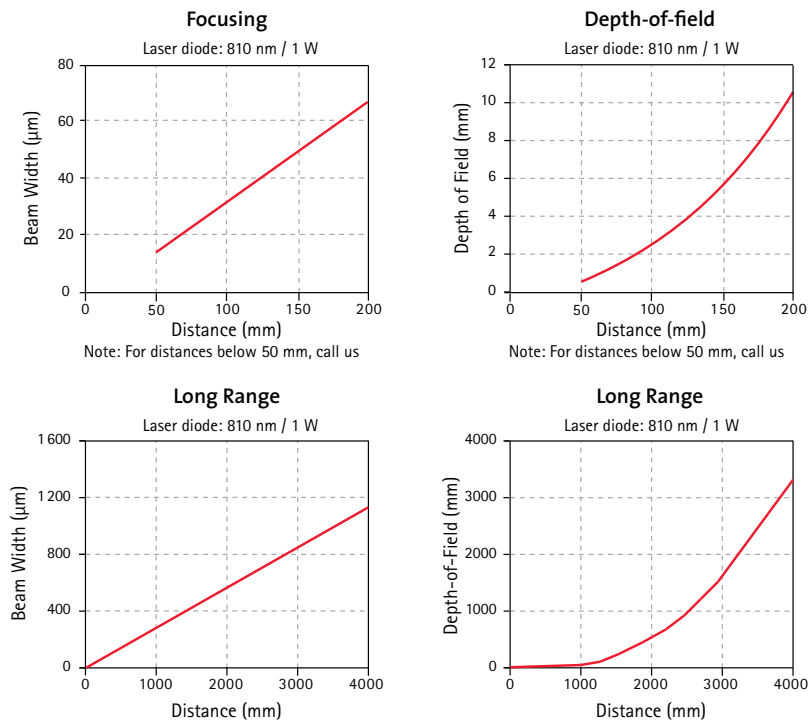

<sup>1</sup> Typical profile.

# Lasiris PowerLine

## Compact, High Power Laser Line Generator

|                                     |                                                          |                                                                                                               |
|-------------------------------------|----------------------------------------------------------|---------------------------------------------------------------------------------------------------------------|
| <b>Optical Specifications</b>       | Wavelength $\pm 10$ (nm)                                 | 375, 405, 445, 670, 810, custom                                                                               |
|                                     | Output Power (mW)                                        | 150 to 2000                                                                                                   |
|                                     | Fan Angles ( $^{\circ}$ )                                | 1 to 120                                                                                                      |
|                                     | Intensity Distribution                                   | Uniform (non-Gaussian) lengthwise, Gaussian widthwise                                                         |
|                                     | Line Thickness (focus)                                   | User adjustable                                                                                               |
|                                     | Pointing Stability ( $\mu\text{rad}/^{\circ}\text{C}$ )  | 5                                                                                                             |
|                                     | Bore Sighting (mrad)                                     | <3 (collimated)                                                                                               |
| <b>Environmental Specifications</b> | Operating Temperature ( $^{\circ}\text{C}$ )             | -20 to +55 for most models                                                                                    |
|                                     | Storage Temperature ( $^{\circ}\text{C}$ )               | -40 to +70                                                                                                    |
|                                     | Wavelength Drift                                         | Maximum $\pm 1$ nm over entire operating temperature range                                                    |
| <b>Electrical Specifications</b>    | Power Supply Voltage (VDC)                               | 5 $\pm$ 0.5<br>Options: 12/24 $\pm$ 0.5 or 110/240 VAC external adapter                                       |
|                                     | Power Supply Current (A)                                 | 3 at ambient temperature<br>4 maximum                                                                         |
|                                     | Built-In Protection                                      | Entire product: ESD, over-voltage, reverse polarity of power supply<br>Laser diode: overheating, over-current |
|                                     | Laser Diode Operating Temperature ( $^{\circ}\text{C}$ ) | 25 $\pm$ 0.5 (adjustable in factory)                                                                          |
|                                     | Maximum Beam Power                                       | User adjustable (trim potentiometer on the back panel)                                                        |
|                                     | Beam Modulation                                          | External, though a DB-9 connector on the back panel                                                           |
|                                     | Monitoring                                               | Laser temperature, laser current, PD current through DB-9                                                     |


### Options

#### Power Options

Power Adjustment Potentiometer

Laser power can be easily changed by adjusting an optional built-in potentiometer with a small screwdriver. Code "P".

Power Adjustment by External Signal

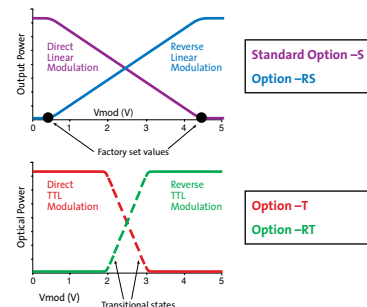

Standard modulation

S (synchro) or RS (reverse synchro): DC to 10 kHz

S: Input voltage = 0VDC  $\rightarrow$  laser "ON"

RS: Input voltage = 5VDC  $\rightarrow$  laser "ON"

T (TTL) or RT (Reverse TTL): DC to 10 kHz

Transition 2 to 3 VDC

T: Input voltage = 0VDC  $\rightarrow$  laser "ON"

RT: Input voltage = 5V  $\rightarrow$  laser "ON"

"Fast" Modulation

FS (fast synchro), FRS (fast reverse synchro) or

FT (fast TTL), FRT (fast reverse TTL)

DC to 300 kHz

Maximum Frequency (kHz)

10

Rise/Fall Time ( $\mu\text{s}$ )

<10

No fan for lower electronic noise (up to 40 $^{\circ}$  C without over-heating)

Base plate for efficient heat dissipation



### Available Patterns

Single Line

Parallel Lines



### DB-9 Pinout and Back Panel Diagram

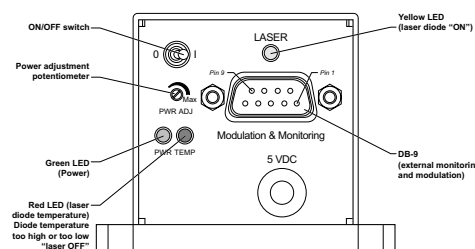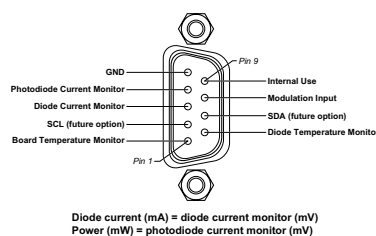

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories

# Lasiris Green PowerLine

Thermoelectrically Cooled Green Laser


Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers



## Mechanical Specifications

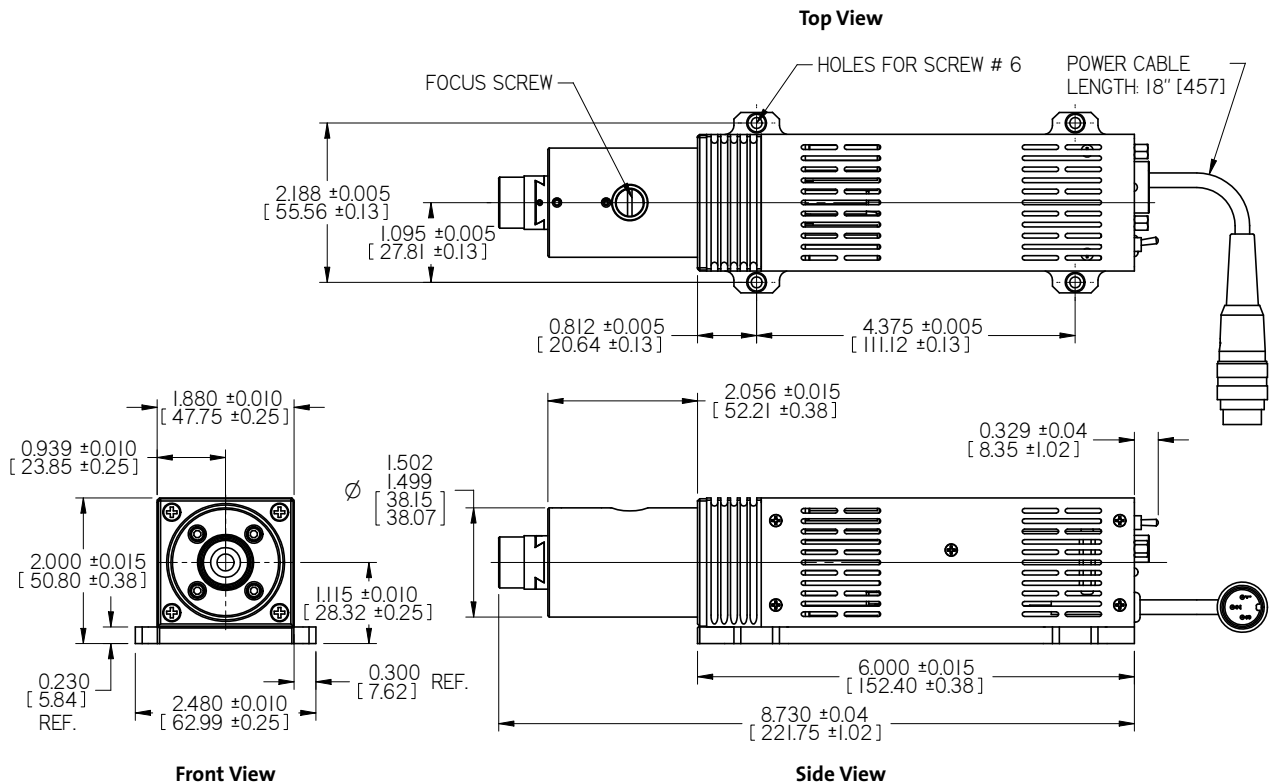

# Lasiris Green PowerLine

## Thermoelectrically Cooled Green Laser

### Key Features

- High power, high visibility, high contrast laser
- External focusing
- Uniform, non-Gaussian intensity distribution along the line
- ESD, over-temperature, and reverse-polarity protection

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

### Ordering Information

Green PowerLine lasers are covered under a one-year warranty (parts and labor). **MODEL-(Fast/Slow axis)(Pattern)-(Interbeam Angle)-("532" for Wavelength)(Electronic Option)-Power-Fan Angle-Custom Options.** E.g., GPL- 503L(1.2°)-532-100-5°. Contact us for more details.

|                              | 501L or 501D                                                                 | 501H      | 503L or 503D       | 505L or 505D      | 509L or 509D      |
|------------------------------|------------------------------------------------------------------------------|-----------|--------------------|-------------------|-------------------|
| Pattern                      | 1 line or 1 dot                                                              | crosshair | 3 lines or 3 dots  | 5 lines or 5 dots | 9 lines or 9 dots |
| Interbeam Angle <sup>1</sup> | -                                                                            | -         | 1.2°, 3.95°, 9.23° | 0.18°, 1.23°      | 0.09°, 0.06°      |
| Diode Power (mW)             | 50, 100, 150, 200                                                            |           |                    |                   |                   |
| Fan Angle                    | 1° <sup>2</sup> , 5°, 10°, 15°, 20°, 30°, 45°, 60°, 75°, custom <sup>3</sup> |           |                    |                   |                   |

|                              | 511L or 511D                                                                 | 515L or 515D        | 519L or 519D        | 533L or 533D        | 599L or 599D        |
|------------------------------|------------------------------------------------------------------------------|---------------------|---------------------|---------------------|---------------------|
| Pattern                      | 11 lines or 11 dots                                                          | 15 lines or 15 dots | 19 lines or 19 dots | 33 lines or 33 dots | 99 lines or 99 dots |
| Interbeam Angle <sup>1</sup> | 1.2°                                                                         | 1.8°                | 0.61°               | 0.07°, 0.30°        | 0.118°              |
| Diode Power (mW)             | 50, 100, 150, 200                                                            |                     |                     |                     |                     |
| Fan Angle                    | 1° <sup>2</sup> , 5°, 10°, 15°, 20°, 30°, 45°, 60°, 75°, custom <sup>3</sup> |                     |                     |                     |                     |

<sup>1</sup> At 532 nm.

<sup>2</sup> Not standard for crosshair projectors.

<sup>3</sup> Custom (please call us).


# Lasiris Green PowerLine

## Thermoelectrically Cooled Green Laser

### Uniform Intensity

Laser line patterns are often generated by cylindrical optics that produce a Gaussian line profile with a bright center and fading ends. Lasiris patented optics spread the light into an evenly illuminated line with sharp ends.

### Line Intensity Profile Along Line Length

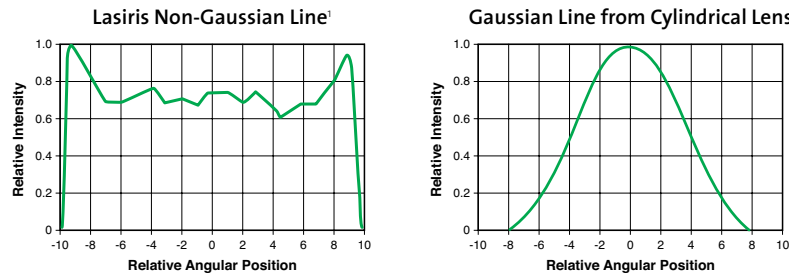


Relative Intensity vs. Angular Position Along Line Length


### Focusing Performance

The following figures show the typical focusing and depth-of-field performance of the Green PowerLine Laser. The focus charts indicate the minimum line thickness (at  $1/e^2$ ) achievable for a specific projection distance. The depth-of-field is defined as twice the distance over which the thickness of the line has increased by a factor of  $\sqrt{2}$ .

### Focusing and Depth-of-Field Performance

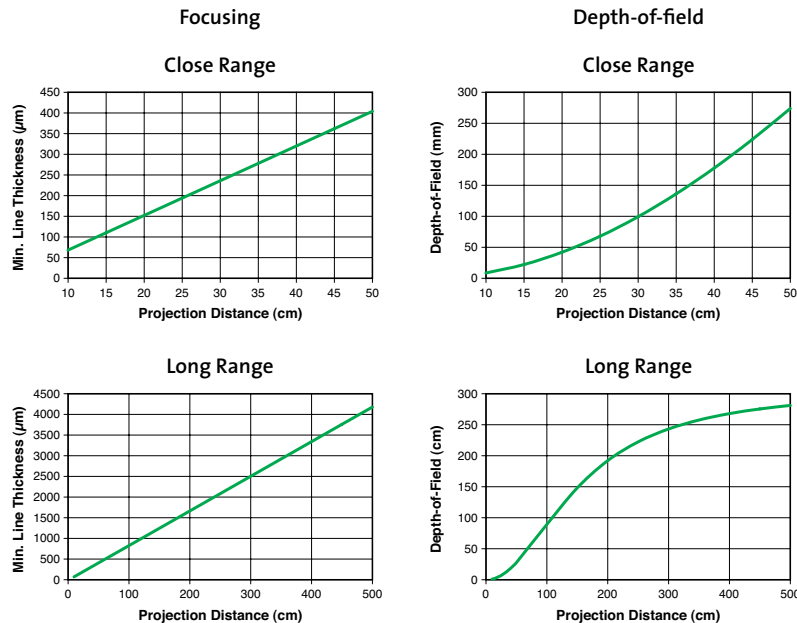

### DB-9 Pinout and Back Panel Diagram

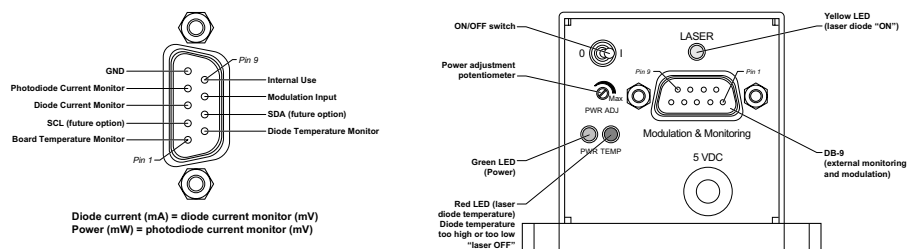

<sup>1</sup> Typical profile.

# Lasiris Green PowerLine

## Thermoelectrically Cooled Green Laser

|                                     |                                        |                                                                                                               |
|-------------------------------------|----------------------------------------|---------------------------------------------------------------------------------------------------------------|
| <b>Optical Specifications</b>       | Wavelength                             | 532                                                                                                           |
|                                     | Output Power (mW)                      | 50,100                                                                                                        |
|                                     | Intensity Uniformity                   | 30% uniformity over 80% of the pattern                                                                        |
|                                     | Mode                                   | TEM <sub>00</sub> (Gaussian profile)                                                                          |
|                                     | Polarization Ratio (typical)           | 4:1                                                                                                           |
|                                     | Line Thickness (focus)                 | See focusing charts                                                                                           |
|                                     | Bore Sighting (mrad)                   | <3                                                                                                            |
| <b>Environmental Specifications</b> | Operating Temperature (°C)             | -20 to +45                                                                                                    |
|                                     | Storage Temperature (°C)               | -40 to +70                                                                                                    |
|                                     | Wavelength Drift                       | Maximum ±1 nm over entire operating temperature range                                                         |
| <b>Electrical Specifications</b>    | Power Supply Voltage (VDC)             | 5 ±0.5<br>Adapter is available to supply the unit from 110/240 VAC line                                       |
|                                     | Power Supply Current (A)               | 3 at ambient temperature<br>4 maximum                                                                         |
|                                     | Built-In Protection                    | Entire product: ESD, over-voltage, reverse polarity of power supply<br>Laser diode: overheating, over-current |
|                                     | Laser Diode Operating Temperature (°C) | 25 ±0.5 (adjustable in factory)                                                                               |
|                                     | Maximum Beam Power                     | User adjustable (trim potentiometer on the back panel)                                                        |
|                                     | Beam Modulation                        | External, though a DB-9 connector on the back panel                                                           |
|                                     | Monitoring                             | Laser temperature, laser current, PD current through DB-9                                                     |


### Options

#### Power Options


##### Power Adjustment Potentiometer

Laser power can be easily changed by adjusting an optional built-in potentiometer with a small screwdriver. Code "P".

##### Power Adjustment by External Signal

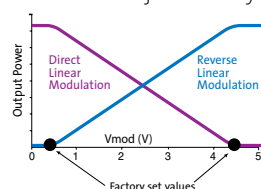

Standard Option -S  
Option -RS

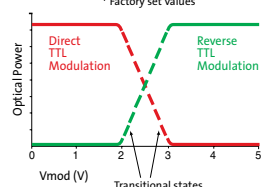

Option -T  
Option -RT

#### Standard modulation

S (synchro) or RS (reverse synchro): DC to 10 kHz

S: Input voltage = 0VDC → laser "ON"

RS: Input voltage = 5VDC → laser "ON"

T (TTL) or RT (Reverse TTL): DC to 10 kHz

Transition 2 to 3 VDC

T: Input voltage = 0VDC → laser "ON"

RT: Input voltage = 5V → laser "ON"

#### "Fast" Modulation

FS (fast synchro), FRS (fast reverse synchro) or

FT (fast TTL), FRT (fast reverse TTL)

DC to 300 kHz

#### Pulsed Modulation Available for both options

Maximum Frequency (kHz)

<100

Rise/Fall Time (μs)

<1

High Polarization Ratio

100:1

Separate Driver

Optional

### Available Patterns

Single Line



Crosshair



Parallel Lines



Single Square



Single Dot



Dot Line



7x7 Dot Matrix



Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

# Lasiris Magnum II

High Power Laser Line Generator


Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers



## Mechanical Specifications

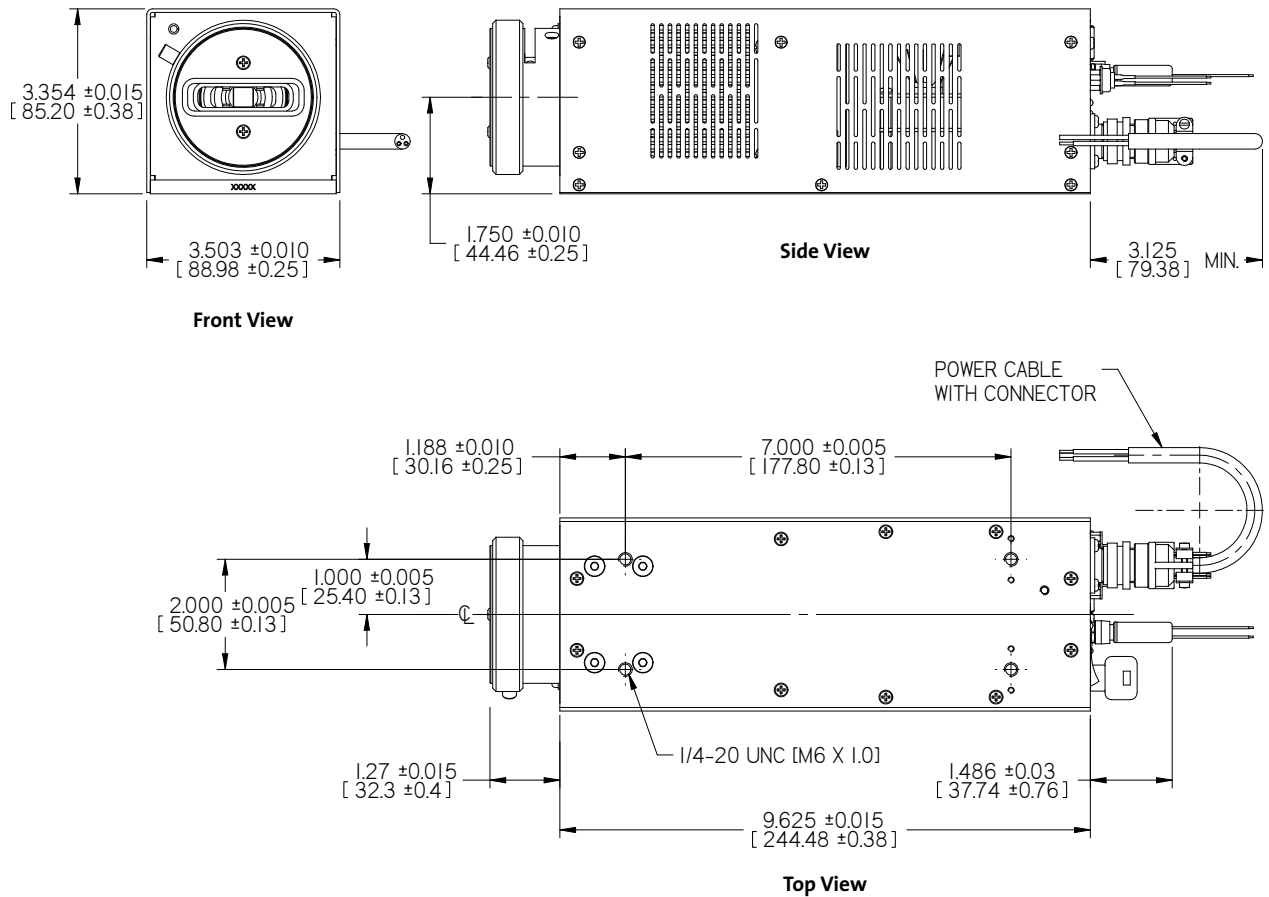

# Lasiris Magnum II

## High Power Laser Line Generator

### Key Features

- Uniform, non-Gaussian intensity distribution along the line
- Power to 7W
- Wide range of powers and fan angles
- High pointing stability
- Focusable
- Protection against over-voltage, reverse polarity of power supply, overheating and ESD
- Rugged, industrial-grade design

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories

### Ordering Information

For all Magnum series, the warranty period shall be four thousand (4000) hours or one (1) year (whichever comes first) for the laser diode chip and one (1) year for the electronic, optical, and mechanical components. To order, use the following code: **MODEL-Wavelength(Electronic Option)-Power-Fan Angle**. Add -100K after the Pulsing Option code for high-speed frequency (100 KHz). See Modulation Graphs for definitions. Add -SD to product code for the separate driver option (e.g., **MAG2-810S-2000-45**). Note that the projected fan angle may be less than the lens fan angle.

|                              | Magnum 500                                          | Magnum 500 | Magnum 1000 | Magnum 1500 | Magnum 2000 |
|------------------------------|-----------------------------------------------------|------------|-------------|-------------|-------------|
| Wavelength <sup>1</sup> (nm) | 635                                                 | 680        | 680         | 665         | 810         |
| Diode Power (mW)             | 500                                                 | 500        | 1000        | 1500        | 2000        |
| Beam Power (mW)              | 400                                                 | 400        | 800         | 1200        | 1600        |
| Electrical Power             | 12 VDC, 3A                                          | 12 VDC, 3A | 12 VDC, 3A  | 12 VDC, 3A  | 12VDC, 5A   |
| Lens Fan Angle <sup>2</sup>  | 10°, 15°, 20°, 30°, 40°, 45°, 55°, 60°, 75°, custom |            |             |             |             |

|                              | Magnum 4000                                         | Magnum 7000 |
|------------------------------|-----------------------------------------------------|-------------|
| Wavelength <sup>1</sup> (nm) | 810                                                 | 808         |
| Diode Power (mW)             | 4000                                                | 7000        |
| Beam Power (mW)              | 3200                                                | 5600        |
| Electrical Power             | 12 VDC, 5A                                          | 12 VDC, 5A  |
| Lens Fan Angle <sup>2</sup>  | 10°, 15°, 20°, 30°, 40°, 45°, 55°, 60°, 75°, custom |             |

<sup>1</sup> ±10 nm.

<sup>2</sup> 10° and 15° not available on 4000 mW and 7000 mW models.


# Lasiris Magnum II

## High Power Laser Line Generator

### Uniform Intensity

Conventional laser line patterns are often generated by cylindrical optics that produce a Gaussian line profile with a bright center and fading ends. Lasiris patented beam shaping optics spread the light into an evenly illuminated line. The result is a crisp, uniform line with sharp ends.

### Line Intensity Profile Along Line Length

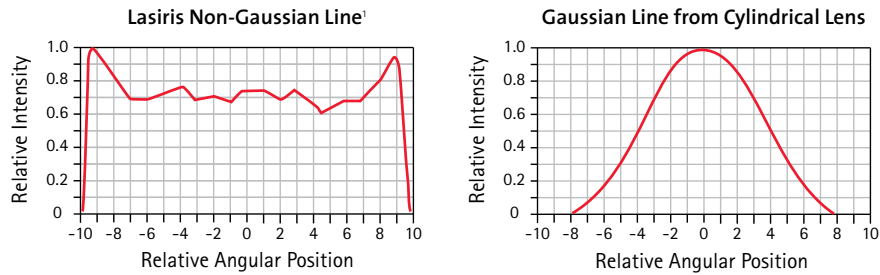


Relative Intensity vs. Angular Position Along Line Length

### Focusing Performance

The following figures show the typical focusing and depth-of-field performance (at  $1/e^2$ ). Lasiris Magnum II lasers are focusable and can be adjusted by the user to produce a focused line at any projection distance. In addition, the line can be collimated so that its thickness remains fairly constant over a long projection distance.

### Focusing and Depth-of-Field Performance

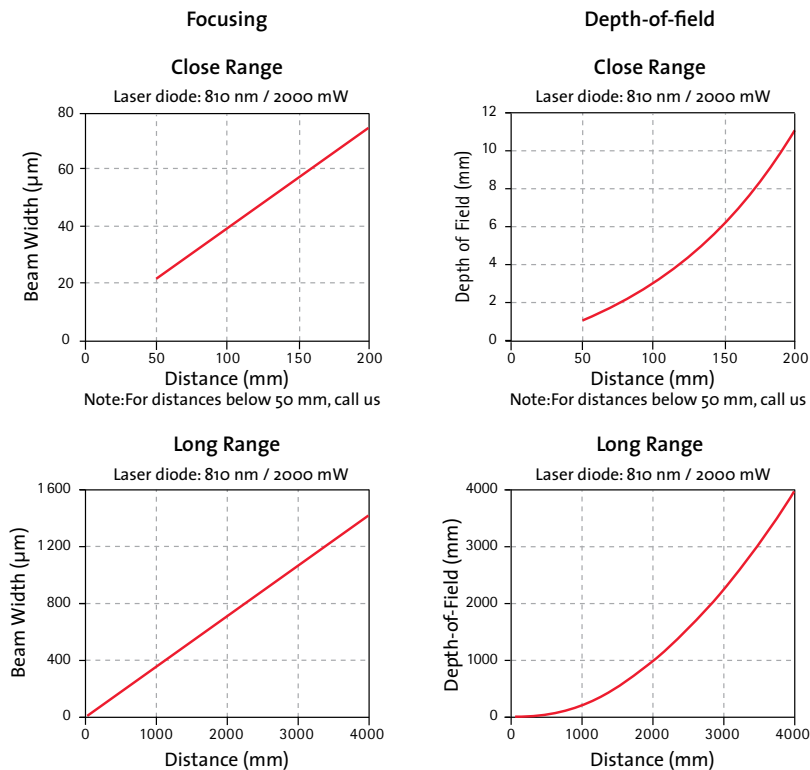


¹ Typical profile.

# Lasiris Magnum II

## High Power Laser Line Generator

|                                     |                                        |                                                                                                                          |
|-------------------------------------|----------------------------------------|--------------------------------------------------------------------------------------------------------------------------|
| <b>Optical Specifications</b>       | Wavelength (nm)                        | 635, 665, 680, 808, 810, custom                                                                                          |
|                                     | Output Power (mW)                      | 500 to 7000                                                                                                              |
|                                     | Intensity Distribution                 | Uniform (non-Gaussian) lengthwise, Gaussian widthwise                                                                    |
|                                     | Fan Angles' (°)                        | 10, 15, 20, 30, 40, 45, 55, 60, 75                                                                                       |
|                                     | Line Thickness (focus)                 | User adjustable                                                                                                          |
|                                     | Noise (%)                              | <1 at ambient temperature (Max. RMS variation)                                                                           |
|                                     | Power Stability (%)                    | <3 at ambient temperature (Max. variation in 24 hrs)                                                                     |
|                                     | Pointing Stability (μrad/°C)           | 10                                                                                                                       |
|                                     | Bore Sighting (mrad)                   | <3 (collimated)                                                                                                          |
| <b>Environmental Specifications</b> | Operating Temperature (°C)             | -35 to +50 for most models                                                                                               |
|                                     | Storage Temperature (°C)               | -40 to 60                                                                                                                |
|                                     | Wavelength Drift                       | Maximum ±1 nm over entire operating temperature range                                                                    |
| <b>Electrical Specifications</b>    | Power Supply Voltage (VDC)             | 12 ±0.5<br>An adapter is available to supply the unit from 110/240 VAC line                                              |
|                                     | Power Supply Current (A)               | 3 to 5 depending on laser power                                                                                          |
|                                     | Built-In Protection                    | Entire product: ESD, over-voltage up to 20V, reverse polarity of power supply.<br>Laser diode: overheating, over-current |
|                                     | Laser Diode Operating Temperature (°C) | 25 ±1 (adjusted in factory)                                                                                              |
|                                     | Maximum Beam Power                     | User adjustable (trim potentiometer on the back panel)                                                                   |
|                                     | Beam Modulation                        | External, through a DB-9 connector on the back panel                                                                     |
|                                     | Input Impedance (kΩ)                   | 1                                                                                                                        |
|                                     | Monitoring                             | Laser temperature, laser current, PD current through DB-9                                                                |
| <b>Options</b>                      | Power Options                          |                                                                                                                          |
|                                     | Power Adjustment Potentiometer         | Laser power can be easily changed by adjusting an optional built-in potentiometer with a small screwdriver. Code "P".    |

### Power Adjustment by External Signal


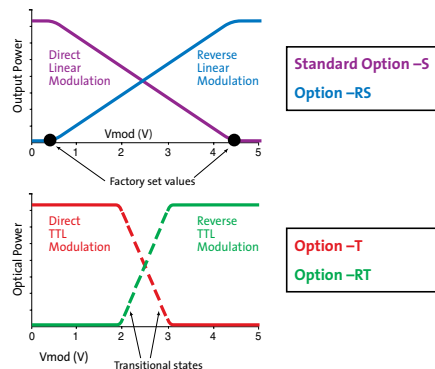


Standard modulation  
S (synchro) or RS (reverse synchro): DC to 10 kHz  
S: Input voltage = 0VDC → laser "ON"  
RS: Input voltage = 5VDC → laser "ON"

T (TTL) or RT (Reverse TTL): DC to 10 kHz  
Transition 2 to 3 VDC  
T: Input voltage = 0VDC → laser "ON"  
RT: Input voltage = 5V → laser "ON"

"Fast" Modulation  
FS (fast synchro), FRS (fast reverse synchro) or  
FT (fast TTL), FRT (fast reverse TTL)  
DC to 300 kHz

### DB-9 Pinout and Back Panel Diagram

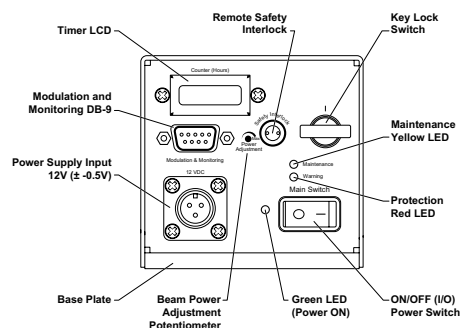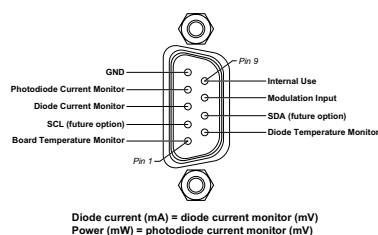

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories

# Lasiris Flat Top Projector

Gaussian to Flat Top Using Refractive Beam Shaping Optics

Diode Laser Systems

Diode Laser Modules


Lasiris Structured Light Lasers



## Mechanical Specifications


### PureBeam/Fiber Option

#### Side View

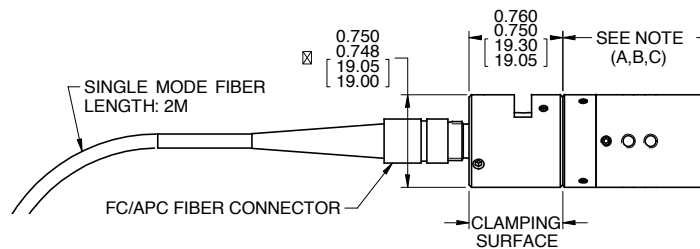


### ColdRay Laser Option

#### Side View

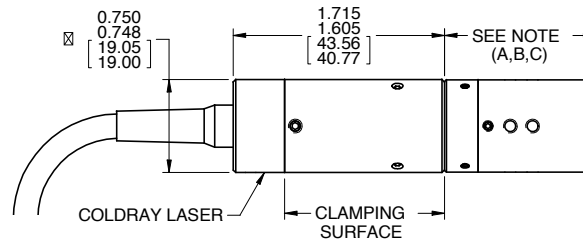


### Free Space Option

M14 x 0.5 THREAD  
X 0.375 [9.50] DEEP.

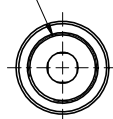

Front View

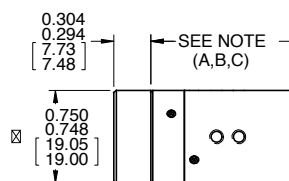

Side View

# Lasiris Flat Top Projector

## Gaussian to Flat Top Using Refractive Beam Shaping Optics

### Key Features

- Very efficient all glass optics
- Achromatic (from UV to near IR)
- Very smooth low ripple uniformity profile with no high frequency noise
- Uniformity coefficient of variation of better than 1% achievable
- High aspect ratios available off the shelf
- High uniformity edge steepness

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

### Ordering Information

To order, use the following code: **BSO-MODEL-Wavelength-POWER-Type-FT Size-WD** (Note: Size in degrees for DFT; size in  $\mu\text{m}$  for FFT; size in mm for WD). E.g., **BSO-FS-660-60-DFT-1x3** (Free Space for 660 nm 60 mW Diverging with a FT size of 1 x 3 degrees). E.g., **BSO-CR-660-60-FFT-300x900-100** (ColdRay for 660 nm 60 mW Focusing with a FT size of 300 x 900  $\mu\text{m}$  at 100 mm).

#### Input Type

|                      |                                                                                                                                                               |
|----------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------|
| ColdRay Laser Option | Coherent TEC direct semiconductor free space laser<br>Order code: CR - Wavelength - Power E.g., CR-640-80)                                                    |
| Free Space Option    | M14 threaded holder. Requires input beam dimension (IBD):<br>0.5 mm up to 4 mm Gaussian at $1/e^2$<br>Order code: FS-IBD (E.g., FS-4)                         |
| Fiber Option         | FC/PC input end, 2 m shielded, NA: 0.12, FC/APC exit end.<br>Requires operating wavelength between 405-830 nm<br>Order code: FIB - Wavelength (E.g., FIB-640) |

#### Flat Top Type

|                                    |                                                                                             |
|------------------------------------|---------------------------------------------------------------------------------------------|
| One Dimensional Diverging Flat Top | Order code: 1DDFT - Divergence (E.g., 1DDFT-0.3°)                                           |
| Two Dimensional Diverging Flat Top | Order code: 2DDFT - Divergence 1 x Divergence 2<br>(E.g., 2DDFT-0.3°x2°)                    |
| Imaged Flat Top                    | Order code: FFT - Divergence 1 x Divergence 2 - Focal Length<br>(E.g., FFT-0.3°x0.2° - F50) |

|                    |                                                                    |
|--------------------|--------------------------------------------------------------------|
| Beam Splitter Type | Order code: Number of beams - Period (E.g., 19B-50 $\mu\text{m}$ ) |
|--------------------|--------------------------------------------------------------------|

# Lasiris Flat Top Projector

## Gaussian to Flat Top Using Refractive Beam Shaping Optics

### Flat Top Options

The flat top is a diverging optical element that provides a uniform intensity profile in the far field. It can be imaged or offered simply as a diverging optical element.

Being cylinder based, it can be readily offered acting on only one axis or crossed pairs creating the desired aspect ratio for your application.

The dimension of the flat top at the focal plane of an imaging lens,  $f$ , is directly proportional to its effective focal length and can be calculated by:

$$FTS \text{ (flat top size)} = 2 * \tan\left(\frac{DIVERGENCE}{2}\right) * f$$

The flat top can be further enhanced by adding Coherent's series of binary phase beam splitters to create a complex array of flat tops. The beam splitters are highly stable binary masters etched directly in fused silica offering very low stray light. They have excellent beam to beam uniformity and have efficiencies in the order of 80%.

|                                         |                                                           |
|-----------------------------------------|-----------------------------------------------------------|
| Divergences <sup>1</sup> ( $\pm 10\%$ ) | 0.2, 0.3, 0.6, 1, 2, 3, 5, 8, 10, 15, 20, 30              |
| Focal Length <sup>2</sup>               | 14, 17.5, 20, 25, 30, 35, 40, 45, 50, 60, 75, 80, 90, 100 |

### Gratings for Beam Splitting Option

Complex flat top patterns can be achieved using our beam splitting option.

The angle between each beam is operating wavelength dependent and can be found via the grating formula:  
 $d \sin \theta = n \lambda$  ( $d$ =period in  $\mu\text{m}$ ,  $\theta$  angle in degrees,  $n$  = diffraction order, 0, +/-1, +/- 2, ....)

|                                        |                         |
|----------------------------------------|-------------------------|
| 2 Beams <sup>3</sup> ( $\mu\text{m}$ ) | 3.33, 7.77, 26.3, 338.6 |
| 3 Beams ( $\mu\text{m}$ )              | 3.33, 7.77, 26.3, 338.6 |
| 5 Beams ( $\mu\text{m}$ )              | 25, 169                 |
| 9 Beams ( $\mu\text{m}$ )              | 338.2, 532              |
| 11 Beams ( $\mu\text{m}$ )             | 25.7                    |
| 15 Beams ( $\mu\text{m}$ )             | 16.4                    |
| 19 Beams ( $\mu\text{m}$ )             | 50                      |
| 33 Beams ( $\mu\text{m}$ )             | 110, 426.1              |
| 99 Beams ( $\mu\text{m}$ )             | 260                     |


<sup>1</sup> Available divergences for one and two dimensional diverging flat tops.

<sup>2</sup> Available focal lengths for the imaged flat tops option ( $\lambda > 400 \text{ nm}$ ). Working distance will be a few mm shorter.


<sup>3</sup> Zero order minimized.

### Flat Top Generators


#### Flat Top 2D

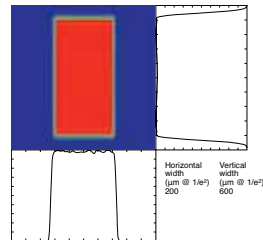


#### Flat Top Matrix 2D

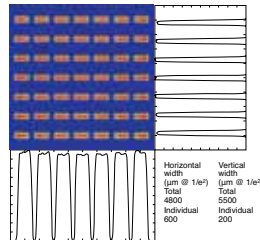

#### Flat Top 1D

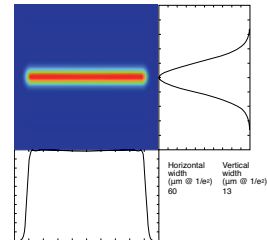

#### Flat Top Matrix 1D

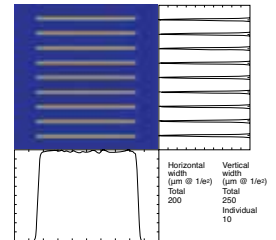


# Lasiris Flat Top Projector

## Gaussian to Flat Top Using Refractive Beam Shaping Optics

### Flow Cytometry Option

Coherent has been very successful in promoting its flat top for flow cytometry applications. Our standard flow cytometry option is available with the following specifications:

|                                                 |                                   |
|-------------------------------------------------|-----------------------------------|
| Working Distance (mm)                           | 45                                |
| Flat Top Dimension ( $\mu\text{m}$ at $1/e^2$ ) | 60                                |
| Gaussian Dimension ( $\mu\text{m}$ at $1/e^2$ ) | 10                                |
| Contained Power (%)                             | >60                               |
| CV Uniformity (%)                               | <1 fiber coupled, <1.5 free space |

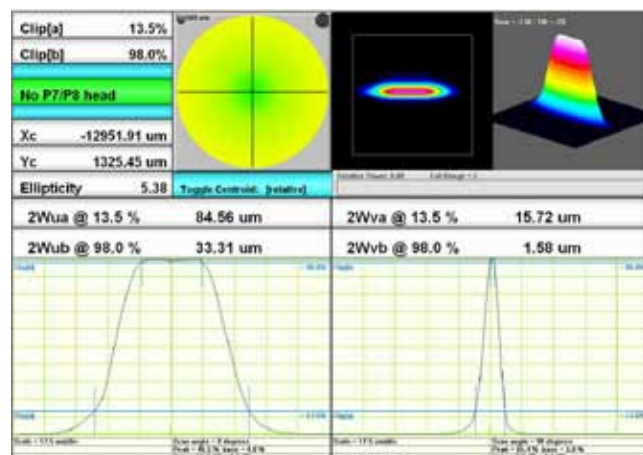

### Input Beam Characteristics

|                              |                             |
|------------------------------|-----------------------------|
| Dimensions <sup>1</sup> (mm) | 0.5 to 4                    |
| Wavelength <sup>2</sup> (nm) | 400 to 1200                 |
| Gaussian Fit (%)             | 97                          |
| Mode                         | TEM <sub>00</sub>           |
| Divergence                   | Collimated/Quasi collimated |

### Output Flat Top Characteristics

|                                                   |                                                   |
|---------------------------------------------------|---------------------------------------------------|
| High Efficiency 1 Dimensional Flat Top Uniformity | Contained power typical >85%                      |
| High Uniformity 1 Dimensional Flat Top Uniformity | CV uniformity typical <5%, down to <1% possible   |
| 2 Dimensional Flat Top                            | CV uniformity < $\pm 20\%$ , contained power >85% |
| Transmission Diverging Flat Top                   | Contained power >95%                              |
| Transmission Focused Flat Top                     | Contained power >90%                              |

<sup>1</sup> Specified to within 5%, larger input beam available.

<sup>2</sup> Other wavelengths available.

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories

# Interference Filters

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers



## Technical Specifications

|                                                       |                                                                                                                   |
|-------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------|
| Transmission at Laser Line (%)                        | 80, minimum (85, typical)                                                                                         |
| Out of Band Transmission (%)                          | 0.001 (50 db attenuation)                                                                                         |
| Spectral Range                                        | X-Ray to 1150 nm                                                                                                  |
| Effective Index                                       | 2.1                                                                                                               |
| Operating Temperature                                 | -20°C to +80°C (-4°F to +176°F)                                                                                   |
| Humidity Resistance                                   | Minimum 40 Cycles (MIL-STD-810E)<br>Method 507 Procedure 1                                                        |
| Dimension (mm)                                        | 25.4 (24.4 mm without ring)                                                                                       |
| Useable Aperture (mm)                                 | 23 ±0.4                                                                                                           |
| Thickness (mm)                                        | 4 ±0.25 (3 mm without ring)                                                                                       |
| Available Band Pass (at FWHM)                         | 20 nm (0 to 15° angle of incidence)<br>40 nm (0 to 25° angle of incidence)<br>70 nm (0 to 35° angle of incidence) |
| Available Center Wavelengths (nm)                     | 532, 635, 660, 670, 680, 685, 780, 810, 830                                                                       |
| Bandwidth Tolerances (nm)                             | 20 ±2<br>40 ±8<br>70 ±14                                                                                          |
| Center Wavelength Tolerances (dependent on bandwidth) | +2/-2<br>+5/-5<br>+10/-10                                                                                         |

Accessories

# Interference Filters

## Key Features

- 80% peak transmission
- Non-hydroscopic
- Wide bandwidth

## Angle of Incidence and Blue Shifting

The center wavelength of the filter will shift toward shorter (blue) wavelengths as the angle of incidence increases; the greater the angle the greater the shift (see % shift vs. angle of incidence graph). This is known as blue shifting. The amount of blue shift may be calculated using the following formula:


$$\lambda_1 = \lambda_0 (1 - \sin^2 \theta / n^2)^{1/2}$$

where  $\theta$  = angle of incidence

$\lambda_1$  = resulting wavelength at angle  $\theta$

$\lambda_0$  = central wavelength at normal incidence

$n$  = effective index of the filter

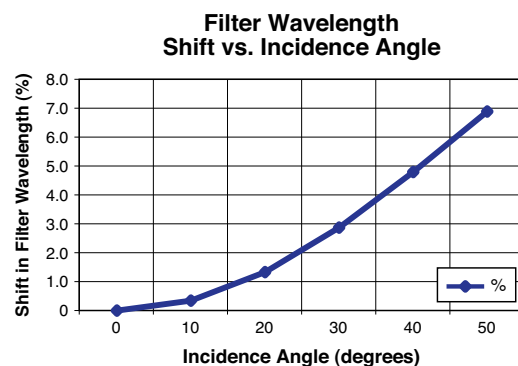

## Laser Diode Tolerances

Another aspect to take into account when selecting a filter is laser diode tolerances. Diodes typically have a wavelength tolerance of  $\pm 10$  nm; this wavelength (measured at room temperature) will increase with temperature increases and decrease with temperature decreases at a rate of 0.25 nm/°C.

## Ordering Information

To order a filter, choose the center wavelength and bandwidth from the Specifications on page 1. The order code is:

**<center wavelength> CW - <band pass> BP - 25.4-80%**

*For example: 635CW-20BP-25.4-80%*

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories

# Diode LabLaser Universal Power Supply



## Key Features

- Fully-featured LabLaser 6 VDC power supply
- Variable laser power control
- Keylock, emission delay and remote interlock

The Diode LabLaser™ Universal Power Supply provides a well-regulated 6 VDC, 250 mA output for driving any of the Coherent Diode LabLasers with the additional benefit of variable power control. To vary the output power of the Diode LabLasers which incorporate the Variable Power and Modulation Control circuitry simply adjust the multi-turn knob on the front panel of the 31-1050 power supply. This varies a voltage on the “third wire” connection in the laser cables and the laser power can be smoothly varied between threshold and the rated power. The supply incorporates CDRH safety features (keylock, power-on indicator, 4-6 second

emission time delay and remote interlock socket), is CE-marked and meets appropriate UL, CSA, and EN requirements.

This compact (120 x 76 x 48 mm) power supply can be connected to any 90-130 or 180-260 VAC 50/60 Hz supply. There is an IEC 320 input socket and a 1.8 m detachable line cord with plug provided. The output to the laser features an easy connect/disconnect 3-terminal jack socket on the front panel which mates to the male plugs of the cables connected to all Diode LabLasers.

|                       |             |                       |                                                                              |
|-----------------------|-------------|-----------------------|------------------------------------------------------------------------------|
| <b>Specifications</b> | 31-1050-000 | LabLaser Power Supply | Power Supply, LabLaser, 6VDC, 250 mA, Adjustable Power Control for MVP       |
|                       | 1057734     | Plug, 3.5 mm Stereo   | Phono Plug to take Laser Diode Module with leads and add a Power Supply Plug |

# Universal Diode Laser Module Power Supply



## Key Features

- Well-regulated 5VDC output
- Autoranging 90 to 260 VAC IEC input
- Output cable with easy connect/disconnect jack

The Universal Diode Laser Module Power Supply uses the latest switchmode technology to provide a well-regulated 5 VDC output suitable for driving any of the Coherent laser diode modules. The supply is CE-marked and meets appropriate UL, CSA, TUV, GS, and FCC requirements.

This compact (75 x 50 x 41 mm) autoranging power module can be connected to any 90-260 VAC 50/60 Hz supply. There is an IEC 320 input socket and a 1.8 m detachable line cord with plug provided. The 1.8 m output cable is terminated with an easy connect/disconnect jack socket which mates to the male plugs connected to all laser modules.

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

Accessories

## Specifications

|             |                     |                                                                              |
|-------------|---------------------|------------------------------------------------------------------------------|
| 31-1001-000 | LDM Power Supply    | Power Supply, Laser Diode Module, Desk Top, 5VDC , 1.5 amp                   |
| 1057734     | Plug, 3.5 mm Stereo | Phono Plug to take Laser Diode Module with leads and add a Power Supply Plug |


# Variable Angle Diode Laser Mounts

Diode Laser  
Systems

Diode Laser  
Modules

Lasiris  
Structured  
Light Lasers

Accessories

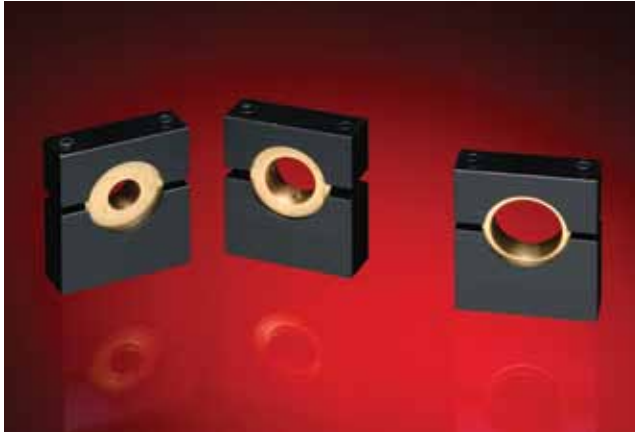


## Key Features

- Patent pending
- Simple tip/tilt and lock in place mechanism
- Very affordable
- Adjustable angle through about  $\pm 15$  degrees

Coherent's variable angle laser diode module mounts are designed to accept three different sizes of diode laser module namely VLM3, VLM2 and LabLasers (9.5 mm, 14.8 mm and 19.0 mm diameter).

This simple mount acts as combined heat sink, electrical insulator and variable-angle mount for 9.5 or 14.8 mm diameter diode modules and 19.0 mm Diode LabLasers. The lasers are held firmly in a spherical-surface brass ring by two small set screws. This ring is clamped between two insulating blocks forming a "ball and socket" arrangement. Two socket-head screws can be loosened to allow easy angular adjustment of the module ( $\pm 15^\circ$  range) and quick, rigid locking in position.

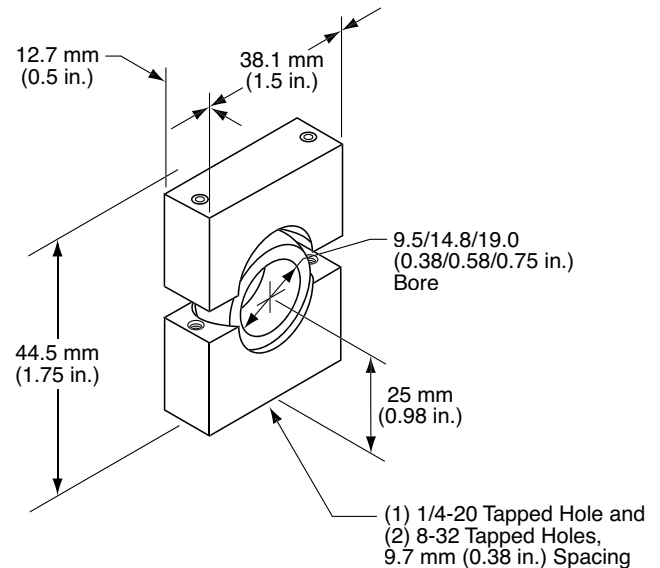

| Specifications |                              |                                                |  |
|----------------|------------------------------|------------------------------------------------|--|
| 0221-437-00    | Mount for VLM2, FVLM2 or LG2 | Variable Angle, 1/4-20 inch Base, 14.8 mm Bore |  |
| 0221-448-00    | Mount for VLM3 or LG3        | Variable Angle, 1/4-20 inch Base, 9.5 mm Bore  |  |
| 0221-449-00    | Mount, LabLaser/Industrial   | Variable Angle, 1/4-20 inch Base, 19 mm Bore   |  |

# FieldMaxII-TO

## Laser Power Meter



FieldMaxII™-TO is an affordable, versatile, easy-to-use digital power meter designed for field service and production test applications. FieldMaxII-TO features a large, easy-to-read backlit LCD and an intuitive user interface offering direct button-driven control for simple operation. To facilitate laser tuning, FieldMaxII-TO

### Key Features

- Measures laser power from nW to kW
- Works with Thermal and Optical sensors
- Large, bright, backlit LCD display
- Simulated analog-like movement for laser tuning
- Intuitive button-driven user interface
- USB 1.1 interface

provides a fast and responsive tuning feature. Statistics are supported, including min., max., and mean. FieldMaxII-TO is compatible with our thermal and optical sensors. It can measure power from nW to kW.

| Specifications | 1098579 | FieldMaxII-TO | Laser Power Meter |
|----------------|---------|---------------|-------------------|
|----------------|---------|---------------|-------------------|

Diode Laser Systems

Diode Laser Modules

Lasiris Structured Light Lasers

# LaserCheck

## Handheld Laser Power Meter



LaserCheck™ is a handheld, inexpensive laser power meter specifically designed to provide power measurements in a small, lightweight, self-contained package that can easily be stored in a pocket or tool kit. At higher powers the silicon sensor is protected from saturation

### Key Features

- Auto-ranging with peak sample-and-hold
- Power measurement from 0.5  $\mu$ W to 1W
- $\pm 5\%$  calibrated accuracy from 400 nm to 1064 nm

with a built-in optical attenuator, and all measurements are displayed on a 3-digit LCD screen. LaserCheck is microprocessor-controlled with wavelength correction, auto-ranging ( $\mu$ W or mW displayed), optical attenuator, power overload warning and automatic shutoff.

| Specifications | 1098293 | LaserCheck | Handheld Laser Power Meter |
|----------------|---------|------------|----------------------------|
|----------------|---------|------------|----------------------------|

Accessories

# High-Sensitivity Thermopile Sensors

To 1W



## Key Features


- Thermally stabilized designs
- Broadband 0.3 to 11  $\mu\text{m}$  spectral range
- High resolution to 10  $\mu\text{W}$
- Large area: 10 and 19 mm diameter
- Fiber-optic connector options

The PS model sensors offer a thermally stabilized, amplified power sensor with a broad spectral response, high sensitivity, and large active area. These sensors are ideal for laser diodes, diode bars and stacks, HeNe,


HeCd and small ion lasers. For optimum performance at low power levels we recommend the PS10Q and PS19Q, which utilize quartz windows. PS10 models include a detachable beam tube to block stray light.

| System Specifications                                              | PS10                    | PS10Q                   | PS19                    | PS19Q                   |
|--------------------------------------------------------------------|-------------------------|-------------------------|-------------------------|-------------------------|
| Wavelength Range ( $\mu\text{m}$ )                                 | 0.3 to 11               | 0.3 to 2                | 0.3 to 11               | 0.3 to 2                |
| Power Range                                                        | 100 $\mu\text{W}$ to 1W | 100 $\mu\text{W}$ to 1W | 100 $\mu\text{W}$ to 1W | 100 $\mu\text{W}$ to 1W |
| Resolution ( $\mu\text{W}$ )                                       | 10                      | 10                      | 10                      | 10                      |
| Max. Thermal Drift                                                 | 40 $\mu\text{W}$        | 20 $\mu\text{W}$        | 400 $\mu\text{W}$       | 20 $\mu\text{W}$        |
| Max. Avg. Power Density ( $\text{kW}/\text{cm}^2$ )                | 0.5                     | 0.5                     | 0.5                     | 0.5                     |
| Max. Pulse Energy Density ( $\text{mJ}/\text{cm}^2$ , 10 ns pulse) | 50                      | 50                      | 50                      | 50                      |
| Response Time (sec.)                                               | 2                       | 2                       | 2                       | 2                       |
| Detector Coating                                                   | Black                   | Black                   | Black                   | Black                   |
| Quartz Filter Window                                               | -                       | Yes                     | -                       | Yes                     |
| Detector Diameter (mm)                                             | 10                      | 10                      | 19                      | 19                      |
| Dimensions (mm)                                                    | $\Phi 84 \times 49$     | $\Phi 84 \times 49$     | $\Phi 84 \times 49$     | $\Phi 84 \times 49$     |
| Calibration Uncertainty                                            | $\pm 1\%$               | $\pm 1\%$               | $\pm 1\%$               | $\pm 1\%$               |
| Calibration Wavelength (nm)                                        | 514                     | 514                     | 514                     | 514                     |
| Cooling Method                                                     | Air                     | Air                     | Air                     | Air                     |
| Connector Type                                                     | PM DB-25                | PM DB-25                | PM DB-25                | PM DB-25                |
| Cable Length (m)                                                   | 2                       | 2                       | 2                       | 2                       |
| Part Number                                                        | 1098350                 | 1098400                 | 1098413                 | 1098341                 |

PS10/PS10Q

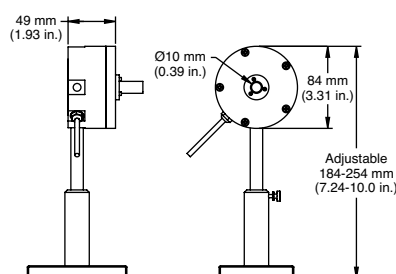

PS19/PS19Q

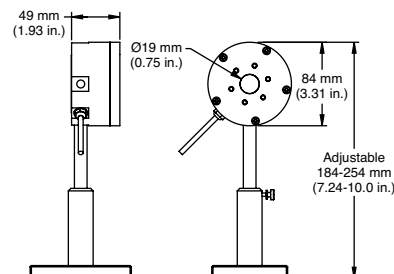

# Visit the Coherent Website

for more information about how Coherent can enable your ultrafast laser application.



[www.Coherent.com](http://www.Coherent.com)

**Use the Product and Application Finders below to quickly get more helpful details.**

---

Extensive Laser Finder

---



---

Powerful Diode Module Finder

---



---

Comprehensive Meter and Sensor Finder

---



---

Enhanced Application Finder

---



---

Brochures, Tech Notes, Technical Illustrations, Videos, Data Sheets

---



#### **Coherent, Inc., Corporate Headquarters**

5100 Patrick Henry Drive  
Santa Clara, CA 95054

#### **United States**

Coherent, Inc., Portland  
27650 SW 95th Avenue  
Wilsonville, OR 97070  
United States  
Toll Free: 1-800-343-4912  
Tel: (408) 764-4042  
Fax: (503) 454-5727  
[www.Coherent.com](http://www.Coherent.com)

#### **China**

Coherent, Inc., Beijing  
Room 215 Yin Gu Building  
No. 9 Beisihuan Xi Road  
Haidian District  
Beijing, 100080, China  
Tel: +86-10-6280-0209  
Fax: +86-10-6280-0129  
[www.Coherent.com.cn](http://www.Coherent.com.cn)

#### **France**

Coherent France  
Domaine Technologique de Saclay  
Bâtiment Azur  
4, Rue René Razel  
91892 Orsay Cedex  
France  
Tel: +33-1-69-85-51-45  
Fax: +33-1-60-19-40-00  
[www.Coherent.fr](http://www.Coherent.fr)

#### **Germany**

Coherent GmbH  
Dieselstraße 5b  
D-64807 Dieburg  
Germany  
Tel: +49-6071-968-333  
Fax: +49-6071-968-499  
[www.Coherent.de](http://www.Coherent.de)

#### **Italy**

Coherent Italia  
Via Borgese, 14  
Milano, 20154  
Italy  
Tel: +39-02-34-530-214  
Fax: +39-02-34-934-165  
[www.Coherent.it](http://www.Coherent.it)

#### **Japan**

Coherent Japan, Inc.  
Toyo MK Building  
7-2-14 Toyo  
Koto-ku, Tokyo 135-0016  
Japan  
Tel: +81 3-5635-8680  
Fax: +81 3-5635-8701  
[www.Coherent.co.jp](http://www.Coherent.co.jp)

#### **Korea**

Coherent Korea  
No. 101, Eagle Town  
278-20, 3 dong  
Seongsu 2 ga  
Seongdong-gu, Seoul  
133-832, South Korea  
Tel: +82 2 460 7900  
[www.Coherent.com](http://www.Coherent.com)

#### **Netherlands**

Coherent B.V.  
Smart Business Park  
Kanaalweg 18A  
3526 KL Utrecht  
Netherlands  
Tel: +31-30-280-6060  
Fax: +31-30-280-6077  
[www.Coherent.nl](http://www.Coherent.nl)

#### **United Kingdom**

Coherent UK Ltd.  
St Thomas' Place  
Cambridgeshire Business Park  
Ely  
CB7 4EX  
United Kingdom  
Tel: +44-1353-658-833  
Fax: +44-1353-659-110  
[www.Coherent.com](http://www.Coherent.com)